

**Quantitative Analysis for Ethanol & Qualitative Analysis for Other Volatiles**

*Analytical Method(s): 1.0*

*Device: Hamilton MICROLAB Liquid Processor/Dilutor Serial Number: ML600HC11378*

**Volatiles Quality Assurance Controls**

**Run Date(s): 07/14/2021 WORKLIST 5108**

**Calibration Date: 07/06/2021**

Control Level	Expiration	Lot #	Target Value	Acceptable Range	Overall Results
Level 1	Jul-23	1907006	0.0764	0.0688-0.0840	0.0712 g/100cc
					0.0741 g/100cc
					0.2063 g/100cc
Level 2	Jul-23	1907007	0.2170	0.1953-0.2387	g/100cc
					g/100cc
					g/100cc
Multi-Component mixture:				Lot # FN007101701	OK
Curve Fit:				Column 1 0.99962	Column2 0.99965

Ethanol Calibration Reference Material						
Calibrator level	Target Value	Acceptable Range	Column 1	Column 2	Precision	Mean
50	0.050	0.045 - 0.055	0.0517	0.0513	0.0004	0.0515
100	0.100	0.090 - 0.110	0.0986	0.0986	0	0.0986
200	0.200	0.180 - 0.220	0.2029	0.2032	0.0003	0.203
300	0.300	0.270 - 0.330	0.2944	0.2948	0.0004	0.2946
400	0.400	0.360 - 0.440				
500	0.500	0.450 - 0.550	0.5022	0.5019	0.0003	0.502

Aqueous Controls			
Control level	Target Value	Acceptable Range	Overall Results
80	0.080	0.076 - 0.084	0.081 g/100cc

**REVIEWED**  
 By Melissa (Nikka) Bradley at 1:55 pm, Jul 15, 2021

MB

The ALCOHOL.GCM method located in  
C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM was copied from  
C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM

7/15/21 GG

I double checked that both ALCOHOL.GCM in the following file paths have the same modification date and time in the computer:

C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM

MB 7/15/21



**Worklist: 5108**


<u>LAB_CASE</u>	<u>ITEM</u>	<u>ITEM_TYPE</u>	<u>DESCRIPTION</u>	
M2021-2949	1	BCK	Alcohol Analysis	
M2021-2950	1	BCK	Alcohol Analysis	
M2021-2955	1	BCK	Alcohol Analysis	
M2021-2956	1	BCK	Alcohol Analysis	
M2021-2956	2	BCK	Alcohol Analysis	
M2021-2956	3	BCK	Alcohol Analysis	
M2021-2991	1	BCK	Alcohol Analysis	
M2021-2992	1	BCK	Alcohol Analysis	
M2021-2993	1	BCK	Alcohol Analysis	
M2021-2994	1	BCK	Alcohol Analysis	
M2021-2995	1	BCK	Alcohol Analysis	
M2021-3002	1	BCK	Alcohol Analysis	
M2021-3027	1	BCK	Alcohol Analysis	
M2021-3028	1	BCK	Alcohol Analysis	
M2021-3032	1	BCK	Alcohol Analysis	
M2021-3036	1	BCK	Alcohol Analysis	
M2021-3054	1	BCK	Alcohol Analysis	
M2021-3055	1	BCK	Alcohol Analysis	
M2021-3070	1	BCK	Alcohol Analysis	



# Calibration Table


Laboratory : MERIDIAN  
 Instrument Name : GC-HS  
 Instrument Serial # : C12595800409 / C12255750548

<<Data File>>  
 Method File : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Batch File : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\CALCURVE\_TEMPLATE.gcb  
 Date Acquired : 7/6/2021 11:00:56 AM  
 Date Created : 7/6/2021 10:56:18 AM  
 Date Modified : 7/6/2021 11:03:58 AM

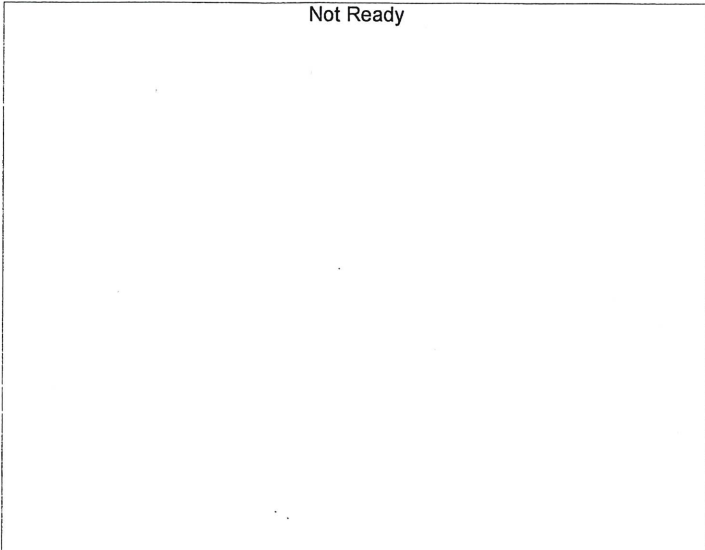

Name : Methanol  
 Detector Name: FID1  
 Function :  $f(x)=0*x+0$   
 R<sup>2</sup> value= 0  
 FitType: Linear  
 ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

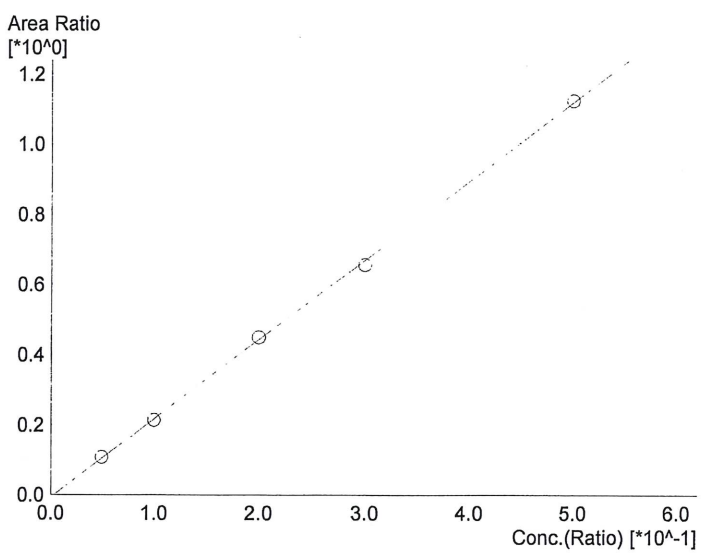

Name : Ethanol  
 Detector Name: FID1  
 Function :  $f(x)=2.25753*x-0.00881545$   
 R<sup>2</sup> value= 0.9996152  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
1	0.050	23858	0.0517
2	0.100	44073	0.0986
3	0.200	96088	0.2029
4	0.300	134135	0.2944
5	0.500	241028	0.5022

*W*




Name : Isopropyl Alcohol  
Detector Name: FID1  
Function :  $f(x)=0*x+0$   
R<sup>2</sup> value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

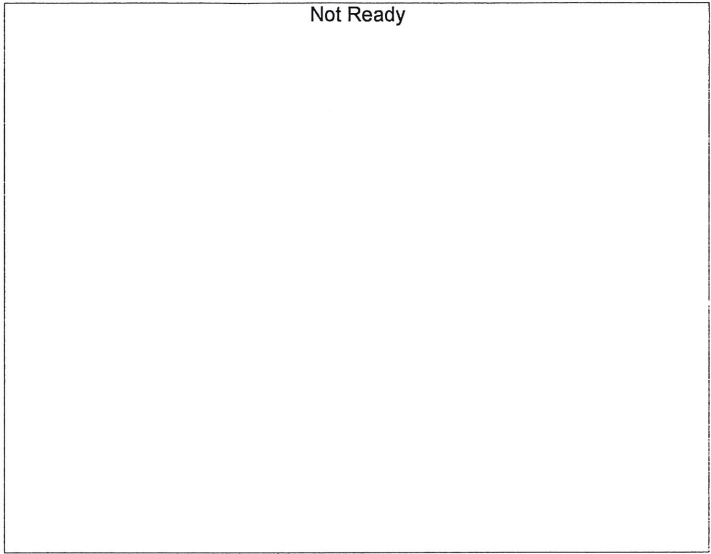


Name : Acetone  
Detector Name: FID1  
Function :  $f(x)=0*x+0$   
R<sup>2</sup> value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

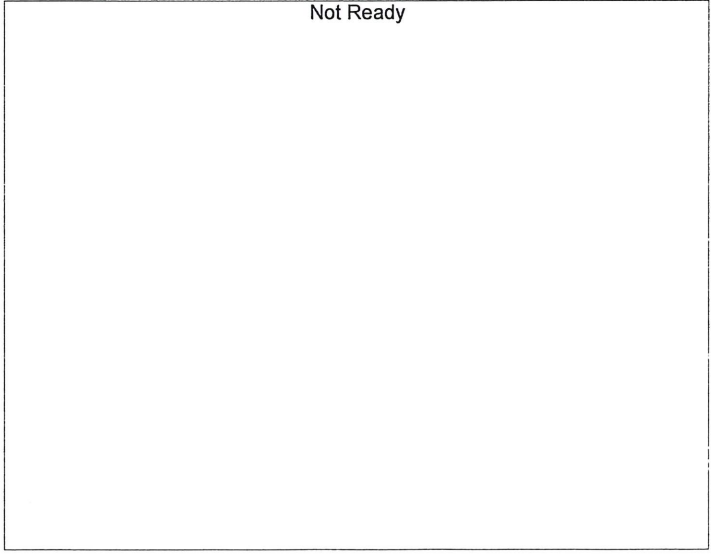


Name : Fluor. Hydrocarbon(s)  
Detector Name: FID1  
Function :  $f(x)=0*x+0$   
R<sup>2</sup> value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

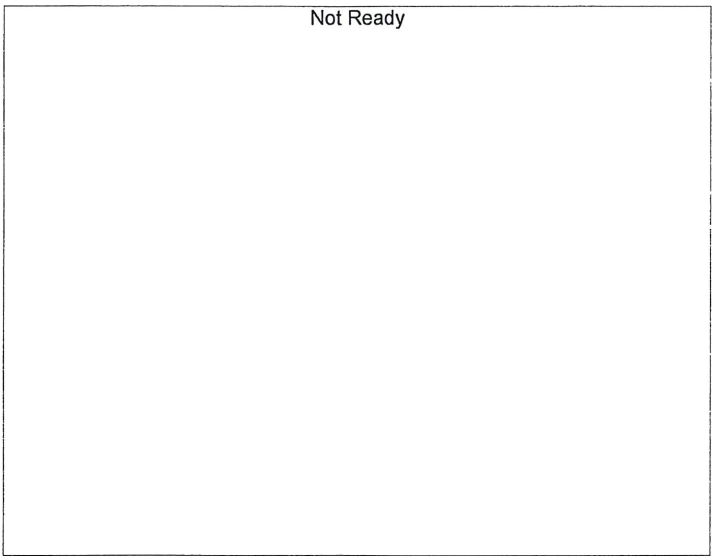


Name : Methanol  
 Detector Name: FID2  
 Function :  $f(x)=0*x+0$   
 R<sup>2</sup> value= 0  
 FitType: Linear  
 ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

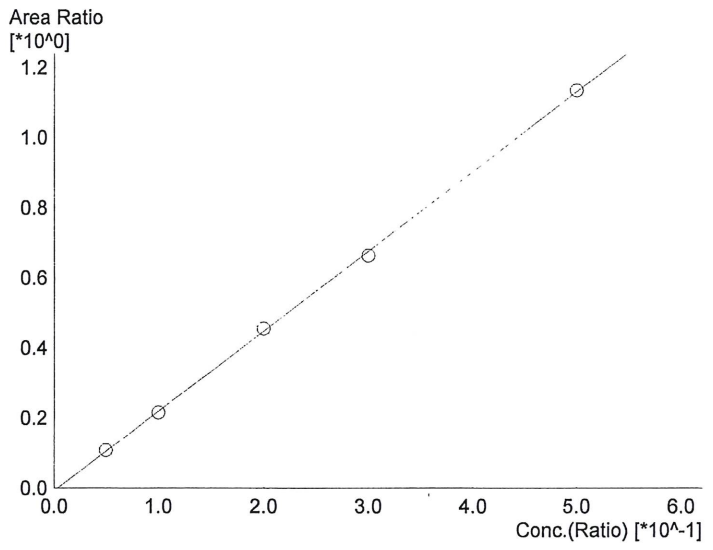

Name : Ethanol  
 Detector Name: FID2  
 Function :  $f(x)=2.27939*x-0.00951992$   
 R<sup>2</sup> value= 0.9996537  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
1	0.050	21210	0.0513
2	0.100	39685	0.0986
3	0.200	86786	0.2032
4	0.300	121227	0.2948
5	0.500	217990	0.5019

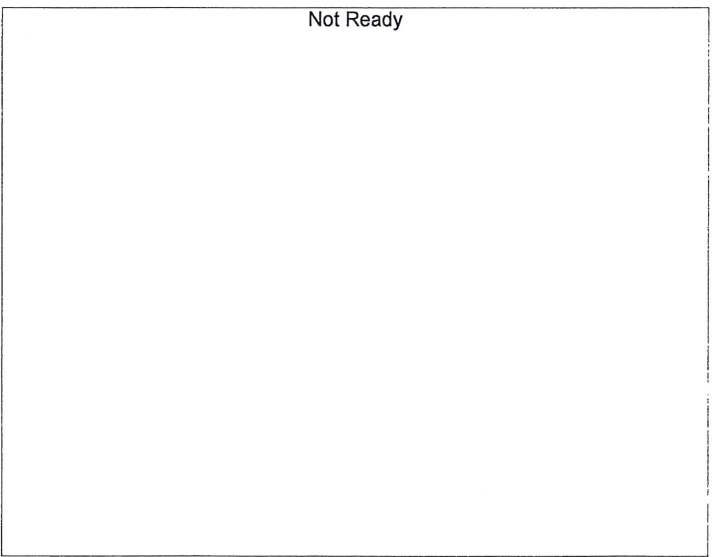

Name : Acetone  
 Detector Name: FID2  
 Function :  $f(x)=0*x+0$   
 R<sup>2</sup> value= 0  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
---	-------	------	------------

*W*

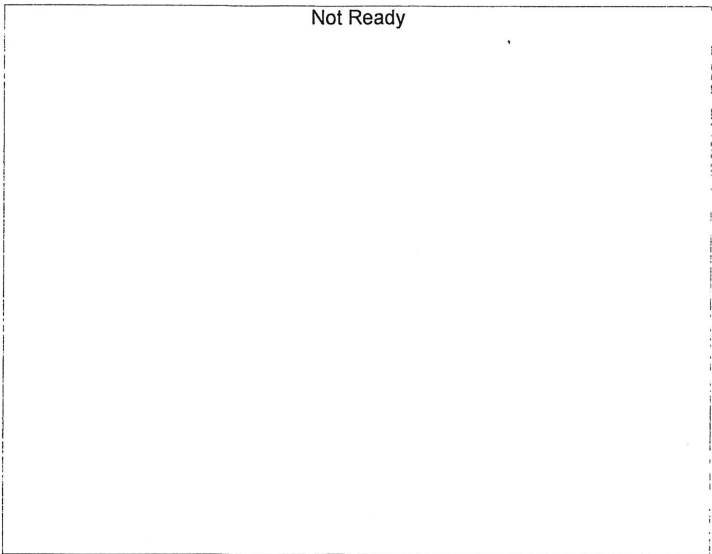

Name : Isopropyl Alcohol  
Detector Name: FID2  
Function :  $f(x)=0*x+0$   
R^2 value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

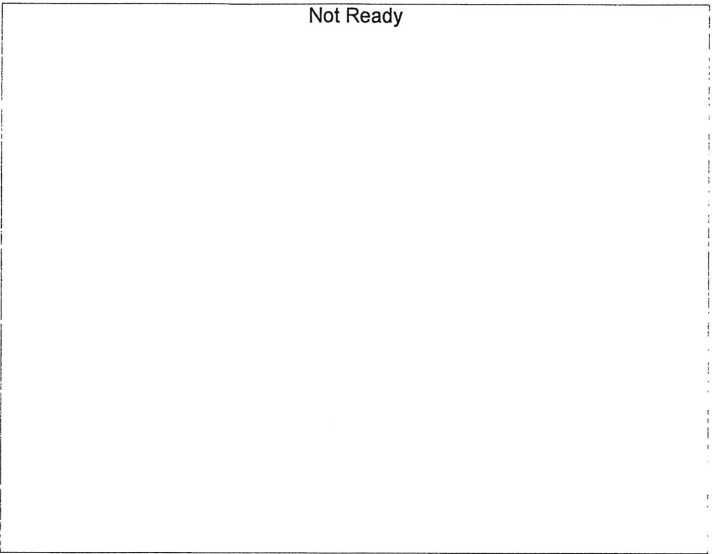

Name : Flour. Hydrocarbon(s)  
Detector Name: FID2  
Function :  $f(x)=0*x+0$   
R^2 value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------

W

# Meridian Blood Alcohol Analysis Batch Table



Shimadzu GC-2030 Serial #C12255750548  
Shimadzu HS-20 Serial #C12595800409  
Lab Solutions Software Ver. 5.99  
Copyright (C) 2008-2020 Shimadzu Corporation

Vial#	Sample Name	Sample Type	Level#	Method File
1	0.050	1:Standard:(I)	1	ALCOHOL.GCM
2	0.100	1:Standard	2	ALCOHOL.GCM
3	0.200	1:Standard	3	ALCOHOL.GCM
4	0.300	1:Standard	4	ALCOHOL.GCM
5	0.500	1:Standard	5	ALCOHOL.GCM
6	INT STD BLNK	0:Unknown	0	ALCOHOL.GCM





Sample Name : 0.050  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 10:29:46 AM  
 Vial # : 1  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

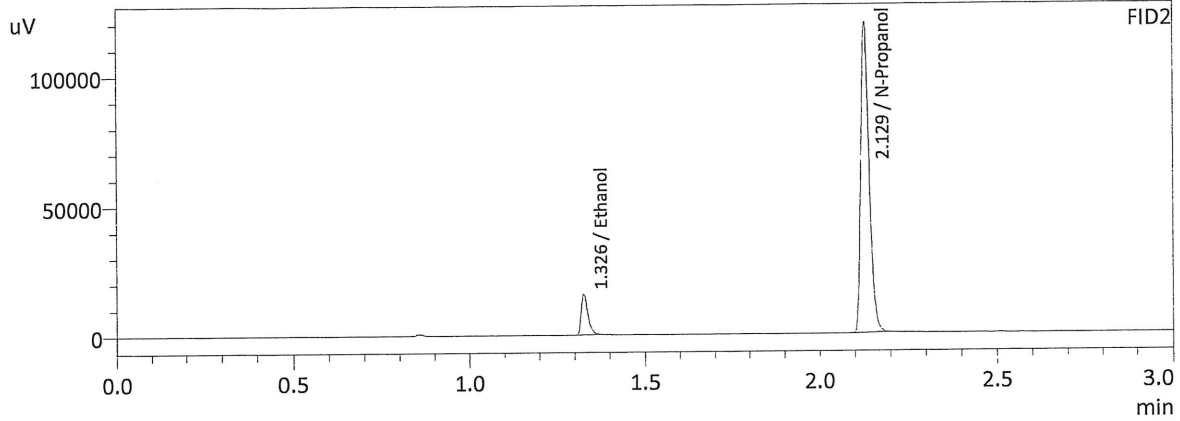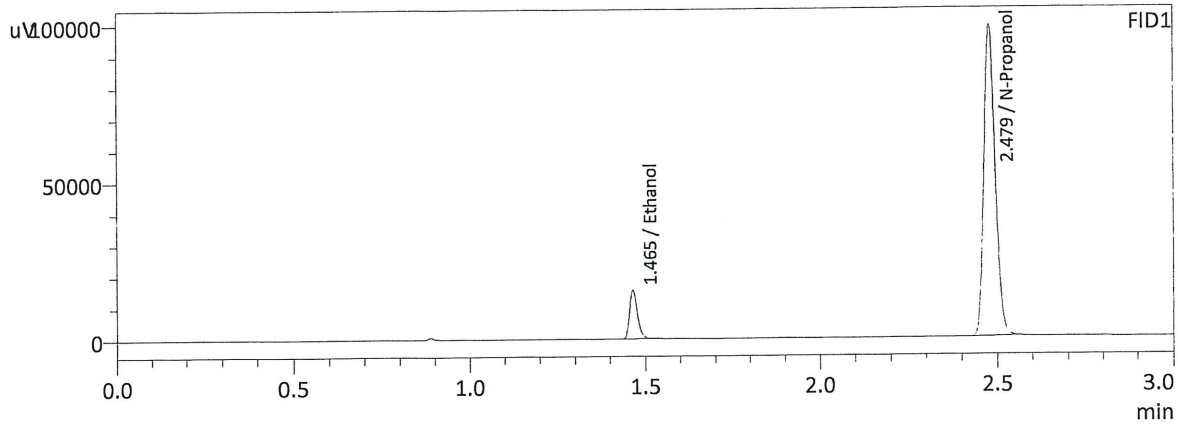

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0517	23858	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	221070	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0513	21210	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	197218	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : 0.100  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 10:37:06 AM  
 Vial # : 2  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

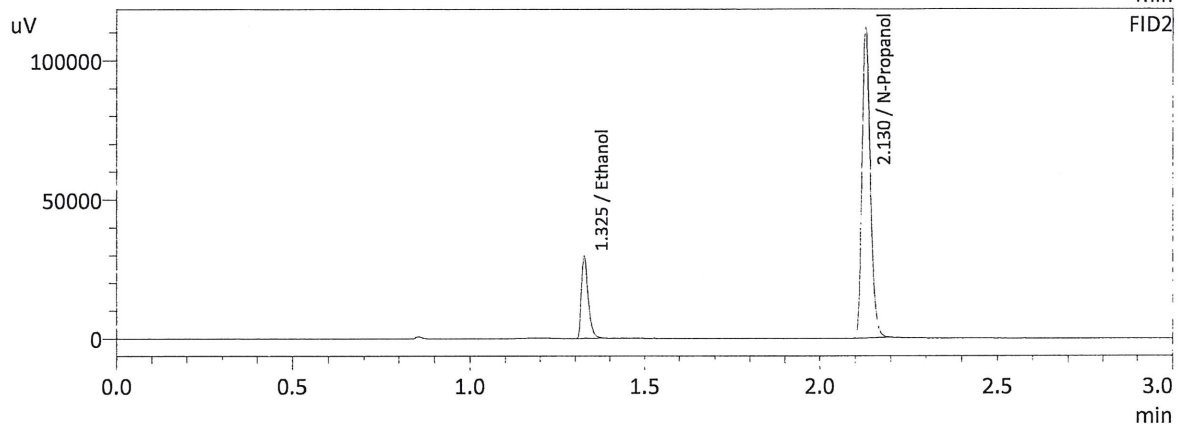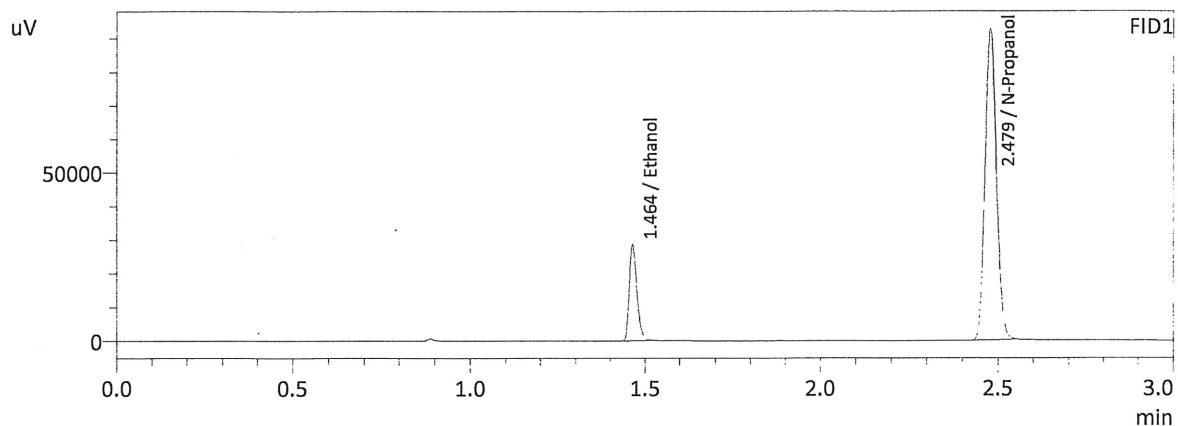

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0986	44073	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	206001	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0986	39685	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	184370	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : 0.200  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 10:44:43 AM  
 Vial # : 3  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

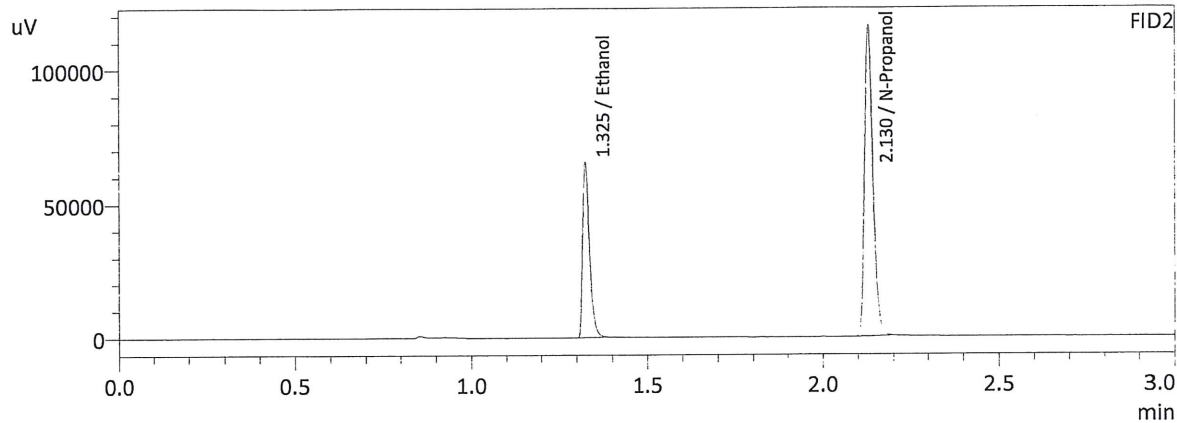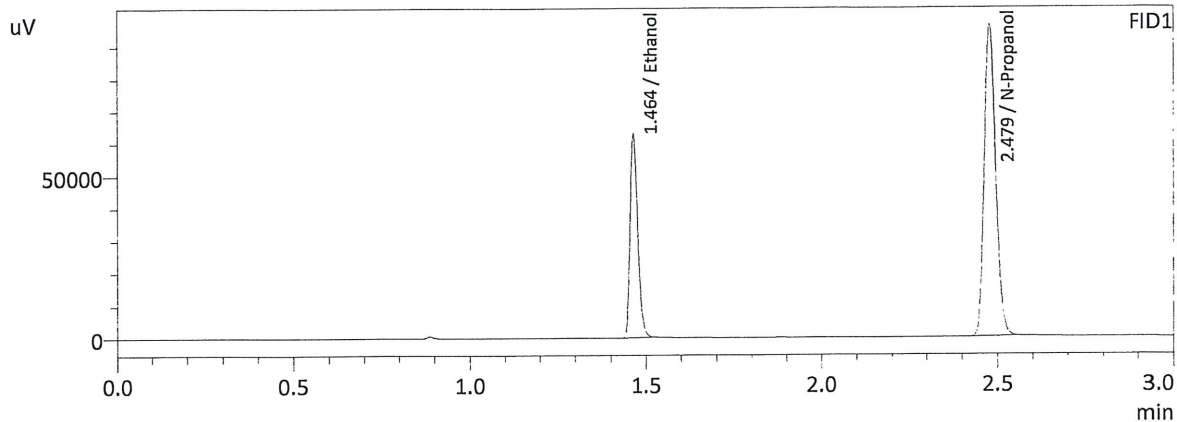

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2029	96088	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	213879	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2032	86786	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	191265	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

W

Sample Name : 0.300  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 10:53:09 AM  
 Vial # : 4  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

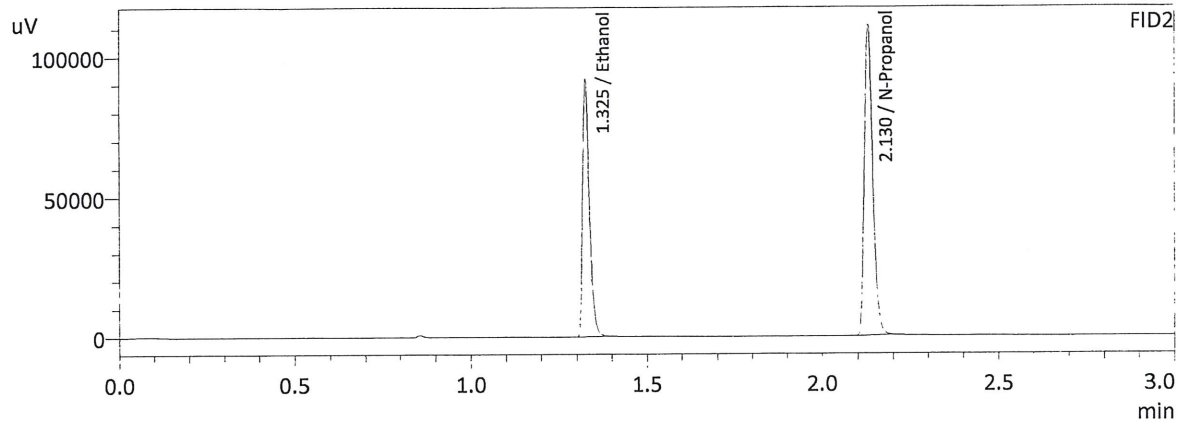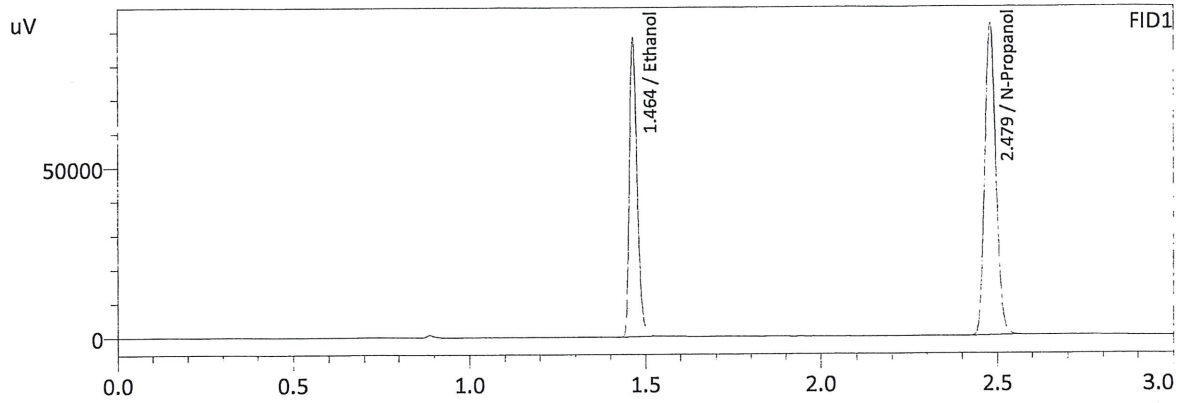

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2944	134135	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	204511	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2948	121227	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	182956	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : 0.500  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 11:00:56 AM  
 Vial # : 5  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

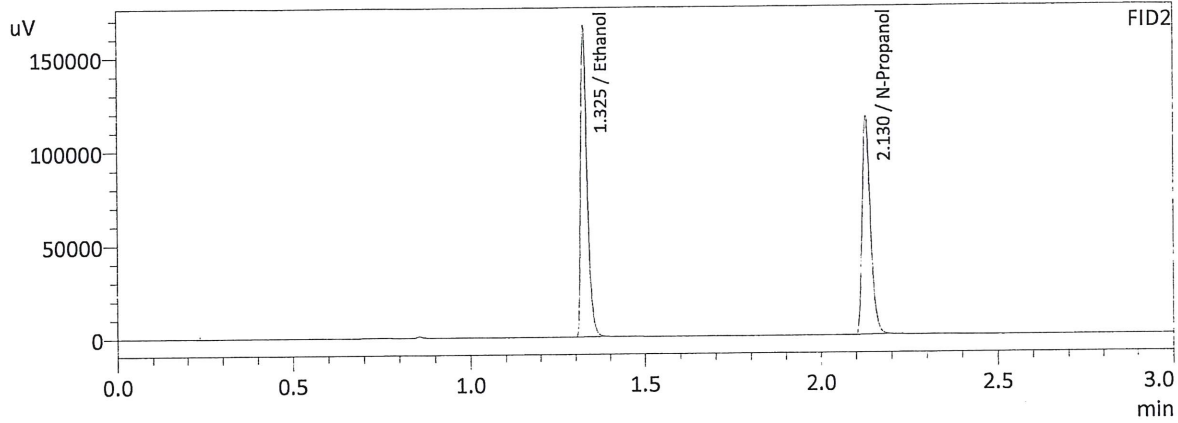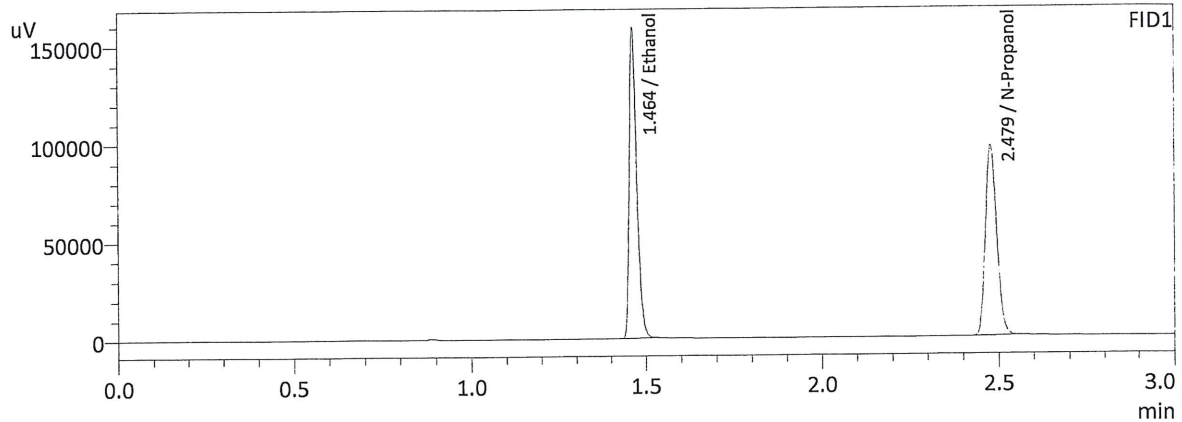

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.5022	241028	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	214234	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.5019	217990	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	192136	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : INT STD BLNK  
 Laboratory : Meridian  
 Injection Date : 7/6/2021 11:09:31 AM  
 Vial # : 6  
 Method Filename : C:\LabSolutions\Data\210706\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

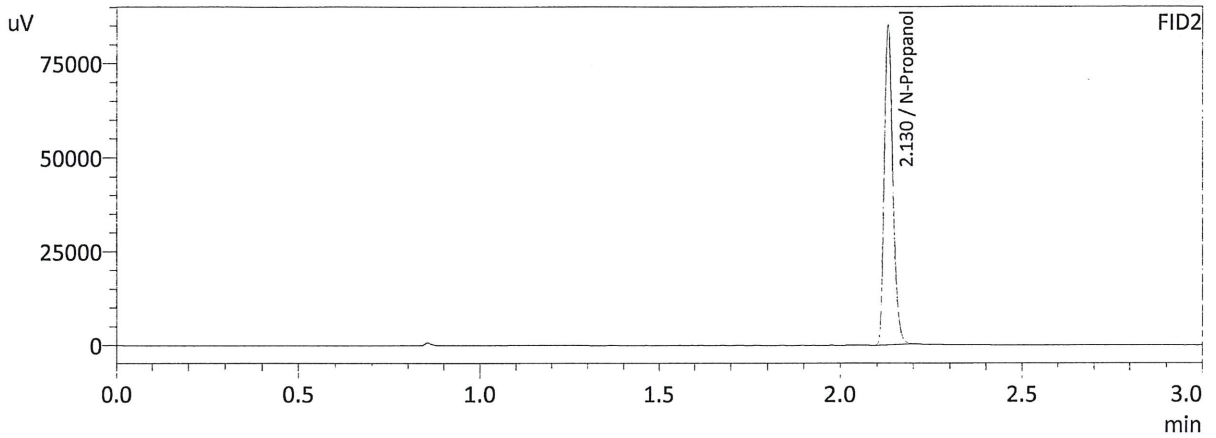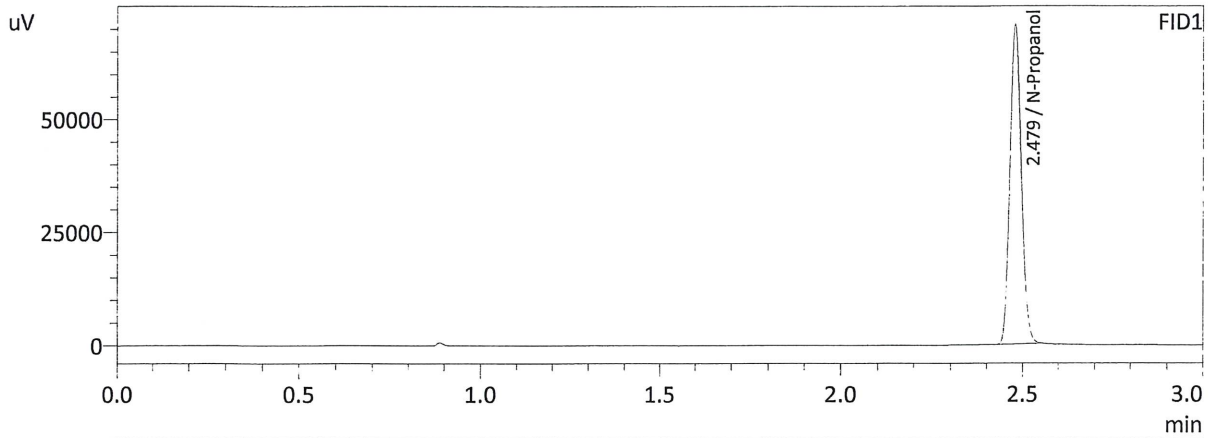

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	157495	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	140616	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc



*W*

# Meridian Blood Alcohol Analysis Batch Table

Shimadzu GC-2030 Serial #C12255750548  
 Shimadzu HS-20 Serial #C12595800409  
 Lab Solutions Software Ver. 5.99  
 Copyright (C) 2008-2020 Shimadzu Corporation

Vial#	Sample Name	Method File
1	INT STD BLK 1	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
2	ED VOLATILES FN 0710	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
3	QC-1-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
4	QC-1-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
5	0.08 QA-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
6	0.08 QA-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
7	M2021-2949-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
8	M2021-2949-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
9	M2021-2950-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
10	M2021-2950-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
11	M2021-2955-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
12	M2021-2955-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
13	M2021-2956-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
14	M2021-2956-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
15	M2021-2956-2-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
16	M2021-2956-2-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
17	M2021-2956-3-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
18	M2021-2956-3-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
19	M2021-2991-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
20	M2021-2991-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
21	M2021-2992-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
22	M2021-2992-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
23	M2021-2993-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
24	M2021-2993-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
25	QC-2-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
26	QC-2-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
27	M2021-2994-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
28	M2021-2994-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
29	M2021-2995-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
30	M2021-2995-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
31	M2021-3002-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
32	M2021-3002-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
33	M2021-3027-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
34	M2021-3027-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
35	M2021-3028-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
36	M2021-3028-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
37	M2021-3032-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
38	M2021-3032-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
39	M2021-3036-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
40	M2021-3036-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
41	M2021-3054-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
42	M2021-3054-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
43	M2021-3055-1-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
44	M2021-3055-1--B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
45	M2021-3070-1--A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
46	M2021-3070-1-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
47	QC1-2-A	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
48	QC1-2-B	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
49	INT STD BLK 2	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
50	DFE 111914OM	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
51	INT STD BLK 3	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
52	TFE 111914	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL
53	INT STD BLNK 4	Solutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL

Sample Name : INT STD BLK 1  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:10:53 PM  
 Vial # : 1  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

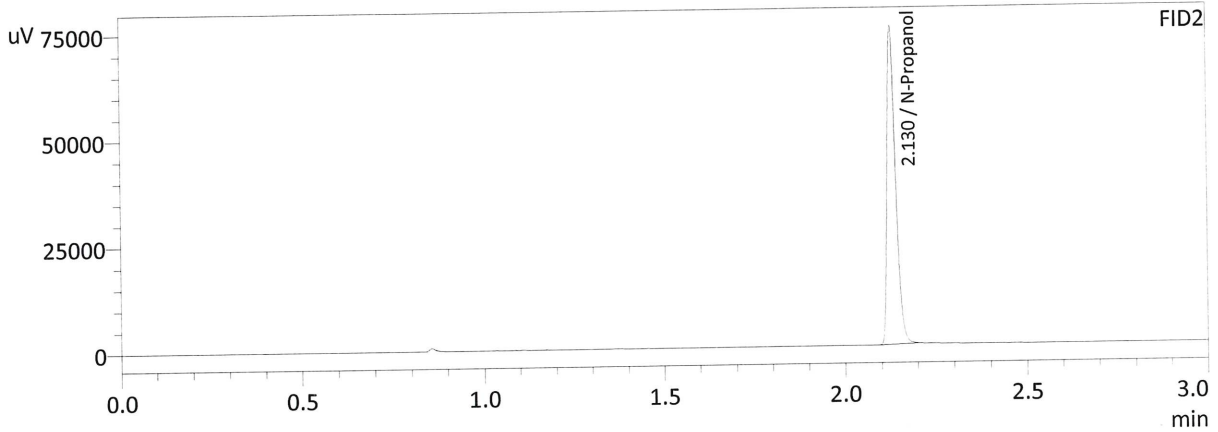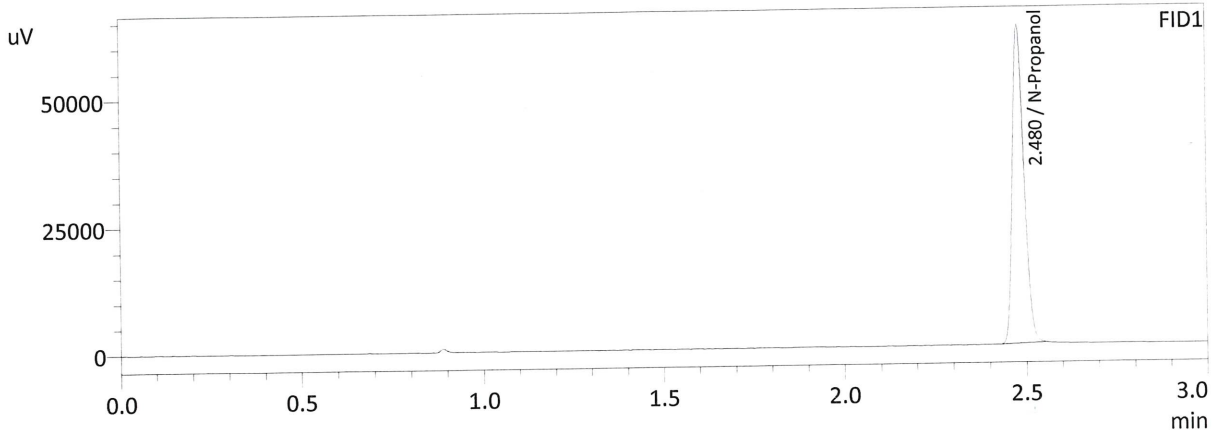

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	139508	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc



FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	124060	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*



Sample Name : MIXED VOLATILES FN 07101701  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:18:12 PM  
 Vial # : 2  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

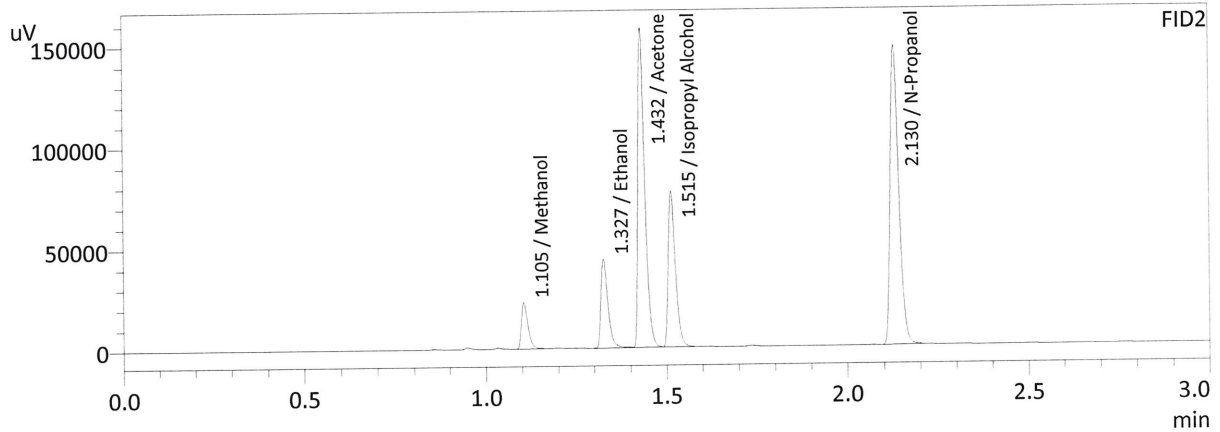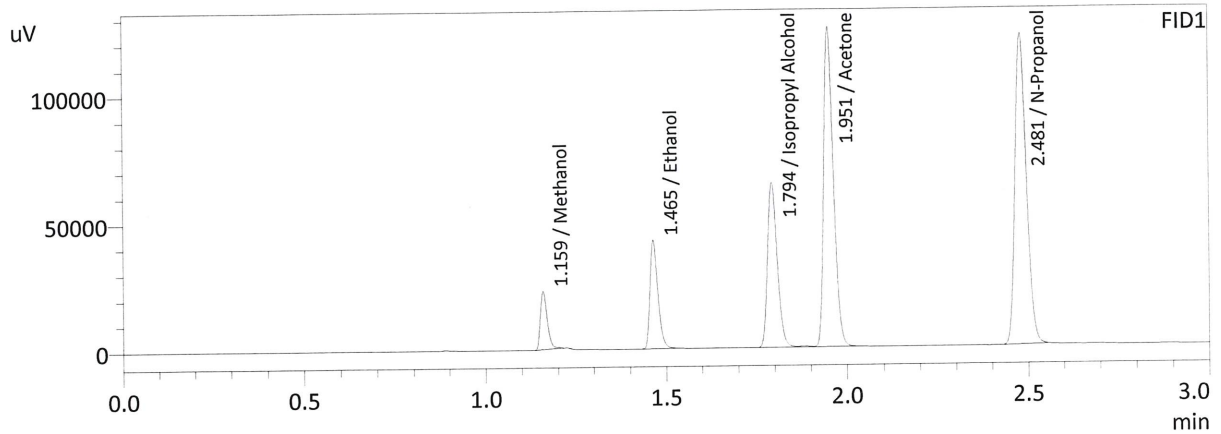

FID1

Name	Conc.	Area	Unit
Methanol	0.0000	30444	g/100cc
Ethanol	0.1098	64707	g/100cc
Isopropyl Alcohol	0.0000	119698	g/100cc
Acetone	0.0000	232127	g/100cc
N-Propanol	0.0000	270650	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	0.0000	28646	g/100cc
Ethanol	0.1115	59154	g/100cc
Acetone	0.0000	211991	g/100cc
Isopropyl Alcohol	0.0000	107820	g/100cc
N-Propanol	0.0000	241679	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

W

**VOLATILES DETERMINATION CASEFILE WORKSHEET**

Laboratory No.: QC 1-1

Analysis Date(s): 07/14/2021

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0720	0.0720	0.0000	0.0720	0.0016	0.0712
(g/100cc)	0.0705	0.0704	0.0001	0.0704		

**Analysis Method**

Refer to Blood Alcohol Method #1

**Instrument Information**

*Instrument information is stored centrally.*

Refer to Instrument Method: Alcohol.m/.gcm, Volatiles.m/.gcm

**Reporting of Results**

Uncertainty of Measurement (UM%): 5.00%


Overall Mean (g/100cc)	Low	High	5% of Mean
0.071	0.067	0.075	0.004

Reported Result	
0.071	

*Calibration and control data are stored centrally.*

*W*

Sample Name : QC-1-1-A  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:25:33 PM  
 Vial # : 3  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

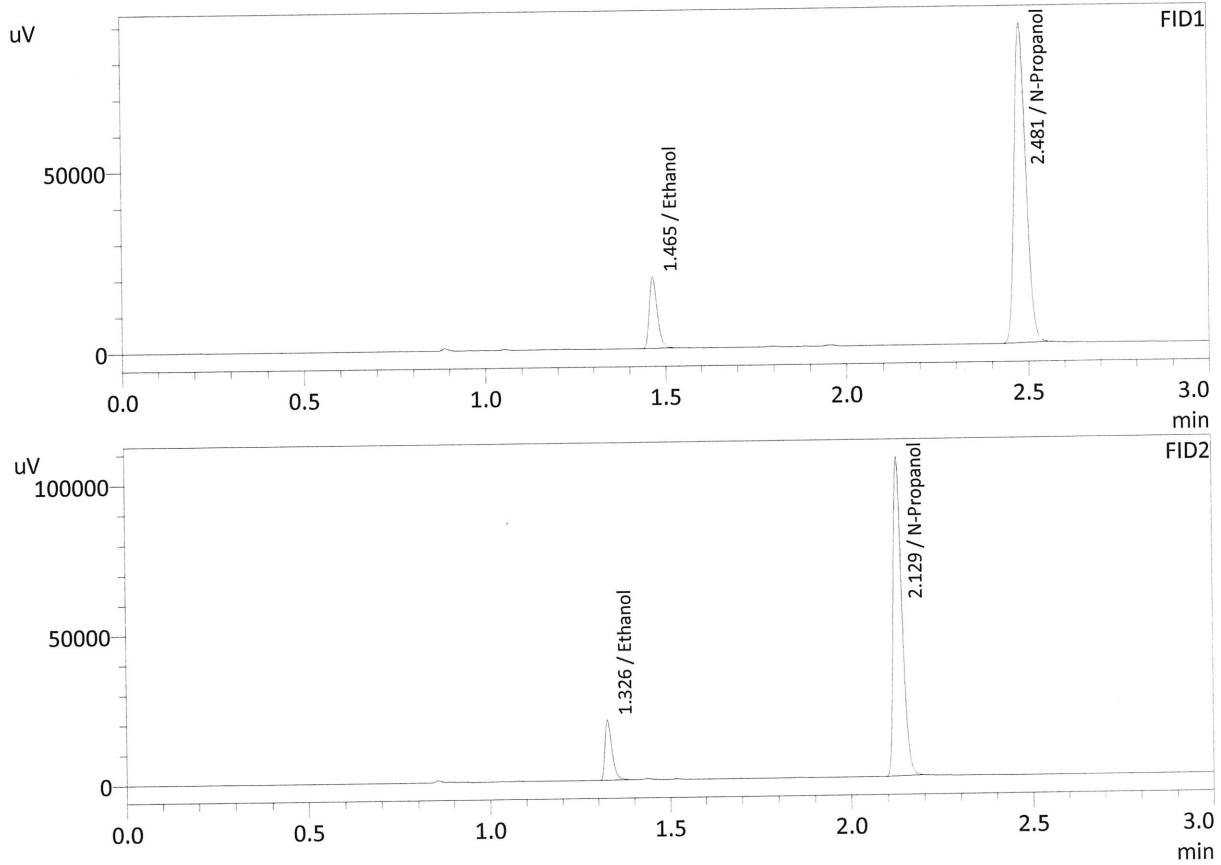

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0720	30210	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	196253	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0720	27096	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	175072	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

W

Sample Name : QC-1-1-B  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:34:26 PM  
 Vial # : 4  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

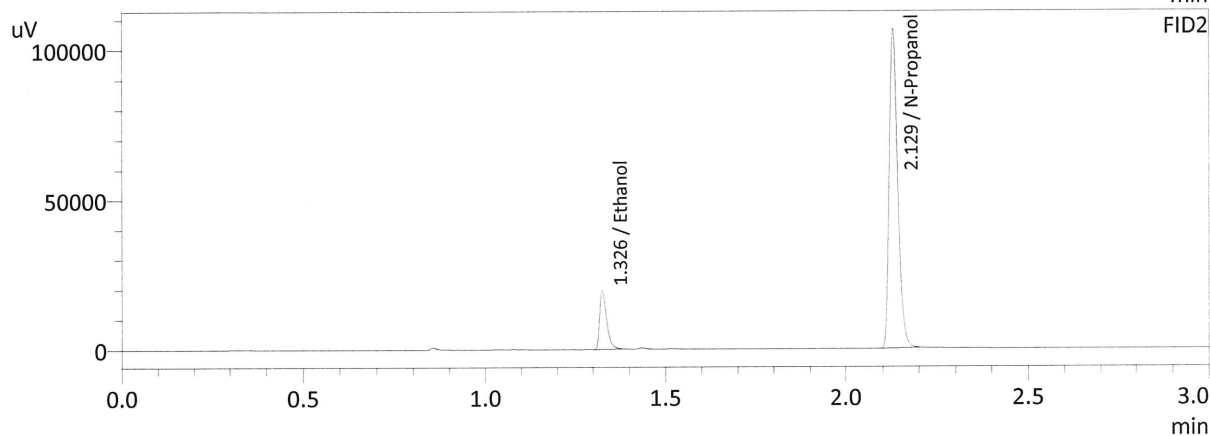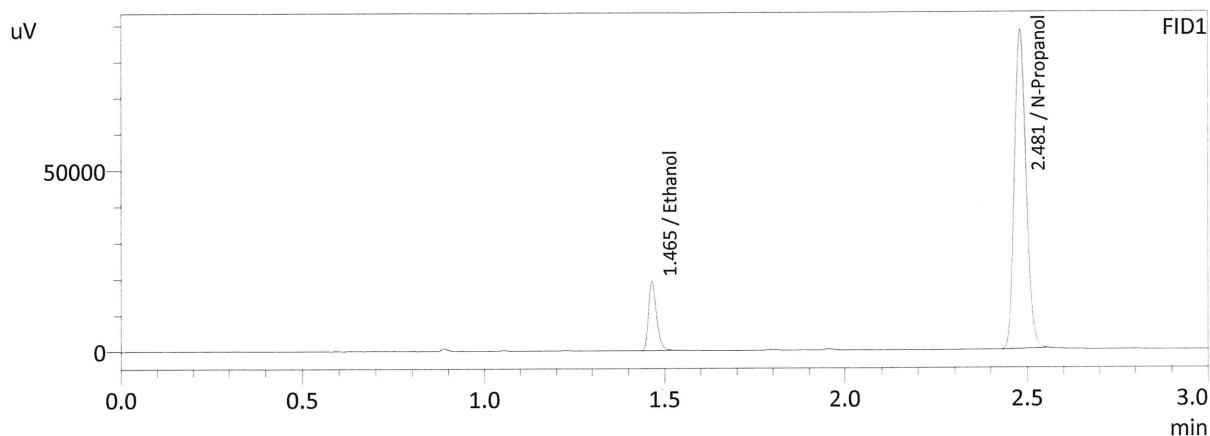

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0705	29508	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	196228	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0704	26476	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	175240	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

W

**VOLATILES DETERMINATION CASEFILE WORKSHEET**

Laboratory No.: QC 1-2

Analysis Date(s): 07/14/2021

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0735	0.0738	0.0003	0.0736	0.0009	0.0741
(g/100cc)	0.0744	0.0747	0.0003	0.0745		

**Analysis Method**

Refer to Blood Alcohol Method #1

**Instrument Information**

*Instrument information is stored centrally.*

Refer to Instrument Method: Alcohol.m/.gcm, Volatiles.m/.gcm

**Reporting of Results**



Uncertainty of Measurement (UM%): 5.00%

Overall Mean (g/100cc)	Low	High	5% of Mean
0.074	0.070	0.078	0.004

Reported Result	
0.074	

*Calibration and control data are stored centrally.*

Sample Name : QC1-2-A  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 8:20:26 PM  
 Vial # : 47  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

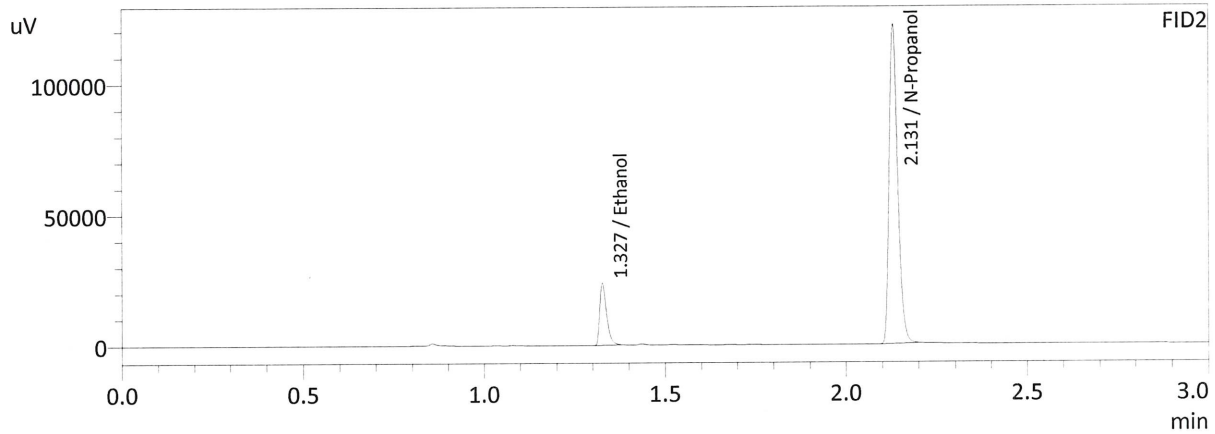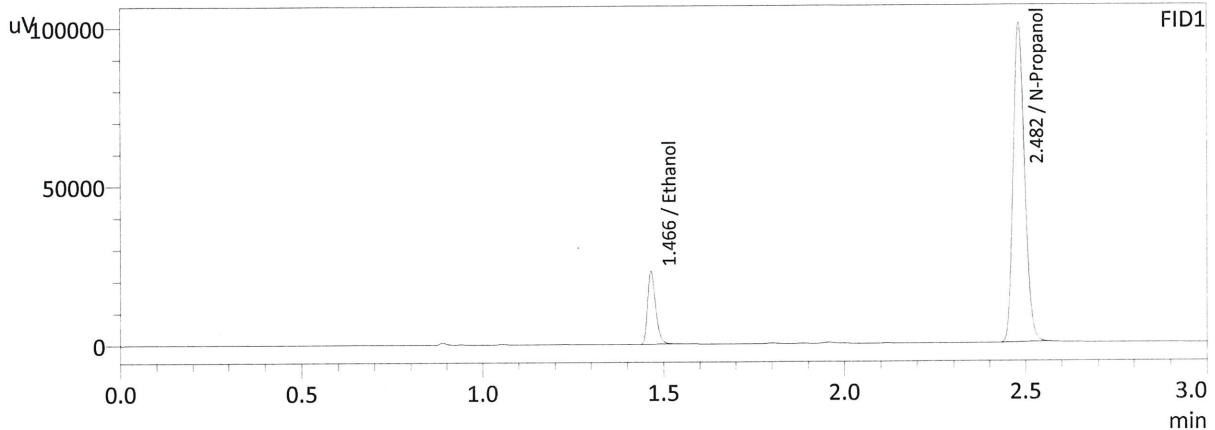

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0735	35204	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	223931	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0738	31790	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	200229	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : QC1-2-B  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 8:29:10 PM  
 Vial # : 48  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

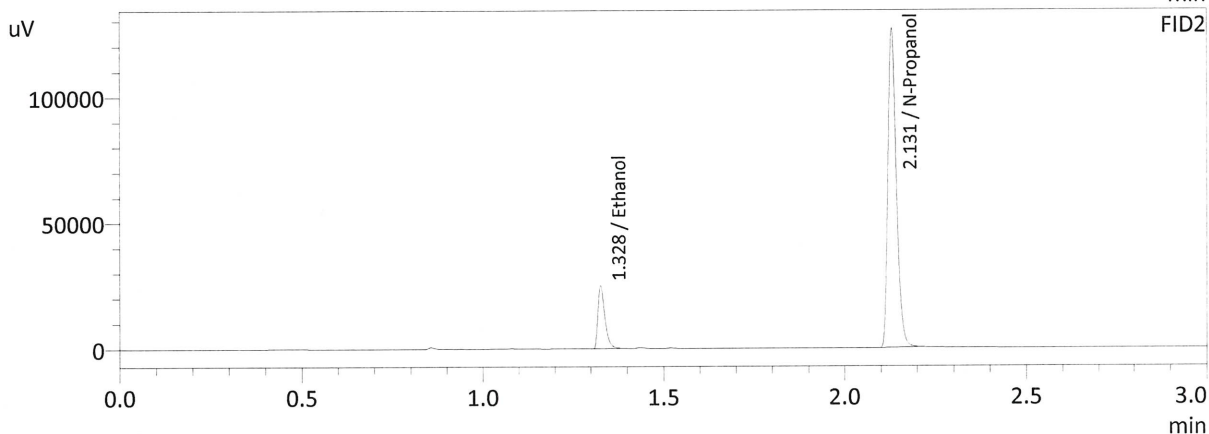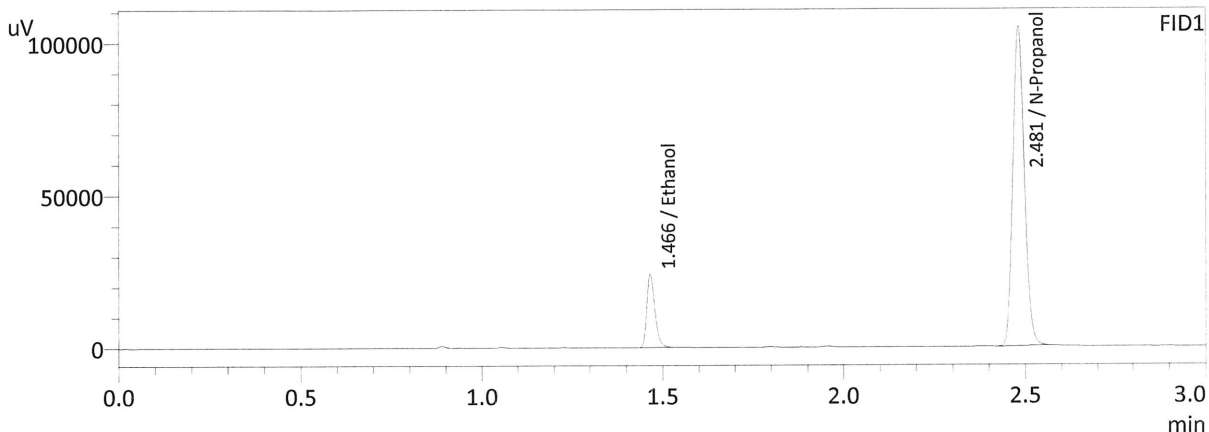

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0744	37090	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	232822	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0747	33479	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	208067	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

W

**VOLATILES DETERMINATION CASEFILE WORKSHEET**

Laboratory No.: QC 2-1

Analysis Date(s): 07/14/2021

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.2077	0.2092	0.0015	0.2084	0.0042	0.2063
(g/100cc)	0.2038	0.2047	0.0009	0.2042		

**Analysis Method**

Refer to Blood Alcohol Method #1

**Instrument Information**

*Instrument information is stored centrally.*

Refer to Instrument Method: Alcohol.m/.gcm, Volatiles.m/.gcm

**Reporting of Results**

**Uncertainty of Measurement (UM%): 5.00%**



Overall Mean (g/100cc)	Low	High	5% of Mean
0.206	0.195	0.217	0.011

Reported Result	
0.206	

*Calibration and control data are stored centrally.*



Sample Name : QC-2-1-A  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 5:23:41 PM  
 Vial # : 25  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

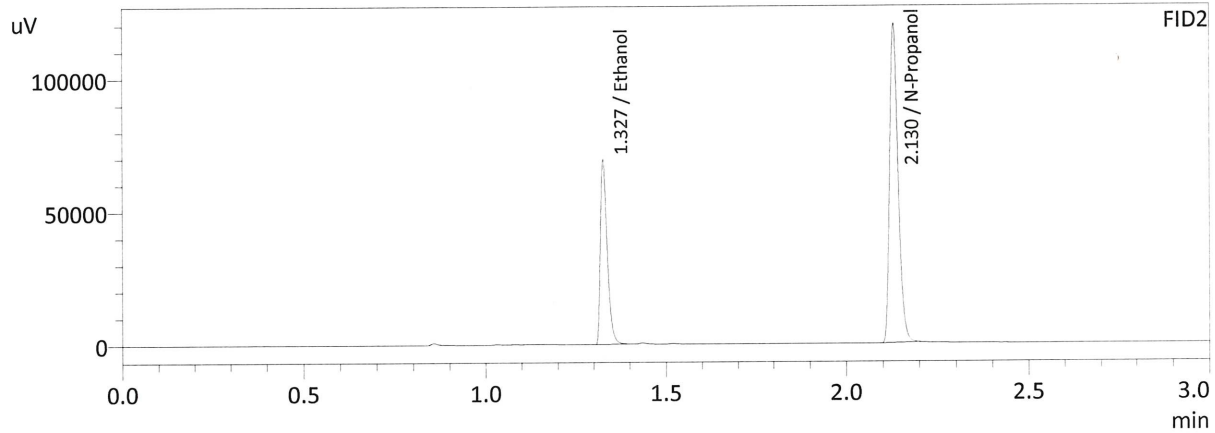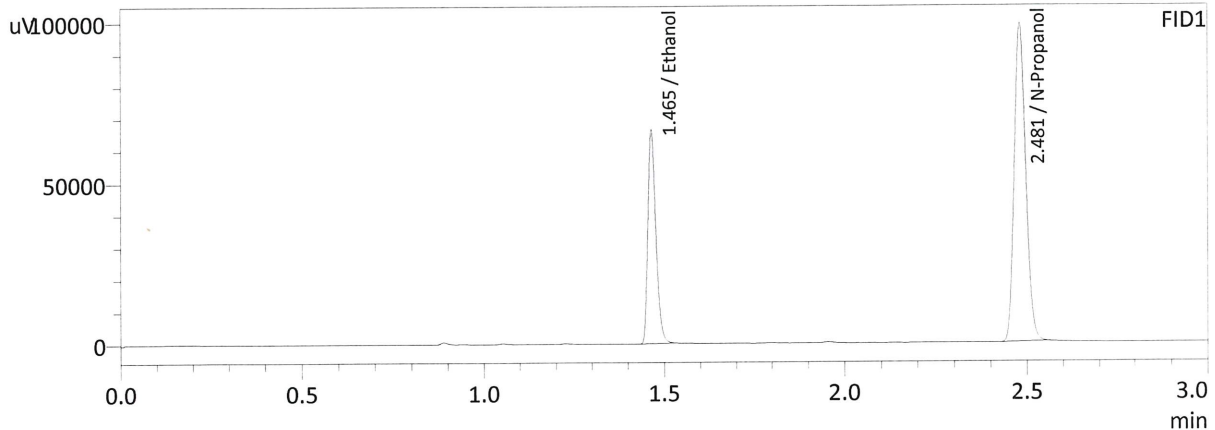

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2077	101322	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	220196	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2092	92051	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	196880	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : QC-2-1-B  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 5:30:58 PM  
 Vial # : 26  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

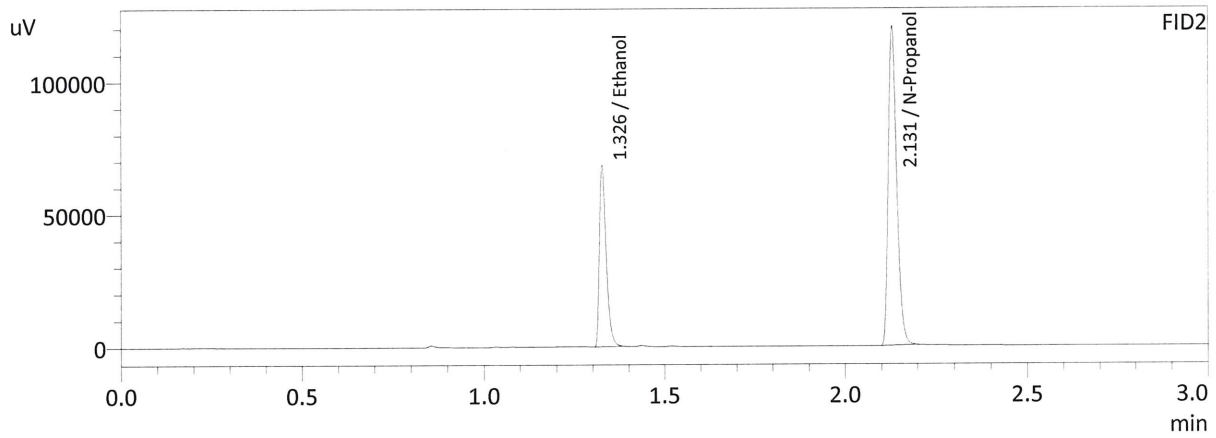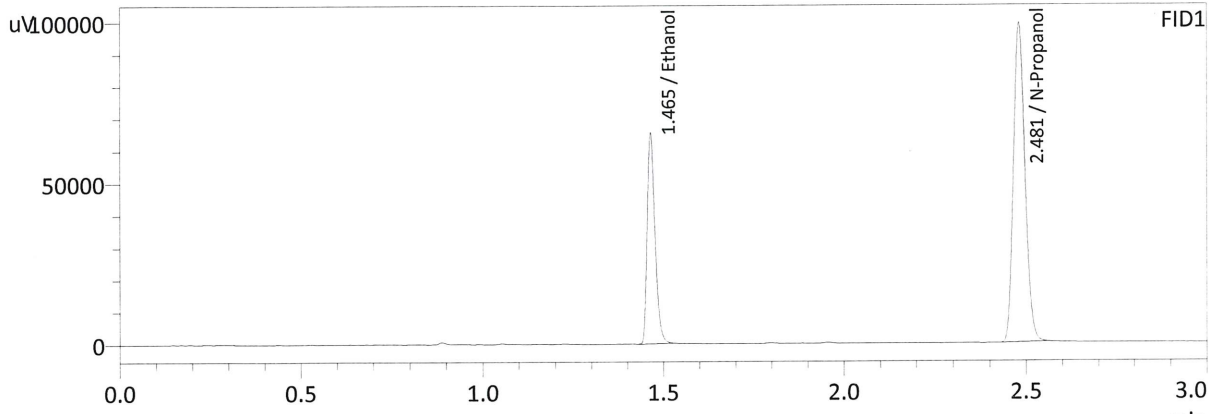

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2038	99790	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	221028	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2047	90346	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	197616	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

**VOLATILES DETERMINATION CASEFILE WORKSHEET**

Laboratory No.: QA 0.08

Analysis Date(s): 07/14/2021

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0816	0.0817	0.0001	0.0816	0.0002	0.0815
(g/100cc)	0.0813	0.0816	0.0003	0.0814		

**Analysis Method**

Refer to Blood Alcohol Method #1

**Instrument Information**

*Instrument information is stored centrally.*

Refer to Instrument Method: Alcohol.m/.gcm, Volatiles.m/.gcm

**Reporting of Results**



**Uncertainty of Measurement (UM%): 5.00%**

Overall Mean (g/100cc)	Low	High	5% of Mean
0.081	0.076	0.086	0.005

Reported Result	
0.081	

*Calibration and control data are stored centrally.*

Sample Name : 0.08 QA-A  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:42:09 PM  
 Vial # : 5  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

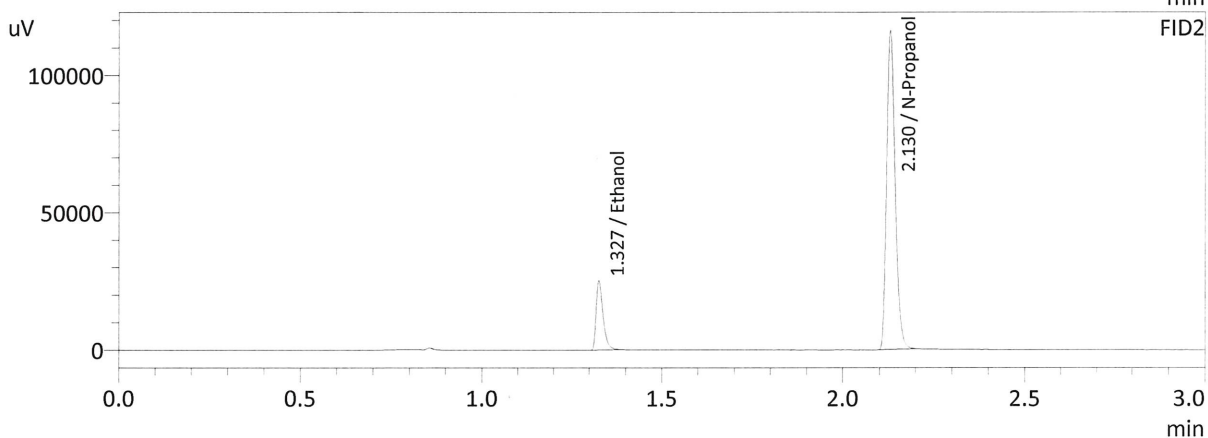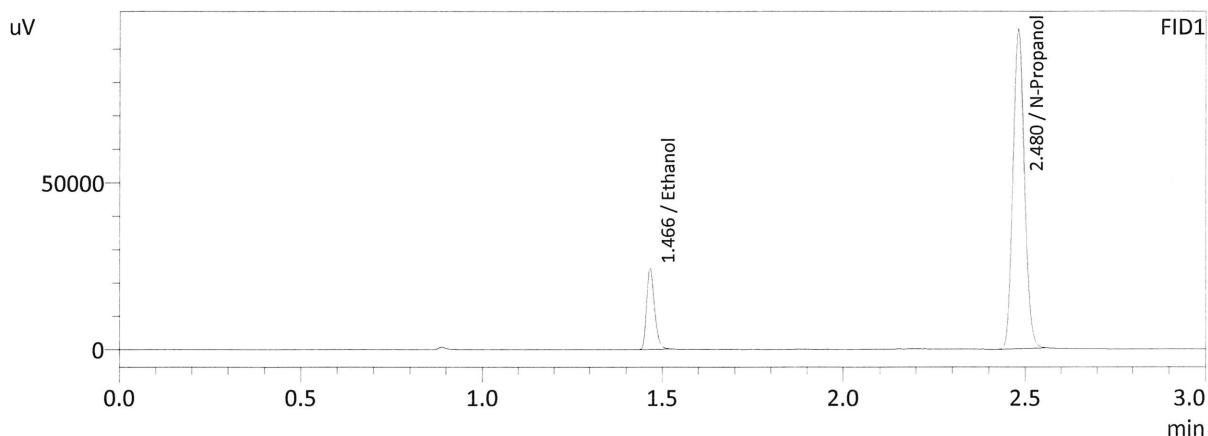

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0816	37448	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	213320	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0817	33725	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	190797	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : 0.08 QA-B  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 2:50:31 PM  
 Vial # : 6  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

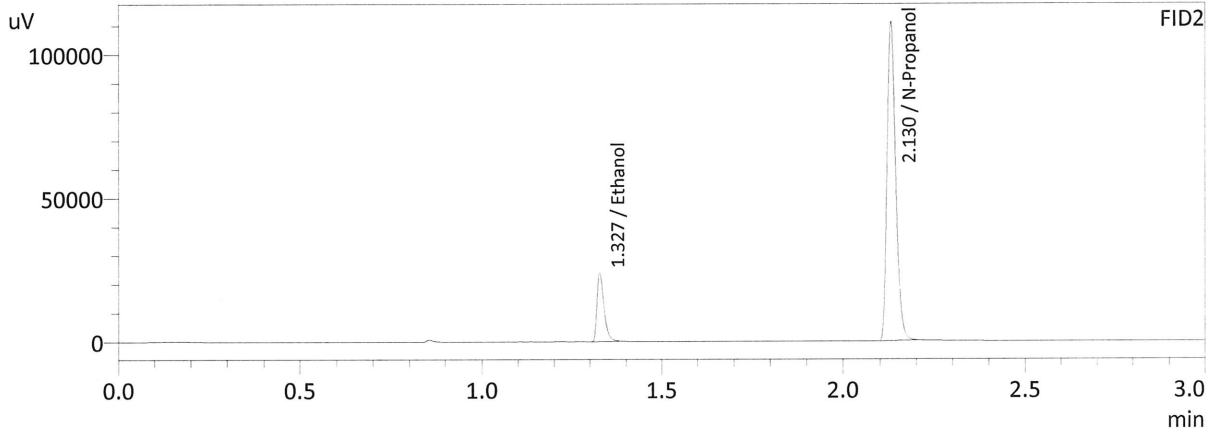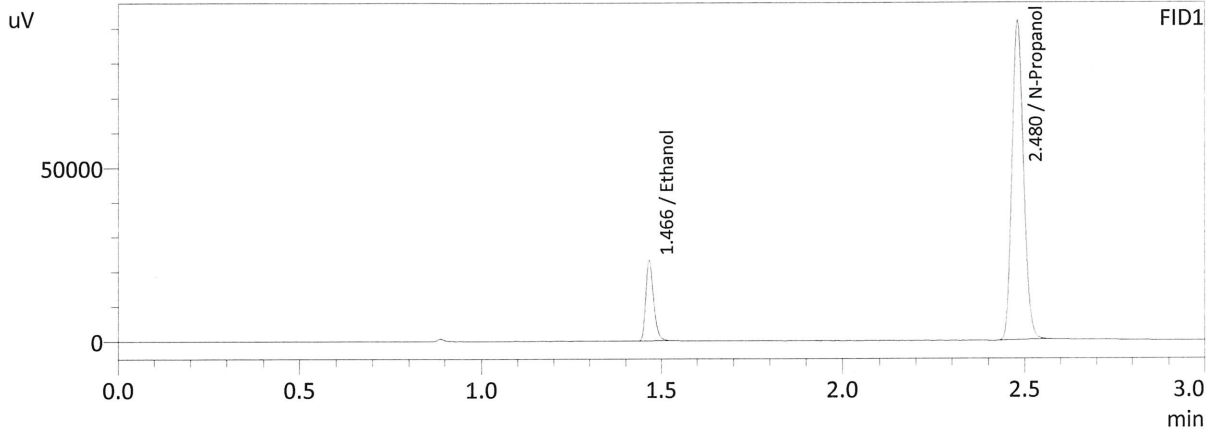

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0813	35782	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	204547	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0816	32254	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	182724	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*66*

Sample Name : INT STD BLK 2  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 8:36:38 PM  
 Vial # : 49  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

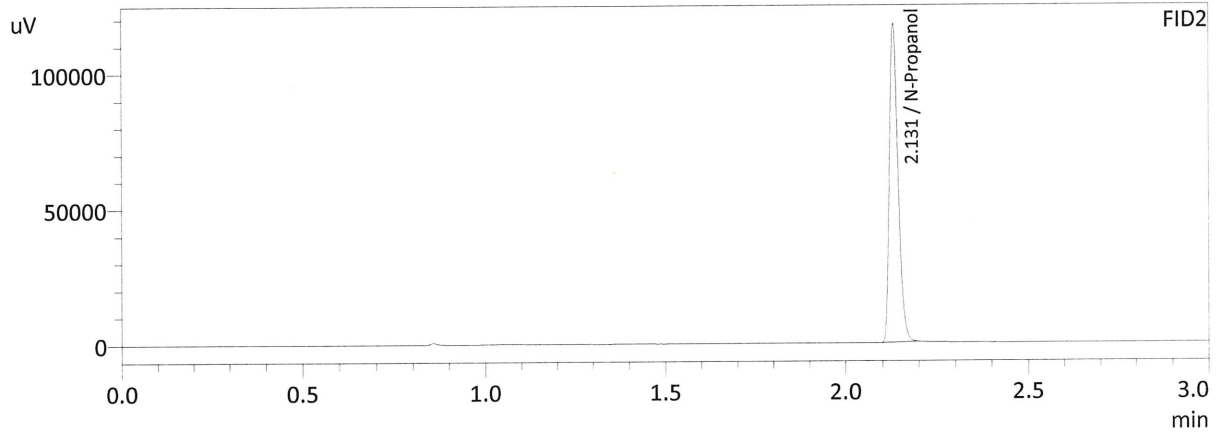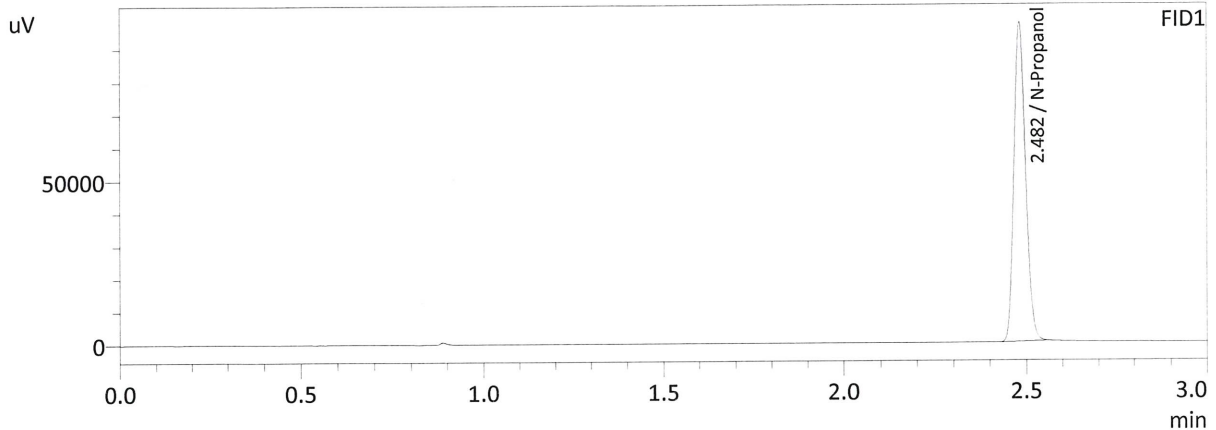

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	216523	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	193321	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

Sample Name : DFE 111914OM  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 8:45:10 PM  
 Vial # : 50  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

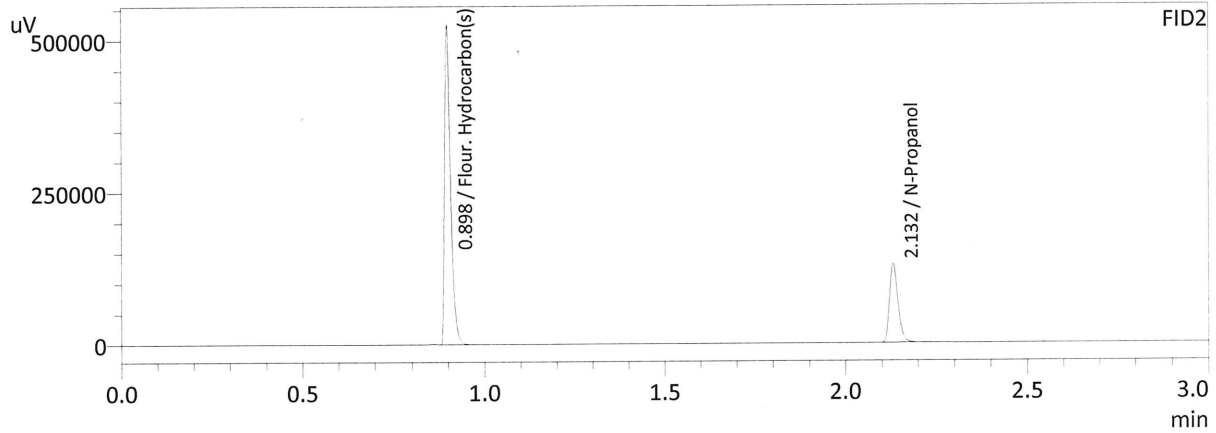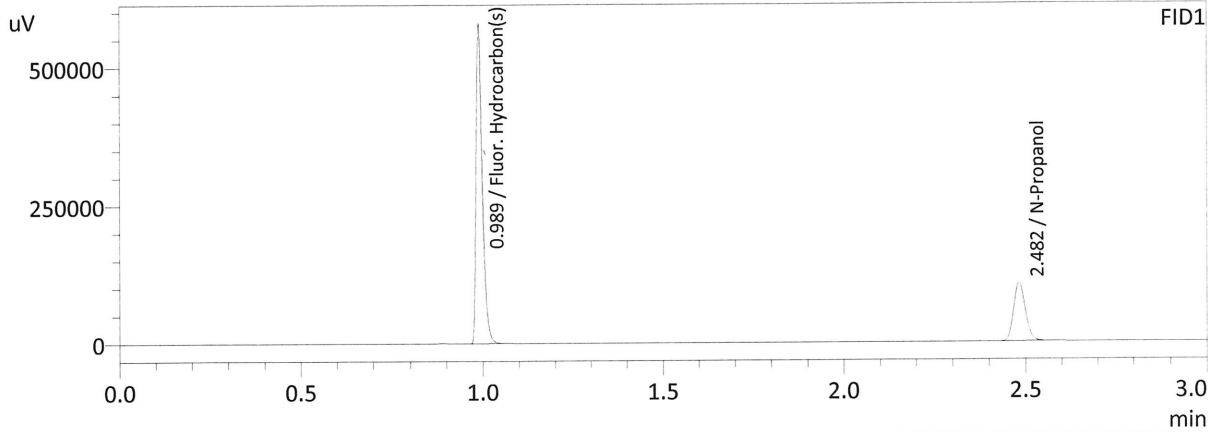

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	237050	g/100cc
Fluor. Hydrocarbon(s)	0.0000	688410	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	212324	g/100cc
Fluor. Hydrocarbon(s)	0.0000	605787	g/100cc

*W*

Sample Name : INT STD BLK 3  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 8:53:55 PM  
 Vial # : 51  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

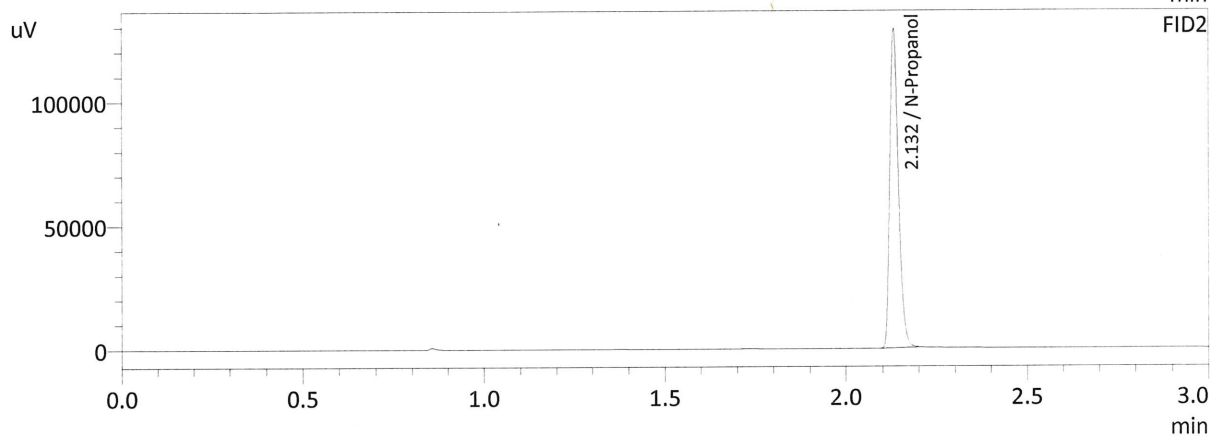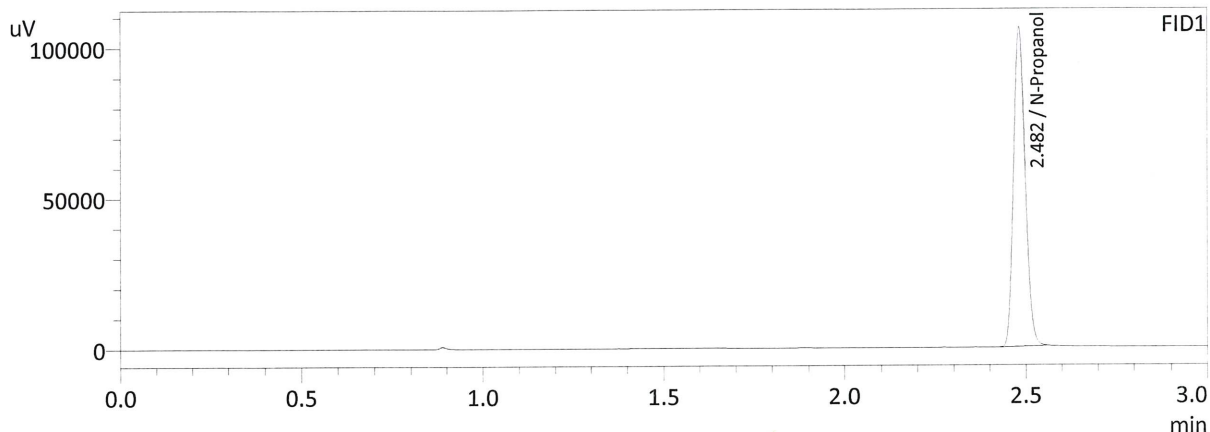

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	236024	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc



FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	211620	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*



Sample Name : TFE 111914  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 9:01:15 PM  
 Vial # : 52  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

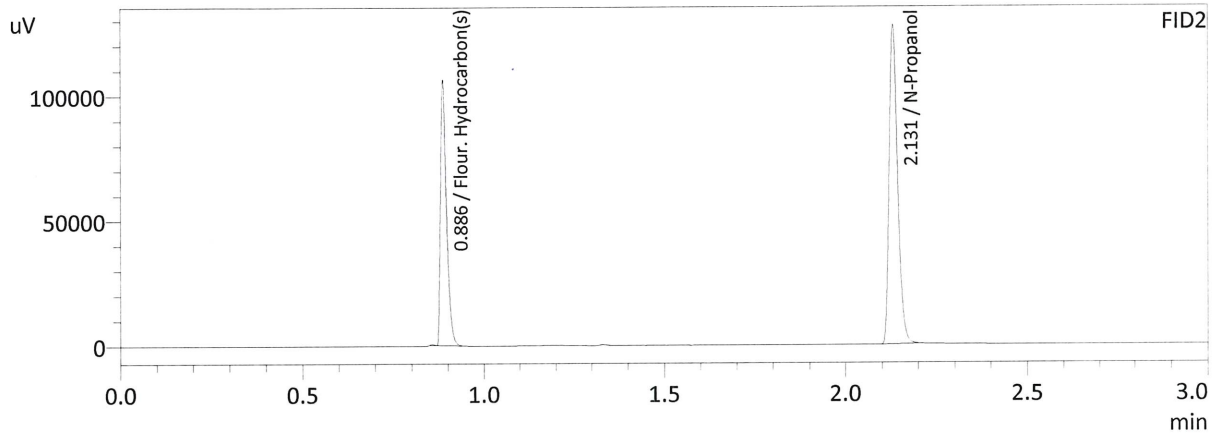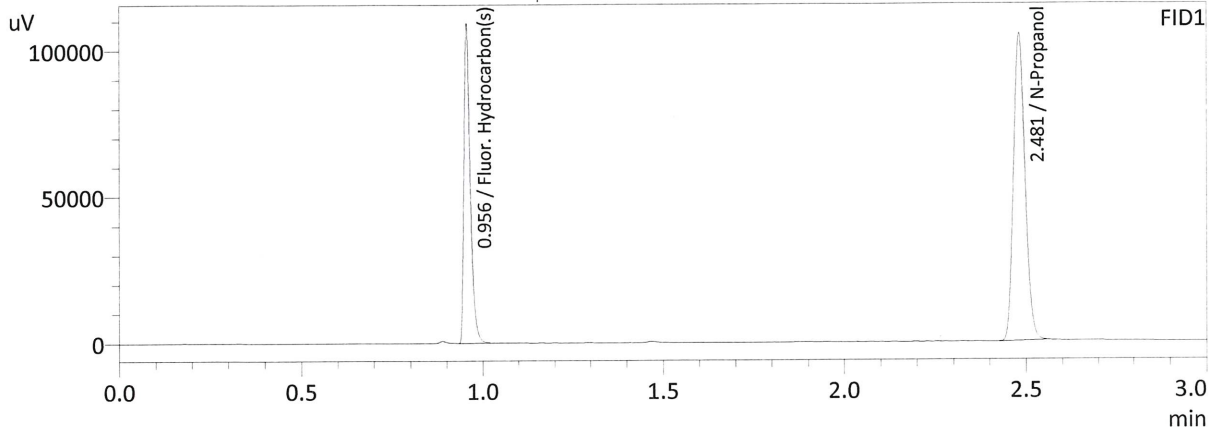

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	234080	g/100cc
Flour. Hydrocarbon(s)	0.0000	132525	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	209867	g/100cc
Flour. Hydrocarbon(s)	0.0000	122953	g/100cc

*W*

Sample Name : INT STD BLNK 4  
 Laboratory : Meridian  
 Injection Date : 7/14/2021 9:09:03 PM  
 Vial # : 53  
 Method Filename : C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM  
 Instrument #GC/HS : C12255750548 / C12595800409

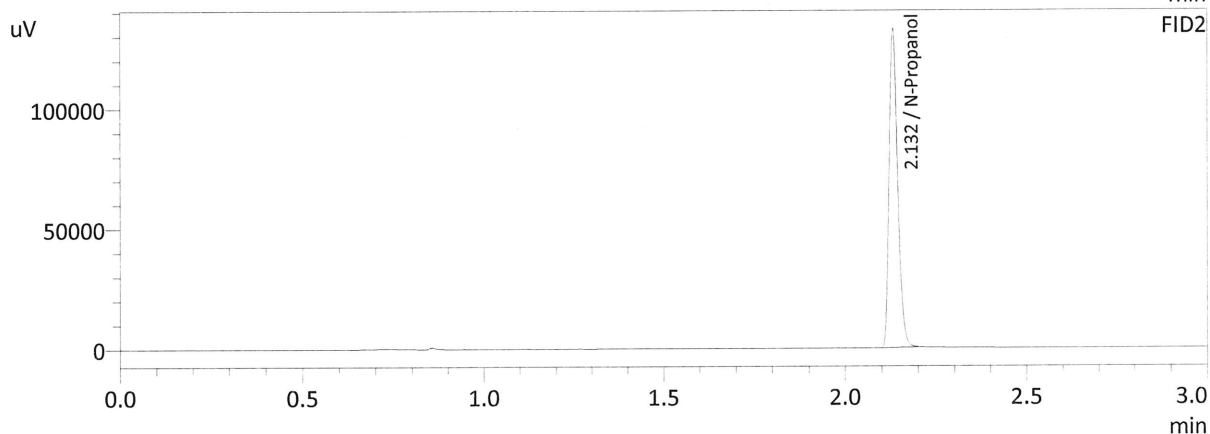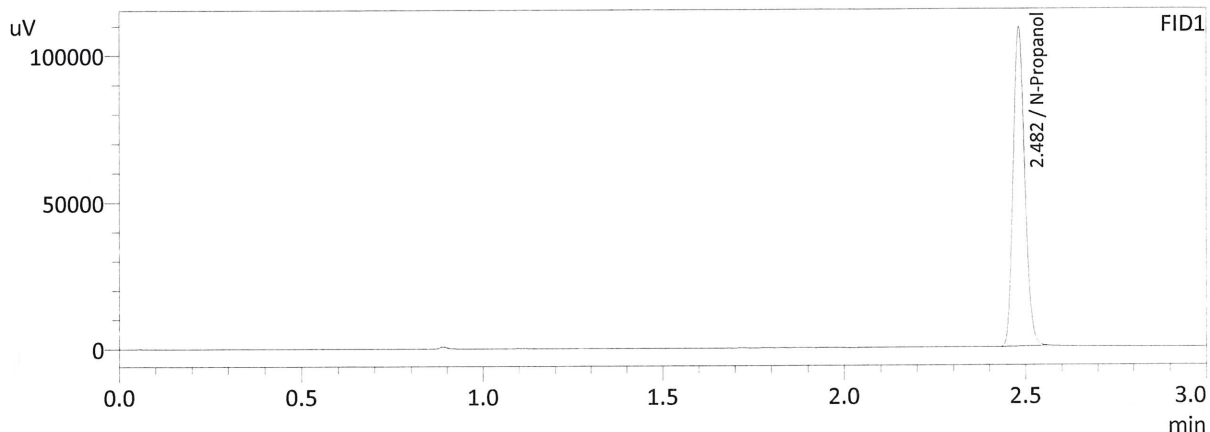

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	242329	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	217431	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

*W*

**Bradley, Nikka**

---

**From:** Bradley, Nikka  
**Sent:** Thursday, July 15, 2021 12:42 PM  
**To:** Giso, Galina  
**Subject:** Calibration Template for Worklist 5108

Will you please create a calibration template for  
C:\LabSolutions\Data\210714\TEMPLATE\CALIBRATION\ALCOHOL.GCM.

At the moment Worklist 5108 looks like calibration table included does not match the method path of the samples.

Nikka Bradley  
Forensic Scientist II  
Idaho State Police Forensic Services  
700 S. Stratford Dr. Ste 125  
Meridian, Idaho 83642  
208.884.7171  
[nikka.bradley@isp.idaho.gov](mailto:nikka.bradley@isp.idaho.gov)  
[www.isp.idaho.gov](http://www.isp.idaho.gov)

CONFIDENTIALITY NOTICE: This e-mail is intended only for the personal and confidential use of the individual(s) named as recipients (or the employee or agent responsible to deliver it to the intended recipient) and is covered by the Electronic Communications Privacy Act, 18 U.S.C. §§ 2510-2521. It may contain information that is privileged, confidential and/or protected from disclosure under applicable law including, but not limited to, the attorney client privilege and/or work product doctrine. If you are not the intended recipient of this transmission, please notify the sender immediately by telephone. Do not deliver, distribute or copy this transmission, disclose its contents or take any action in reliance on the information it contains.