







99

**Worklist: 6659**

<u>LAB CASE</u>	<u>ITEM</u>	<u>ITEM TYPE</u>	<u>DESCRIPTION</u>	
C2024-0008	1	BCK	Alcohol Analysis	
C2024-0011	2	BCK	Alcohol Analysis	
C2024-0014	1	BCK	Alcohol Analysis	
C2024-0033	2	BCK	Alcohol Analysis	
C2024-0035	1	BCK	Alcohol Analysis	
C2024-0070	1	BCK	Alcohol Analysis	
C2024-0072	1	BCK	Alcohol Analysis	
C2024-0081	1	BCK	Alcohol Analysis	
C2024-0087	1	UCK	Alcohol Analysis	
C2024-0088	1	BCK	Alcohol Analysis	
C2024-0111	1	BCK	Alcohol Analysis	
C2024-0128	1	BCK	Alcohol Analysis	
C2024-0131	1	BCK	Alcohol Analysis	
C2024-0132	1	BCK	Alcohol Analysis	

99

# Region 1 CDA Blood Alcohol Analysis Batch Table

Shimadzu GC-2030 Serial #C12255850700  
 Shimadzu HS-20 Serial #C12595700181  
 Lab Solutions DB Software Ver. 6.111  
 Copyright (C) 2008-2020 Shimadzu Corporation

Vial#	Sample Name	Sample Type	Level#	Method File
78	INT STD BLK 5	0:Unknown	0	ALCOHOL Long.gcm
79	INT STD BLK 6	0:Unknown	0	ALCOHOL Long.gcm
80	INT STD BLK 7	0:Unknown	0	ALCOHOL Long.gcm
81	INT STD BLK 8	0:Unknown	0	ALCOHOL Long.gcm
82	INT STD BLK 9	0:Unknown	0	ALCOHOL Long.gcm
83	INT STD BLK 10	0:Unknown	0	ALCOHOL Long.gcm
1	INT STD BLK 1	0:Unknown	0	ALCOHOL Long.gcm
2	0.050	1:Standard:(R)	1	ALCOHOL Long.gcm
3	0.100	1:Standard:(R)	2	ALCOHOL Long.gcm
4	0.200	1:Standard:(R)	3	ALCOHOL Long.gcm
5	0.400	1:Standard:(R)	4	ALCOHOL Long.gcm
6	0.500	1:Standard:(R)	5	ALCOHOL Long.gcm
7	INT STD BLK 2	0:Unknown	0	ALCOHOL Long.gcm
8	MULTI-COMP MIX	1:Standard:(R)	6	ALCOHOL Long.gcm
9	INT STD BLK 3	0:Unknown	0	ALCOHOL Long.gcm
10	QC-1-1	0:Unknown	0	ALCOHOL Long.gcm
11	QC-1-1-B	0:Unknown	0	ALCOHOL Long.gcm
12	0.08 QA	0:Unknown	0	ALCOHOL Long.gcm
13	0.08 QA - B	0:Unknown	0	ALCOHOL Long.gcm
14	C2024-0008-1	0:Unknown	0	ALCOHOL Long.gcm
15	C2024-0008-1-B	0:Unknown	0	ALCOHOL Long.gcm
16	C2024-0011-2	0:Unknown	0	ALCOHOL Long.gcm
17	C2024-0011-2-B	0:Unknown	0	ALCOHOL Long.gcm
18	C2024-0014-1	0:Unknown	0	ALCOHOL Long.gcm
19	C2024-0014-1-B	0:Unknown	0	ALCOHOL Long.gcm
20	C2024-0033-2	0:Unknown	0	ALCOHOL Long.gcm
21	C2024-0033-2-B	0:Unknown	0	ALCOHOL Long.gcm
22	C2024-0035-1	0:Unknown	0	ALCOHOL Long.gcm
23	C2024-0035-1-B	0:Unknown	0	ALCOHOL Long.gcm
24	C2024-0070-1	0:Unknown	0	ALCOHOL Long.gcm
25	C2024-0070-1-B	0:Unknown	0	ALCOHOL Long.gcm
26	C2024-0072-1	0:Unknown	0	ALCOHOL Long.gcm
27	C2024-0072-1-B	0:Unknown	0	ALCOHOL Long.gcm
28	C2024-0081-1	0:Unknown	0	ALCOHOL Long.gcm
29	C2024-0081-1-B	0:Unknown	0	ALCOHOL Long.gcm
30	C2024-0087-1	0:Unknown	0	ALCOHOL Long.gcm
31	C2024-0087-1-B	0:Unknown	0	ALCOHOL Long.gcm
32	QC-2-1	0:Unknown	0	ALCOHOL Long.gcm
33	QC-2-1-B	0:Unknown	0	ALCOHOL Long.gcm
34	C2024-0088-1	0:Unknown	0	ALCOHOL Long.gcm
35	C2024-0088-1-B	0:Unknown	0	ALCOHOL Long.gcm
36	C2024-0111-1	0:Unknown	0	ALCOHOL Long.gcm
37	C2024-0111-1-B	0:Unknown	0	ALCOHOL Long.gcm
38	C2024-0128-1	0:Unknown	0	ALCOHOL Long.gcm
39	C2024-0128-1-B	0:Unknown	0	ALCOHOL Long.gcm
40	C2024-0131-1	0:Unknown	0	ALCOHOL Long.gcm
41	C2024-0131-1-B	0:Unknown	0	ALCOHOL Long.gcm
42	C2024-0132-1	0:Unknown	0	ALCOHOL Long.gcm
43	C2024-0132-1-B	0:Unknown	0	ALCOHOL Long.gcm
44	QC-2-2	0:Unknown	0	ALCOHOL Long.gcm
45	QC-2-2-B	0:Unknown	0	ALCOHOL Long.gcm
46	INT STD BLK 4	0:Unknown	0	ALCOHOL Long.gcm

99

### Quantitative Analysis for Ethanol & Qualitative Analysis for Other Volatiles

*Analytical Method(s): 1.0*

*Device: Hamilton MICROLAB Liquid Processor/Dilutor Serial Number: ML600HC11379*

Volatiles Quality Assurance Controls

Run Date(s):

1-23-2024

*Calibration Date: (if different)*

Worklist #

6659

Control level	Expiration	Lot #	Target Value	Acceptable Range	Overall Results	
Level 1	Feb-25	2101199	0.0808	0.0727 - 0.0889	0.0800 g/100cc	
					g/100cc	
					g/100cc	
Level 2	Mar-26	2110181	0.2030	0.1827 - 0.2233	0.1985 g/100cc	
					0.1993 g/100cc	
					g/100cc	
<b>Multi-Component mixture:</b>		<b>Exp:</b>	January 31, 2026	<b>Lot #</b>	FN01212104	OK
<b>Curve Fit:</b>			<b>Column 1</b>	0.99972	<b>Column2</b>	0.99967

#### Ethanol Calibration Reference Material

Calibrator level	Target Value	Acceptable Range	Column 1	Column 2	Precision	Mean
50	0.050	0.045 - 0.055	0.0530	0.0532	0.0002	0.0531
100	0.100	0.090 - 0.110	0.1004	0.1005	0.0001	0.1004
200	0.200	0.180 - 0.220	0.1952	0.1950	0.0002	0.1951
300	0.300	0.270 - 0.330			0	#DIV/0!
400	0.400	0.360 - 0.440	0.3986	0.3980	0.0006	0.3983
500	0.500	0.450 - 0.550	0.5026	0.5031	0.0005	0.5028

#### Aqueous Controls

Control level	Target Value	Acceptable Range	Overall Results
80	0.080	0.076 - 0.084	0.081 g/100cc

**REVIEWED**

*By Rachel Cutler at 10:49 am, Jan 25, 2024*

Revision: 5

Issue Date: 07/05/2022

Issuing Authority: Quality Manager

99

### Internal Standard Monitoring Worksheet

<b>Worklist #:</b>	<b>6659</b>	<b>Run Date(s):</b>	<b>1-23-2024</b>
--------------------	-------------	---------------------	------------------

Internal Standard Solution: Lot# A014463901	Prep Date: 11/13/2023	Exp Date: 5/13/2024
---	-----------------------	---------------------

Sample Name	Column 1 Value	Column 2 Value
0.080	226205	232812
0.080	224391	230046
QC1	225931	232441
QC1	226546	233661
QC1		
QC1		
QC1		
QC1		
QC2	249629	257807
QC2	244048	251530
QC2	264805	270889
QC2	262743	268832
QC2		
QC2		


	Average	(-)20%	(+)20%
Column 1	240537.3	192429.8	288644.7
Column 2	247252.3	197801.8	296702.7

99

# Calibration Table


Laboratory : Coeur d'Alene  
 Instrument Name : BML8F33-Instrument1  
 Instrument Serial #: C12255850700 / C12595700181

<<Data File>>  
 Method File :Default Project - ALCOHOL Long.gcm  
 Batch File :Default Project - 1-23-24.gcb  
 Date Acquired :1/23/2024 2:26:32 PM  
 Date Created :1/23/2024 2:23:56 PM  
 Date Modified :1/23/2024 2:32:33 PM

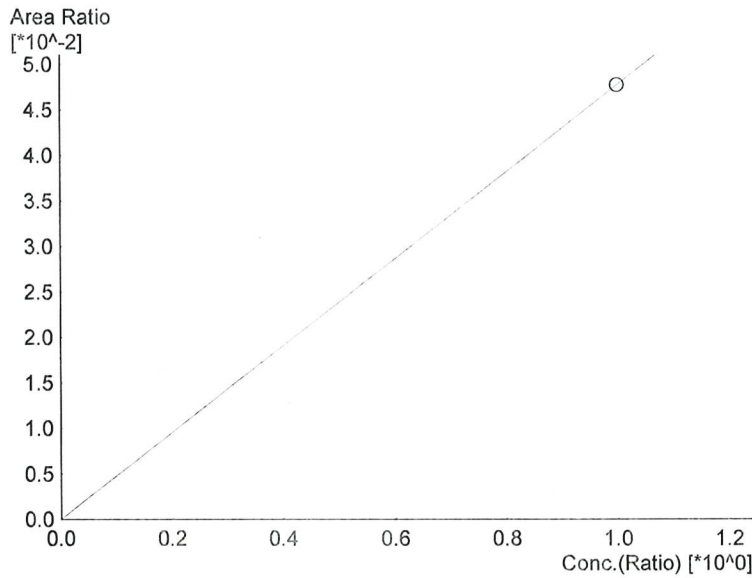

Name : Methanol  
 Detector Name: FID1  
 Function :  $f(x)=0.0477403*x+0$   
 R<sup>2</sup> value= 1.000000  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
6	1.000	10433	1.0000

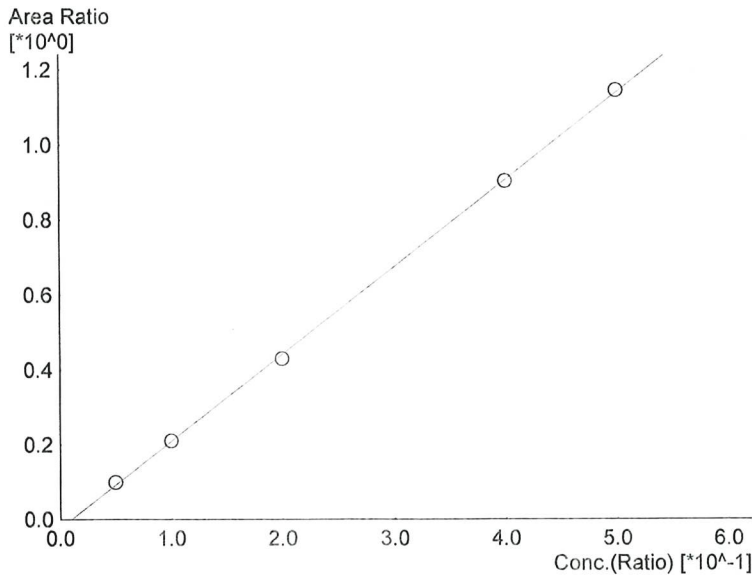

Name : Ethanol  
 Detector Name: FID1  
 Function :  $f(x)=2.32743*x-0.0242292$   
 R<sup>2</sup> value= 0.9997257 ✓  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
1	0.050	20847	0.0530
2	0.100	44581	0.1004
3	0.200	92731	0.1952
4	0.400	198796	0.3986
5	0.500	251428	0.5026

99

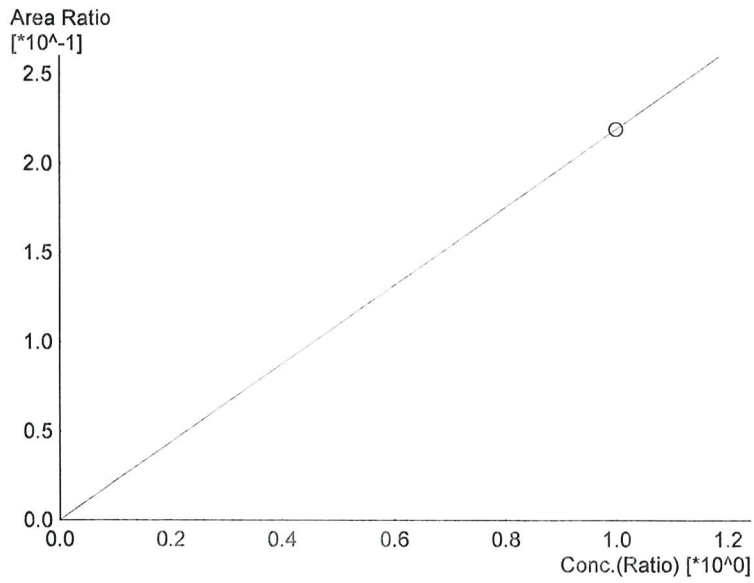

Name : Isopropyl Alcohol  
Detector Name: FID1  
Function :  $f(x)=0.219201*x+0$   
R<sup>2</sup> value= 1.000000  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
6	1.000	47905	1.0000

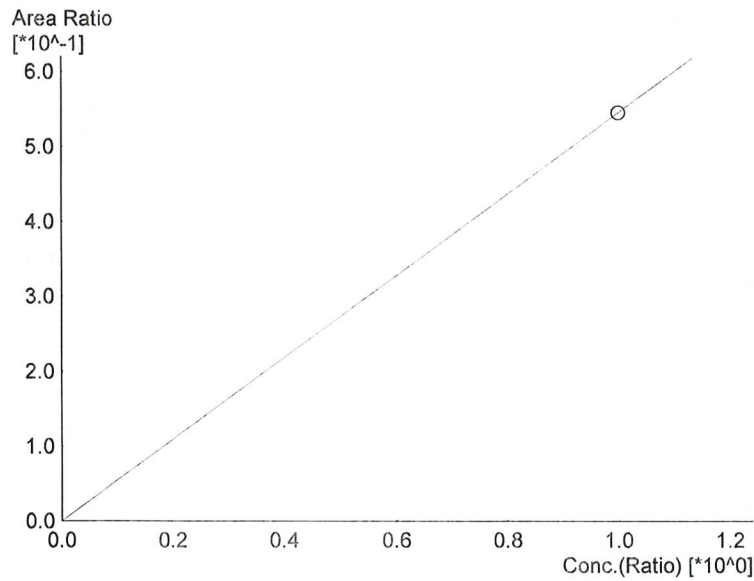

Name : Acetone  
Detector Name: FID1  
Function :  $f(x)=0.545940*x+0$   
R<sup>2</sup> value= 1.000000  
FitType: Linear  
ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
6	1.000	119313	1.0000



Name : Fluor. Hydrocarbon(s)  
Detector Name: FID1  
Function :  $f(x)=0*x+0$   
R<sup>2</sup> value= 0  
FitType: Linear  
ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
---	-------	------	------------

99

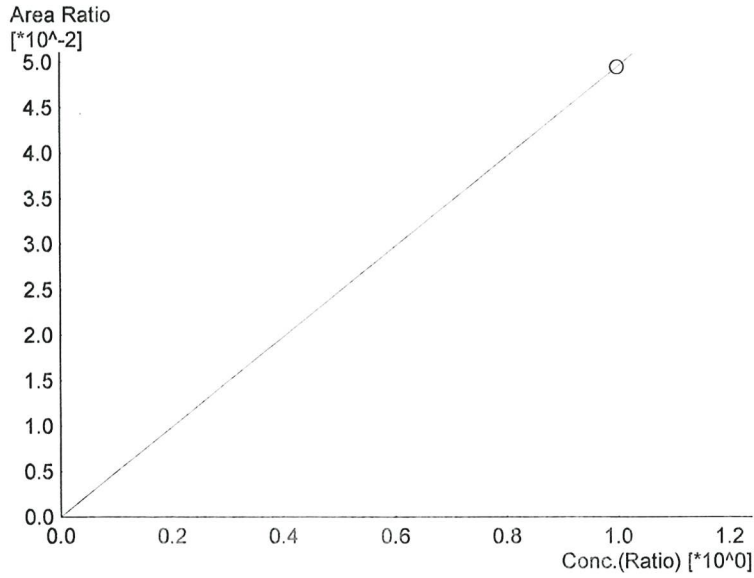


Name : Methanol  
 Detector Name: FID2  
 Function :  $f(x)=0.0494931*x+0$   
 R<sup>2</sup> value= 1.000000  
 FitType: Linear  
 ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
6	1.000	11110	1.0000

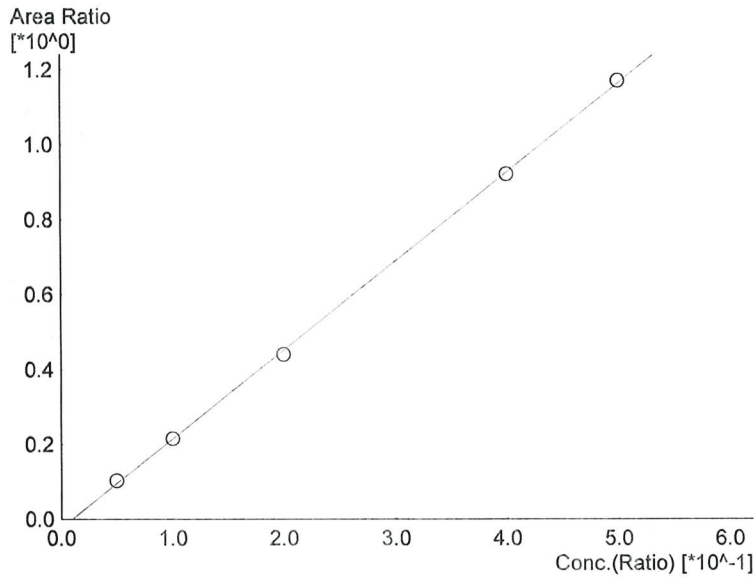

Name : Ethanol  
 Detector Name: FID2  
 Function :  $f(x)=2.37554*x-0.0236133$   
 R<sup>2</sup> value= 0.9996774 ✓  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
1	0.050	22153	0.0532
2	0.100	47079	0.1005
3	0.200	97502	0.1950
4	0.400	208922	0.3980
5	0.500	264541	0.5031

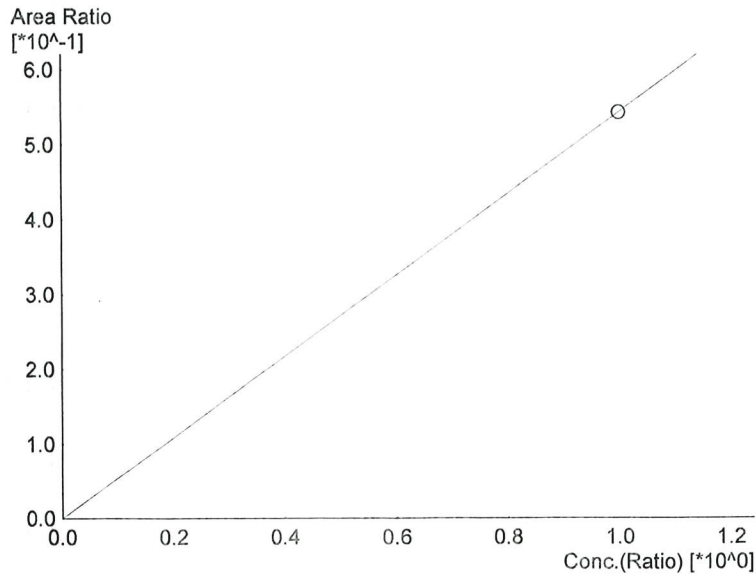

Name : Acetone  
 Detector Name: FID2  
 Function :  $f(x)=0.542993*x+0$   
 R<sup>2</sup> value= 1.000000  
 FitType: Linear  
 ZeroThrough: Not Through


#	Conc.	Area	Std. Conc.
6	1.000	121886	1.0000

99

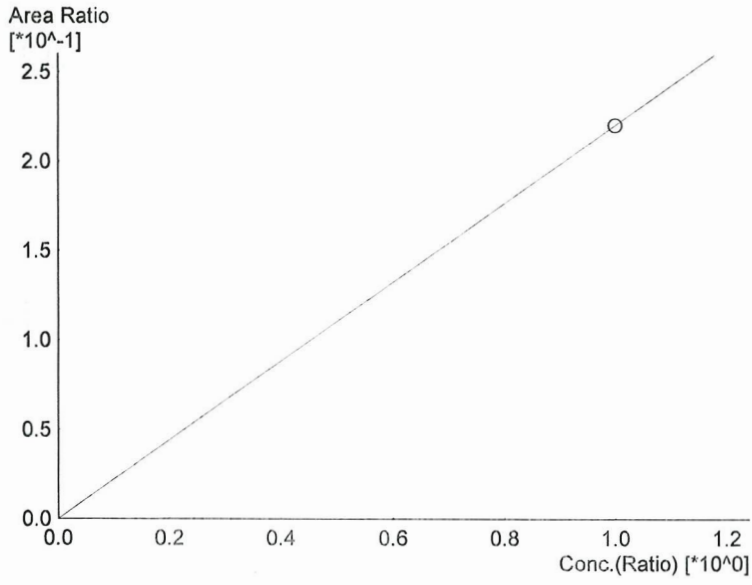

Name : Isopropyl Alcohol  
Detector Name: FID2  
Function :  $f(x)=0.220564*x+0$   
R<sup>2</sup> value= 1.000000  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
6	1.000	49510	1.0000

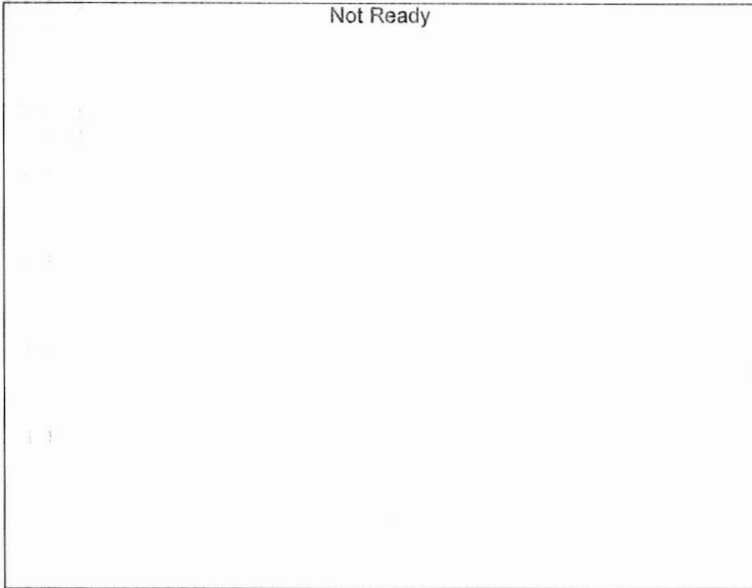



Name : Flour. Hydrocarbon(s)  
Detector Name: FID2  
Function :  $f(x)=0*x+0$   
R<sup>2</sup> value= 0  
FitType: Linear  
ZeroThrough: Not Through

#	Conc.	Area	Std. Conc.
---	-------	------	------------



99

Sample Name : 0.050  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 1:47:46 PM  
 Vial # : 2  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

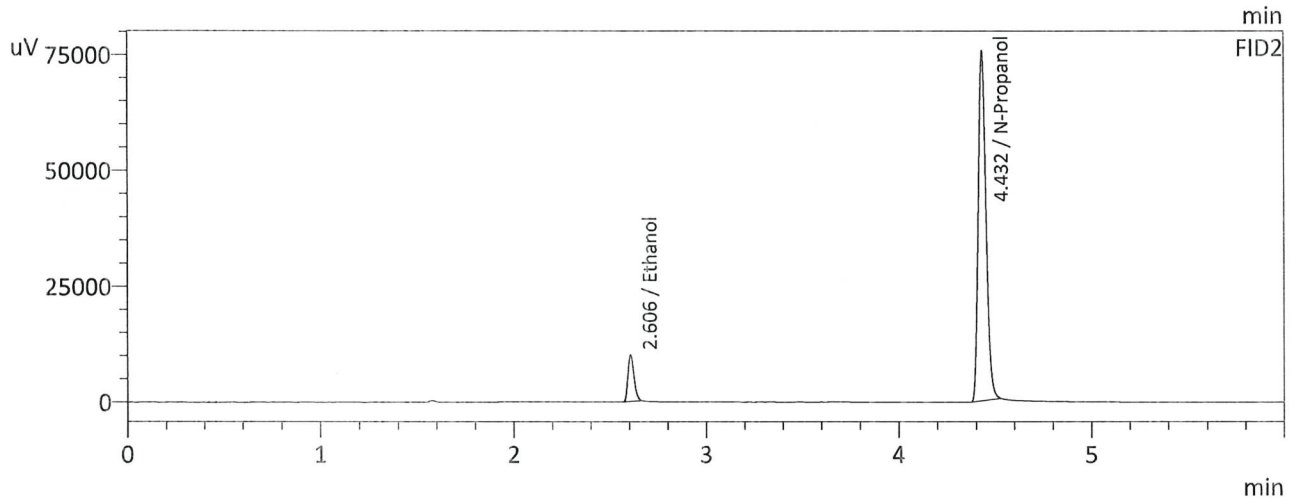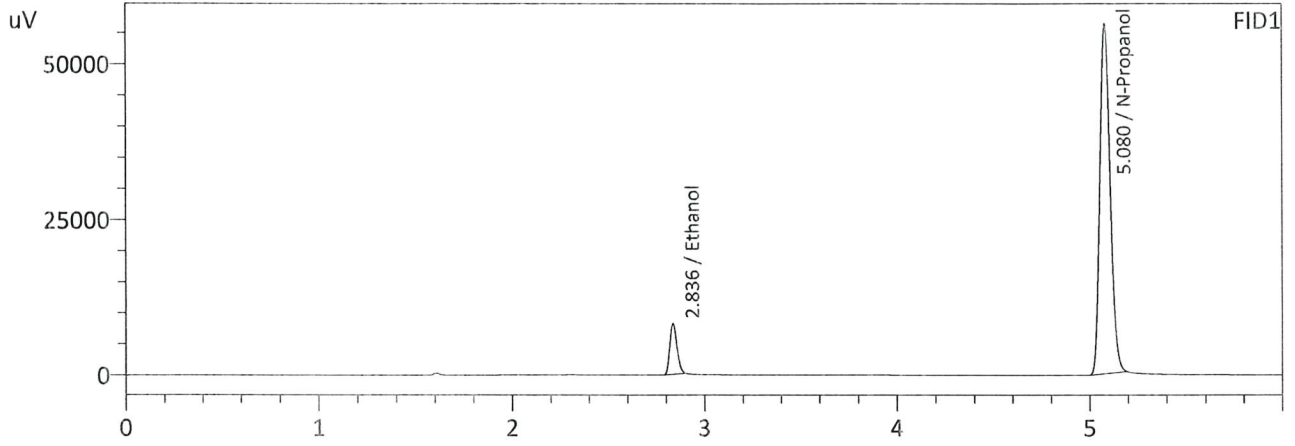

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0530	20847	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	209961	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0532	22153	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	215434	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : 0.100  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 1:58:26 PM  
 Vial # : 3  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

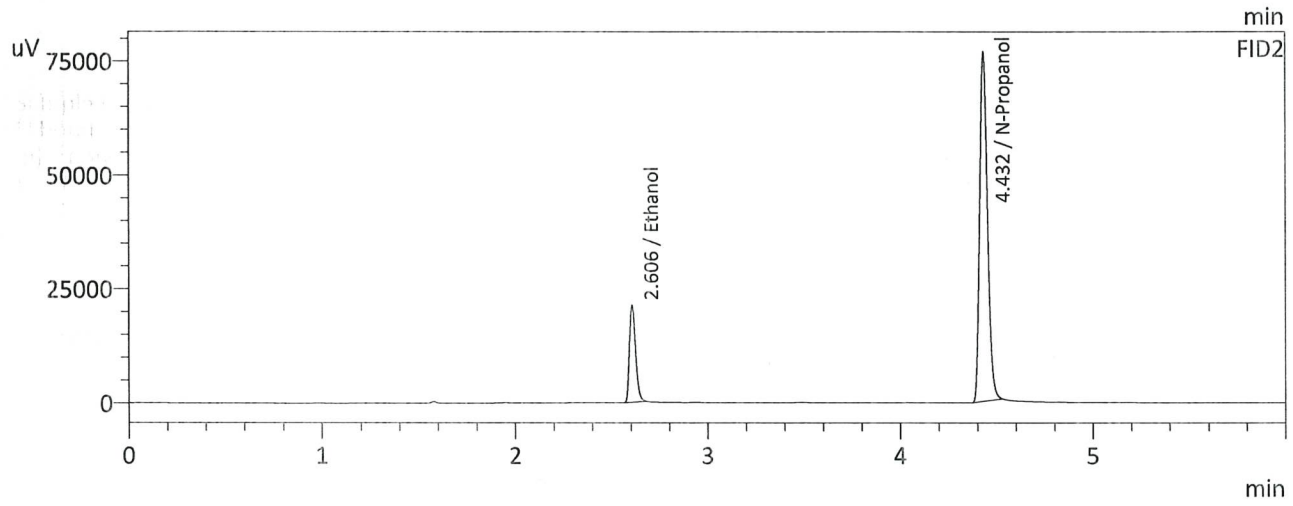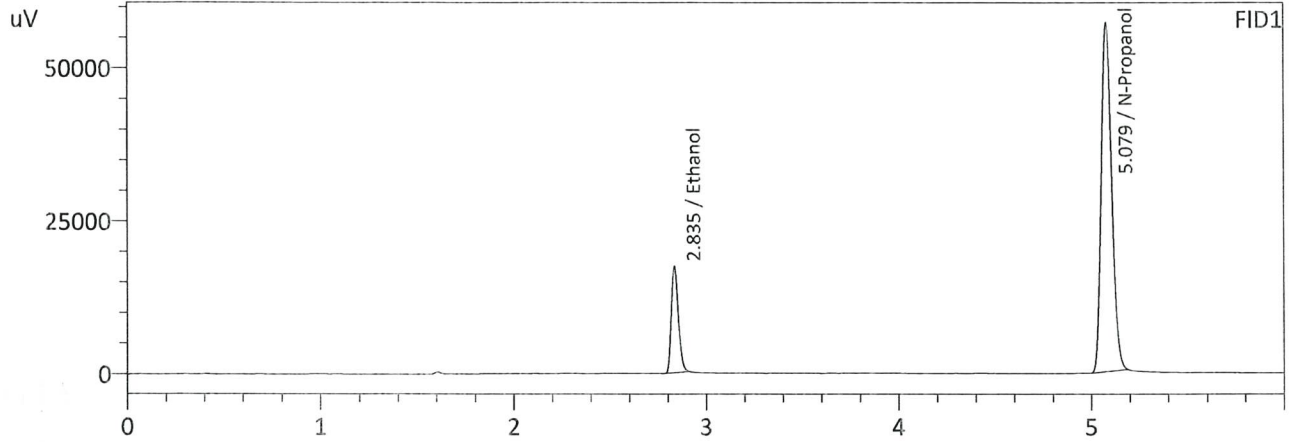

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1004	44581	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	212695	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1005	47079	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	218723	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : 0.200  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:07:07 PM  
 Vial # : 4  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

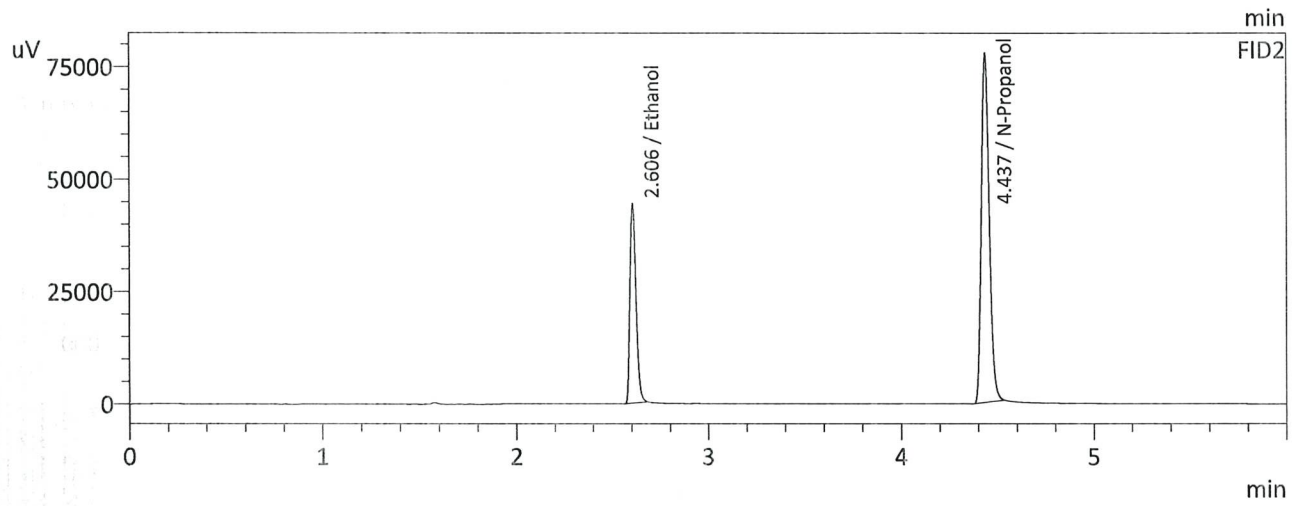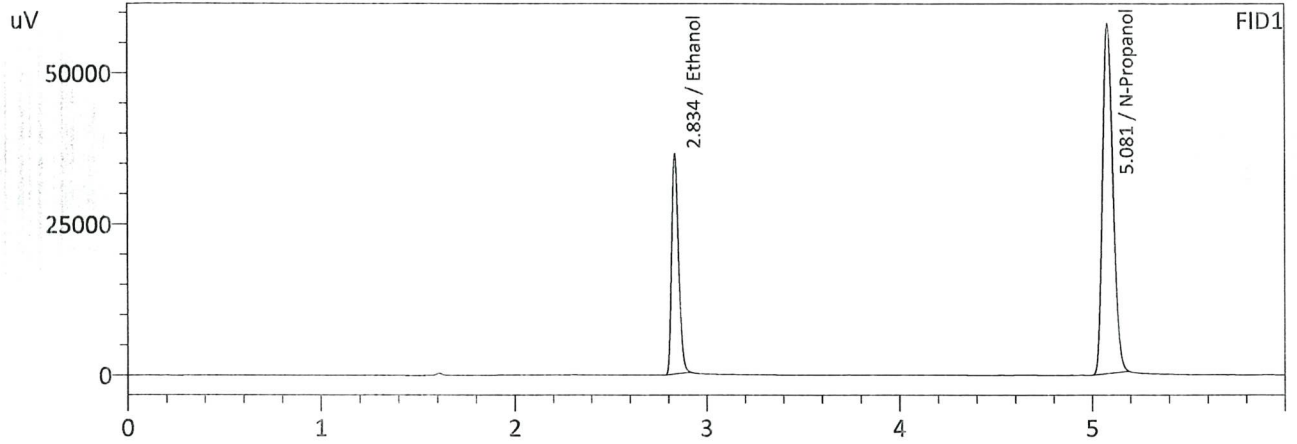

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1952	92731	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	215577	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1950	97502	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	221691	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : 0.400  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:17:52 PM  
 Vial # : 5  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

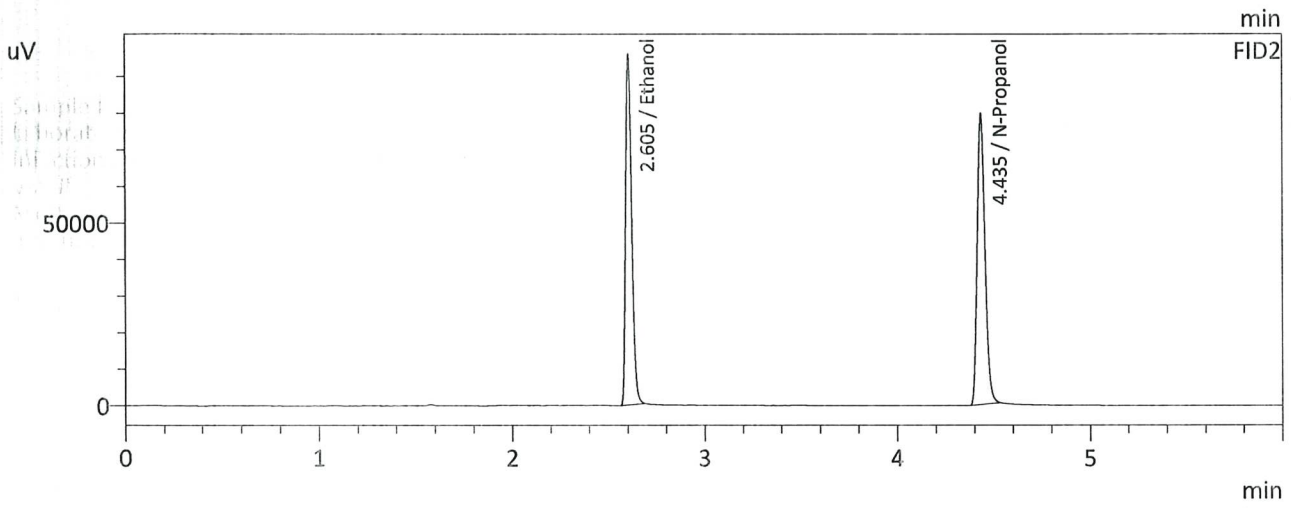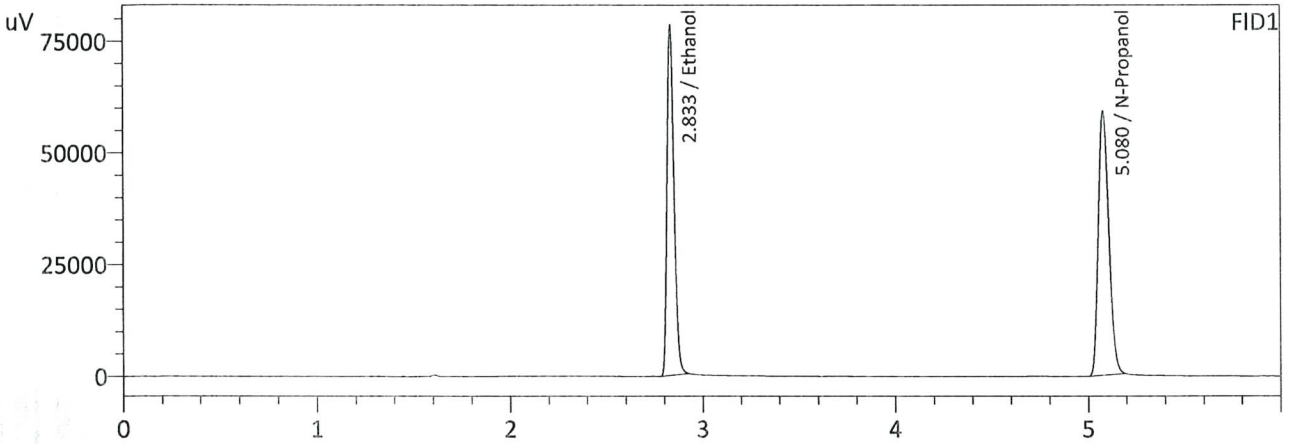

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.3986	198796	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	220022	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.3980	208922	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	226610	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : 0.500  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:26:32 PM  
 Vial # : 6  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

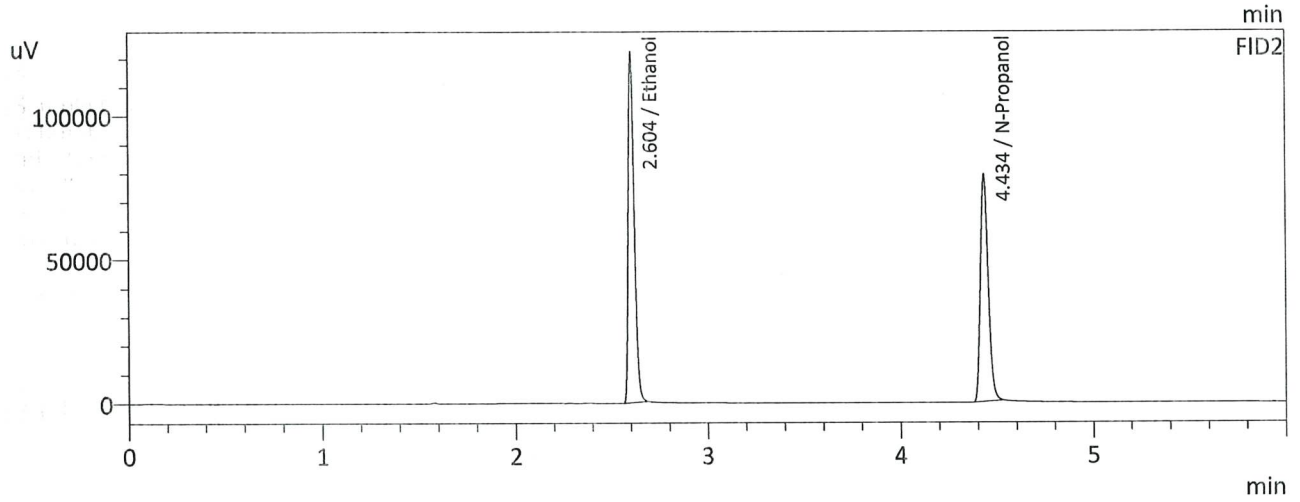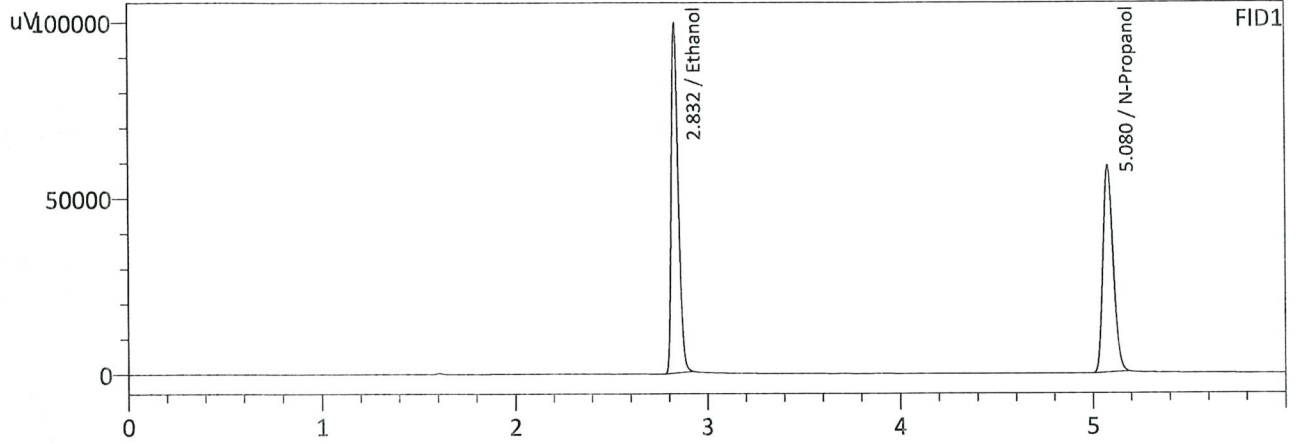

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.5026	251428	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	219479	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.5031	264541	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	225807	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : MULTI-COMP MIX  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:45:57 PM  
 Vial # : 8  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

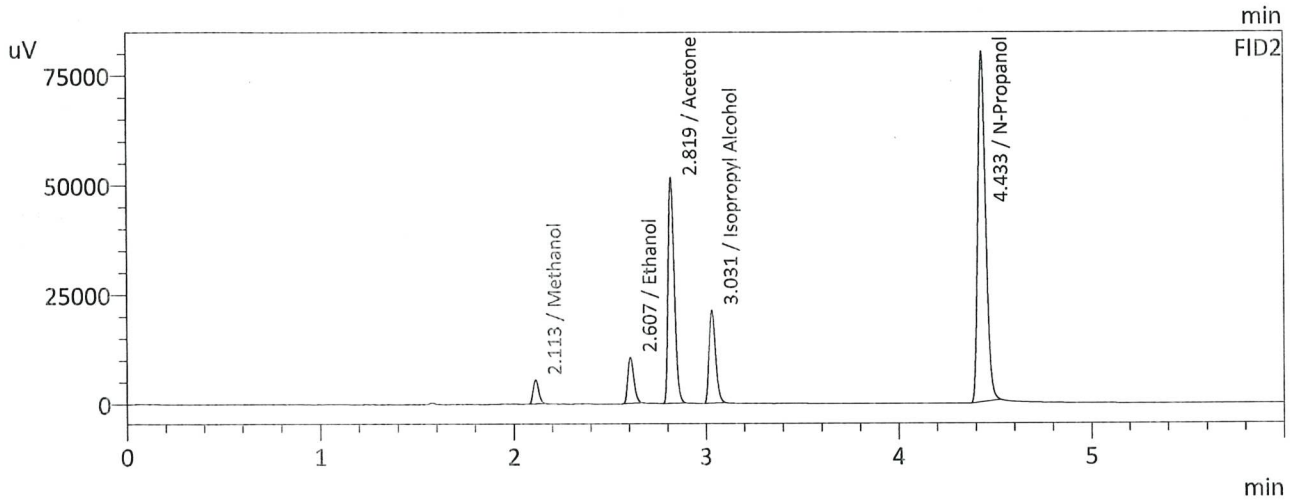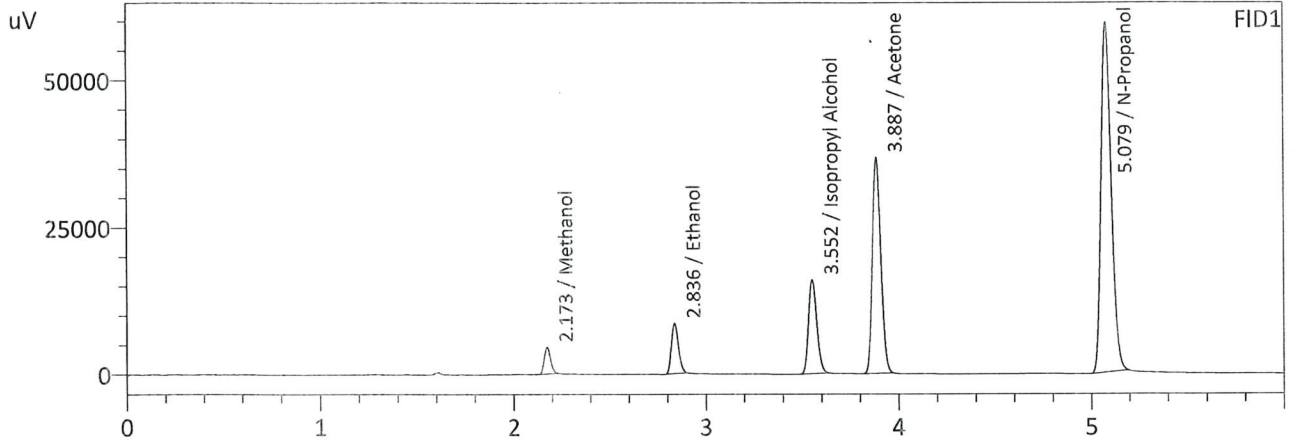

FID1



Name	Conc.	Area	Unit
Methanol	1.0000	10394	g/100cc
Ethanol	0.0527	21818	g/100cc
Isopropyl Alcohol	1.0000	47966	g/100cc
Acetone	1.0000	112421	g/100cc
N-Propanol	0.0000	221401	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	1.0000	11171	g/100cc
Ethanol	0.0531	23356	g/100cc
Acetone	1.0000	114720	g/100cc
Isopropyl Alcohol	1.0000	49626	g/100cc
N-Propanol	0.0000	227416	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : INT STD BLK 1  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 1:39:02 PM  
 Vial # : 1  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

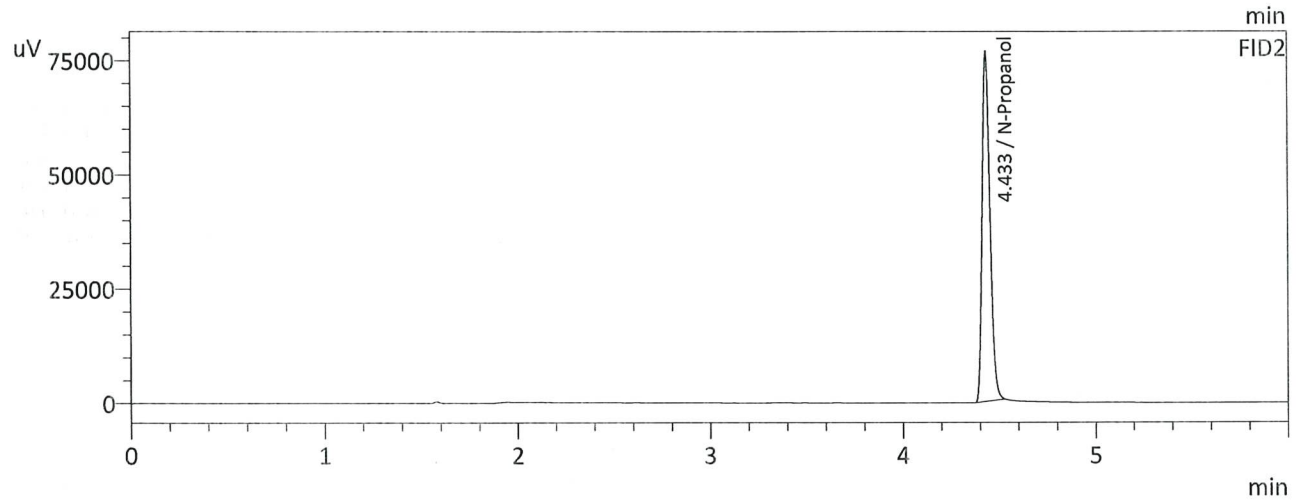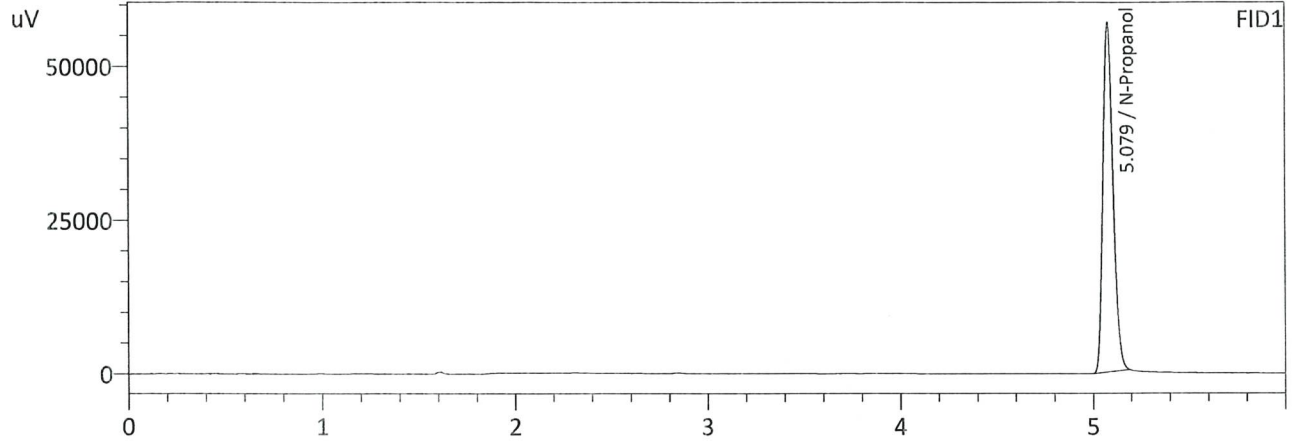

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	212251	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	217909	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : INT STD BLK 2  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:37:17 PM  
 Vial # : 7  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

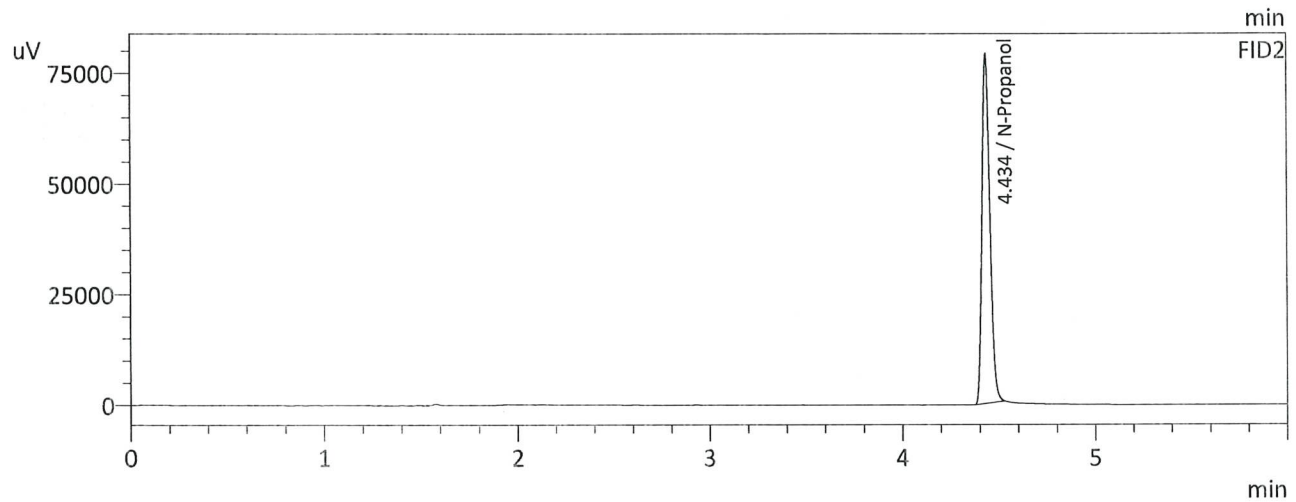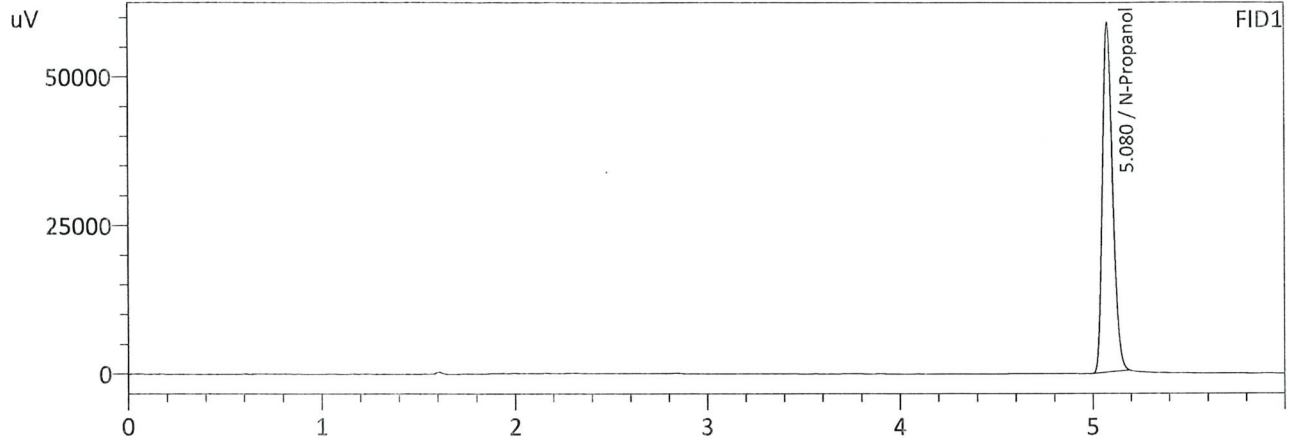

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	219142	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc



FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	225162	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc



99

Sample Name : INT STD BLK 3  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 2:56:41 PM  
 Vial # : 9  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

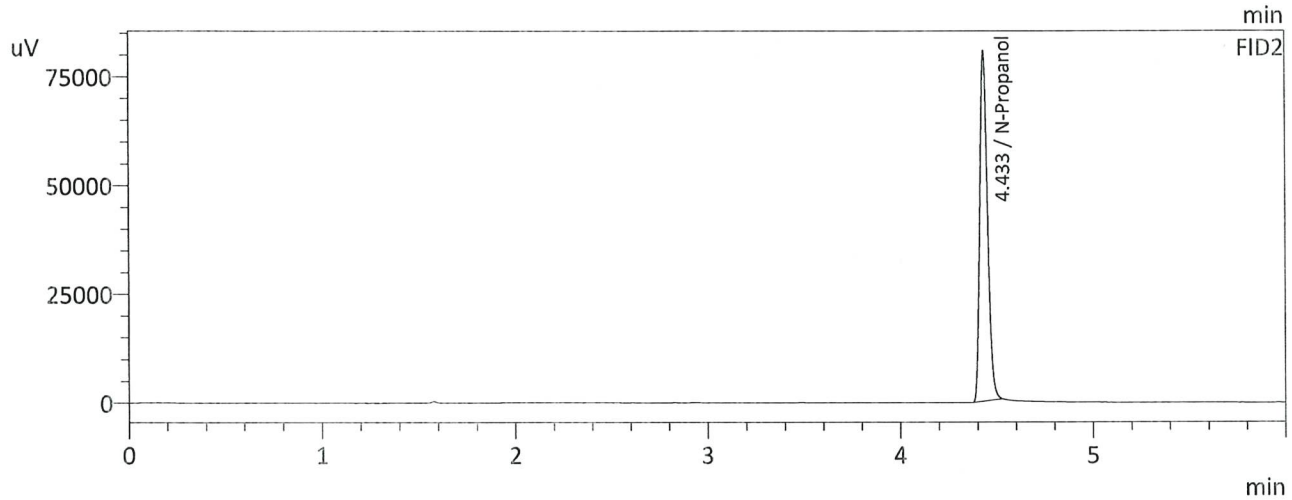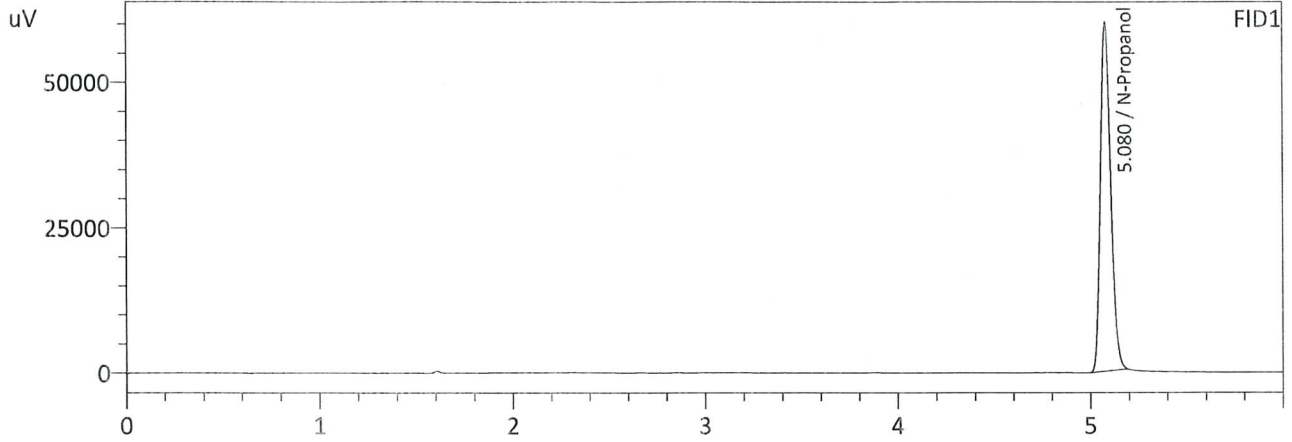

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	223159	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	229225	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : INT STD BLK 4  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 8:54:35 PM  
 Vial # : 46  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

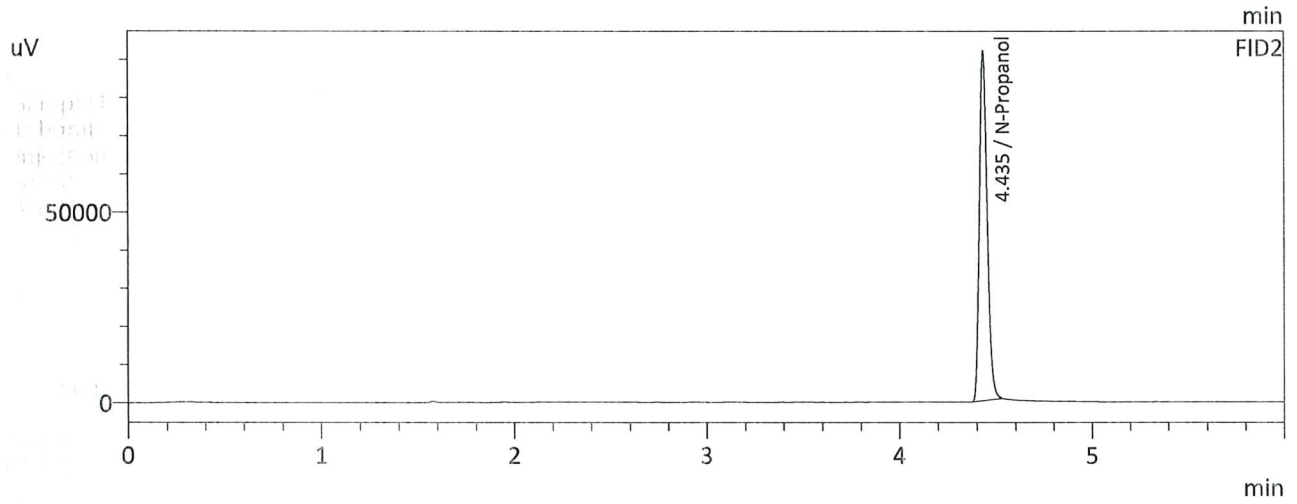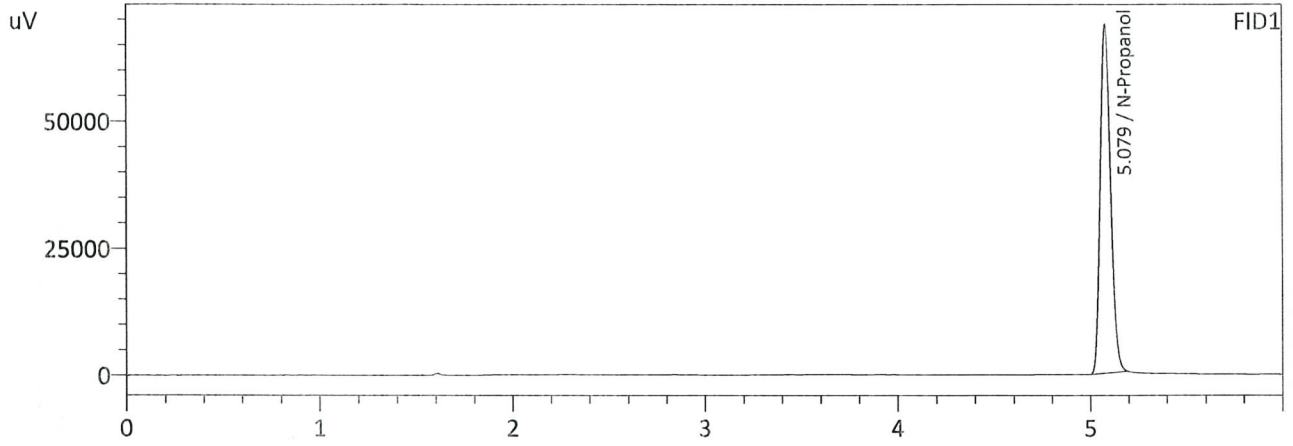

FID1:

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	254275	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2:

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	260198	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

## VOLATILES DETERMINATION CASEFILE WORKSHEET

## VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-1-1		Analysis Date(s): 1/23/2024 3:05:20 PM(-08:00)				
	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0799	0.0798	0.0001	0.0798	0.0004	0.0800
(g/100cc)	0.0802	0.0803	0.0001	0.0802		

## Analysis Method

Refer to Blood Alcohol Method #1

## Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL Long.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.080	0.076	0.084	0.004

	Reported Results	
	0.080	

Calibration and control data are stored centrally.

99

Sample Name : QC-1-1  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 3:05:20 PM  
 Vial # : 10  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

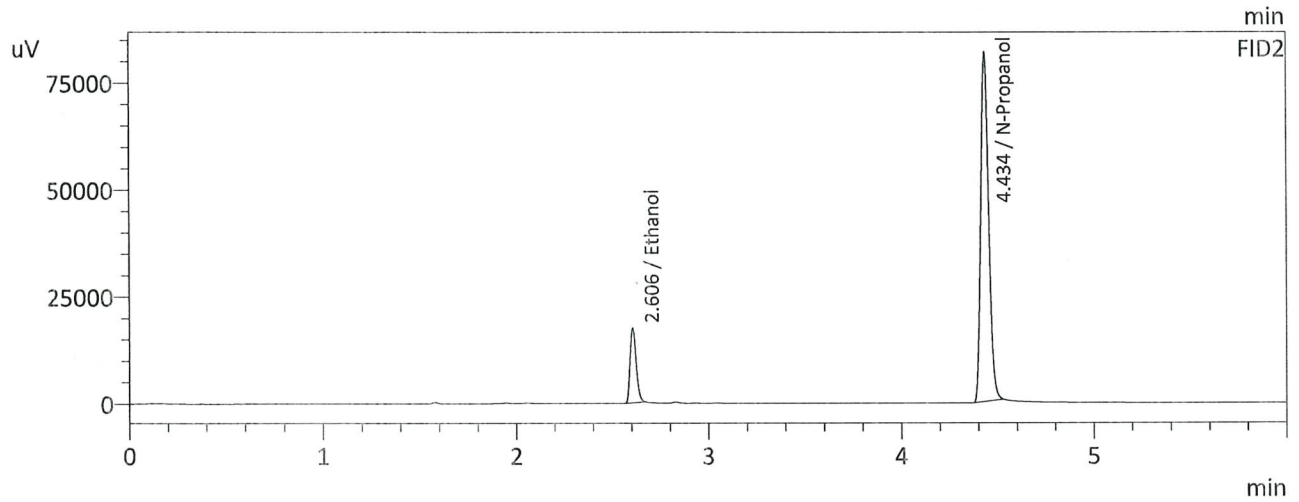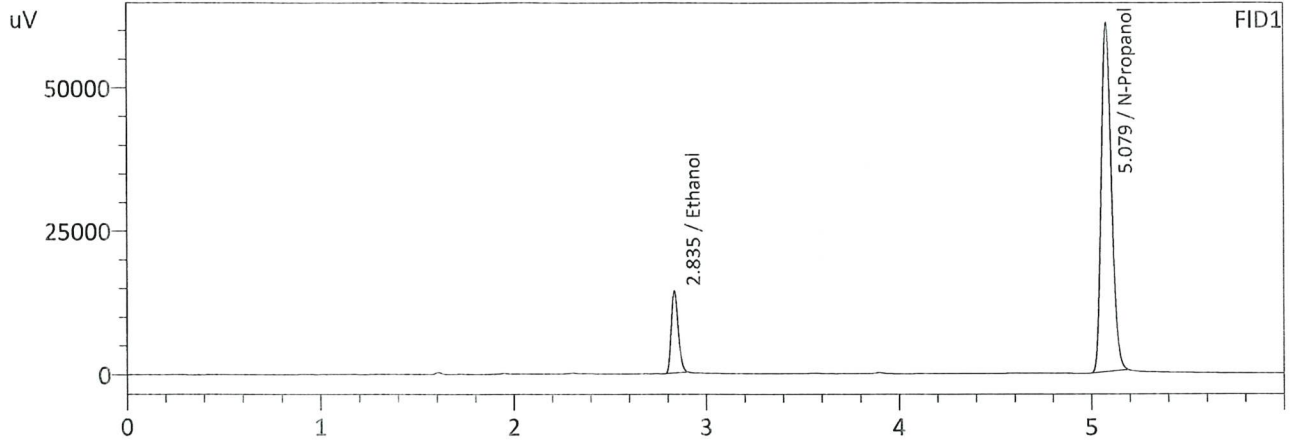

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0799	36578	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	225931	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0798	38628	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	232441	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : QC-1-1-B  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 3:16:03 PM  
 Vial # : 11  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

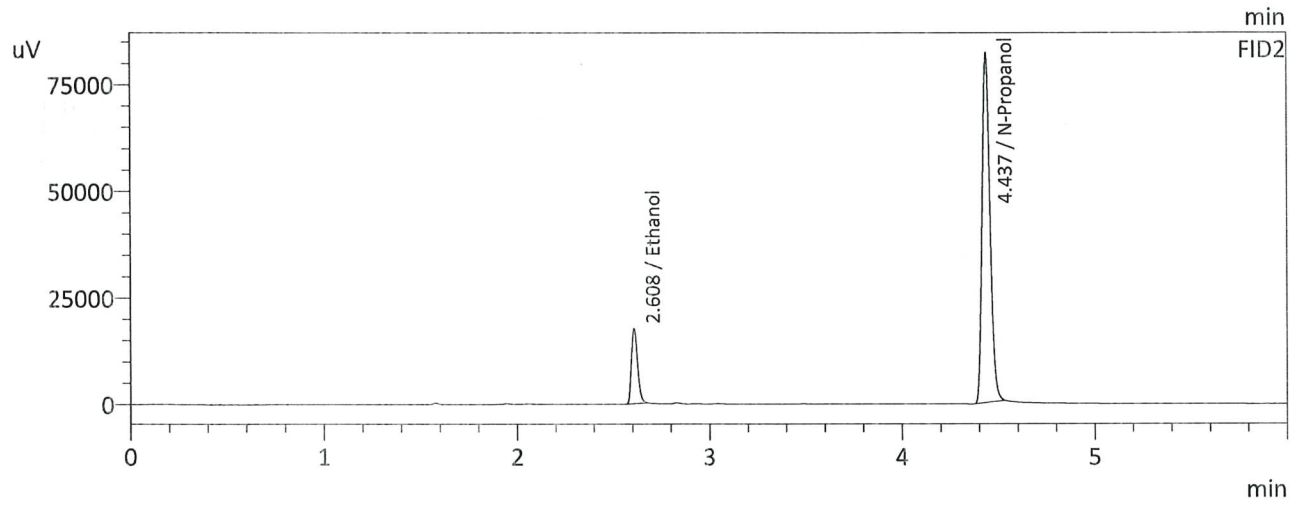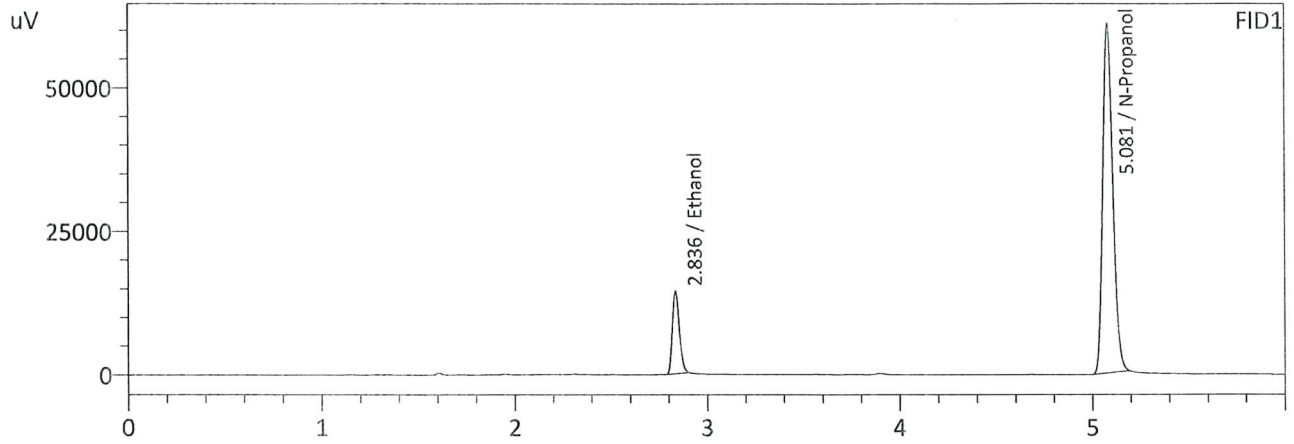

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0802	36808	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	226546	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0803	39100	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	233661	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

## VOLATILES DETERMINATION CASEFILE WORKSHEET

## VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: 0.08 QA		Analysis Date(s): 1/23/2024 3:24:44 PM(-08:00)				
	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0813	0.0818	0.0005	0.0815	0.0004	0.0817
(g/100cc)	0.0815	0.0823	0.0008	0.0819		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL Long.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.081	0.076	0.086	0.005

Reported Results	
0.081	

Calibration and control data are stored centrally.

99

Sample Name : 0.08 QA  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 3:24:44 PM  
 Vial # : 12  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

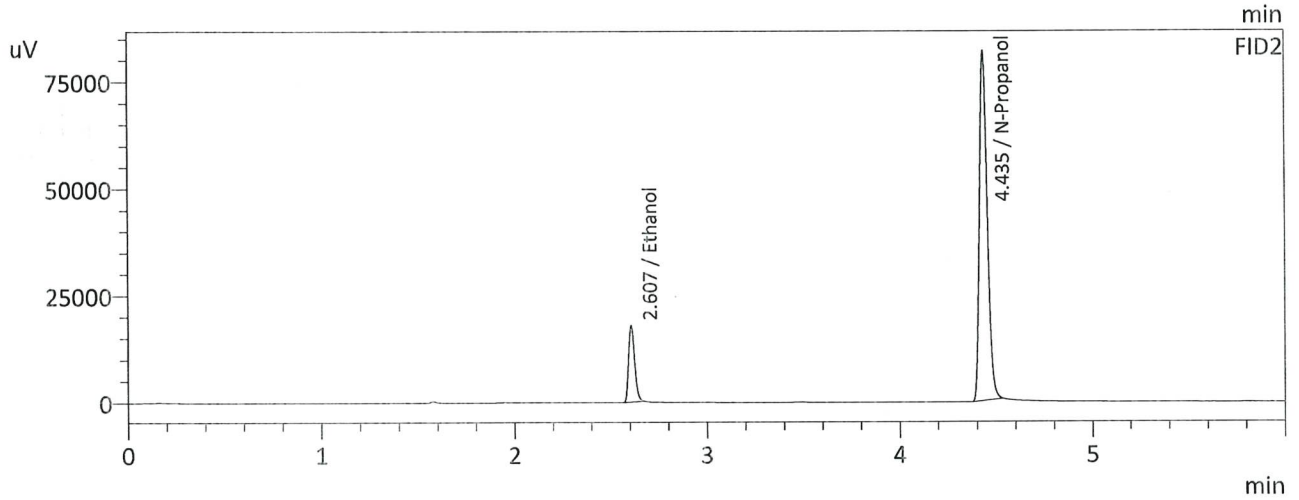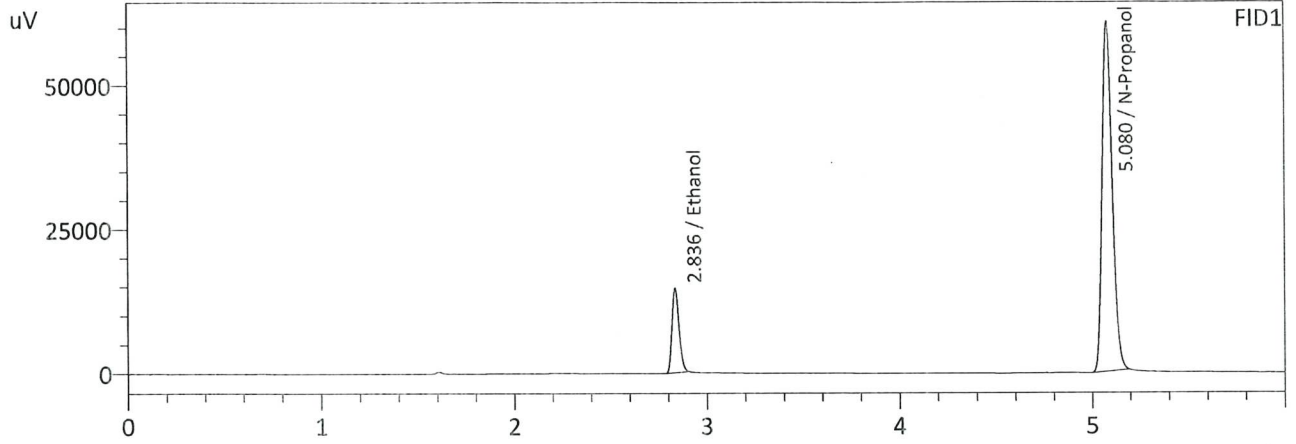

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0813	37354	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	226205	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0818	39774	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	232812	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : 0.08 QA - B  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 3:35:29 PM  
 Vial # : 13  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

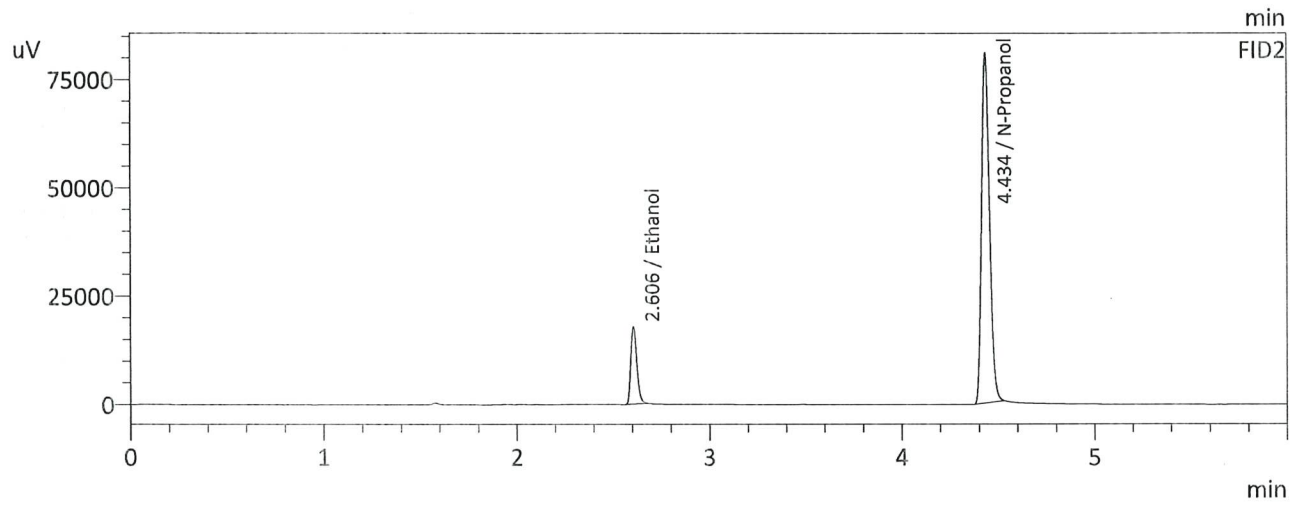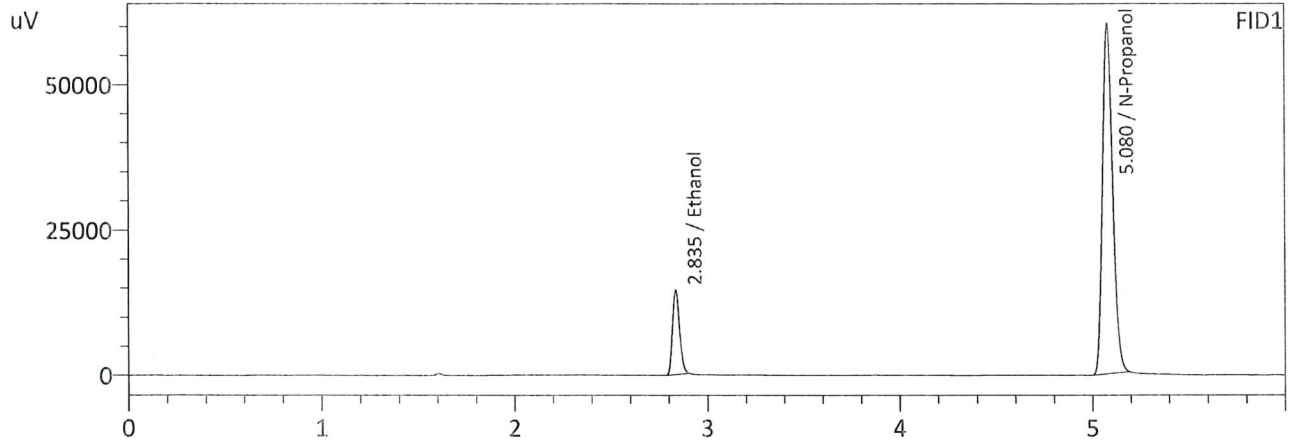

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0815	37177	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	224391	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0823	39555	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	230046	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc



99

## VOLATILES DETERMINATION CASEFILE WORKSHEET

## VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-2-1		Analysis Date(s): 1/23/2024 6:38:47 PM(-08:00)				
	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.1987	0.1985	0.0002	0.1986	0.0002	0.1985
(g/100cc)	0.1988	0.1980	0.0008	0.1984		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.



Refer To Instrument Method: ALCOHOL Long.gcm

Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.198	0.188	0.208	0.010
	Reported Results		
	0.198		

Calibration and control data are stored centrally.

99

Sample Name : QC-2-1-B  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 6:49:30 PM  
 Vial # : 33  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

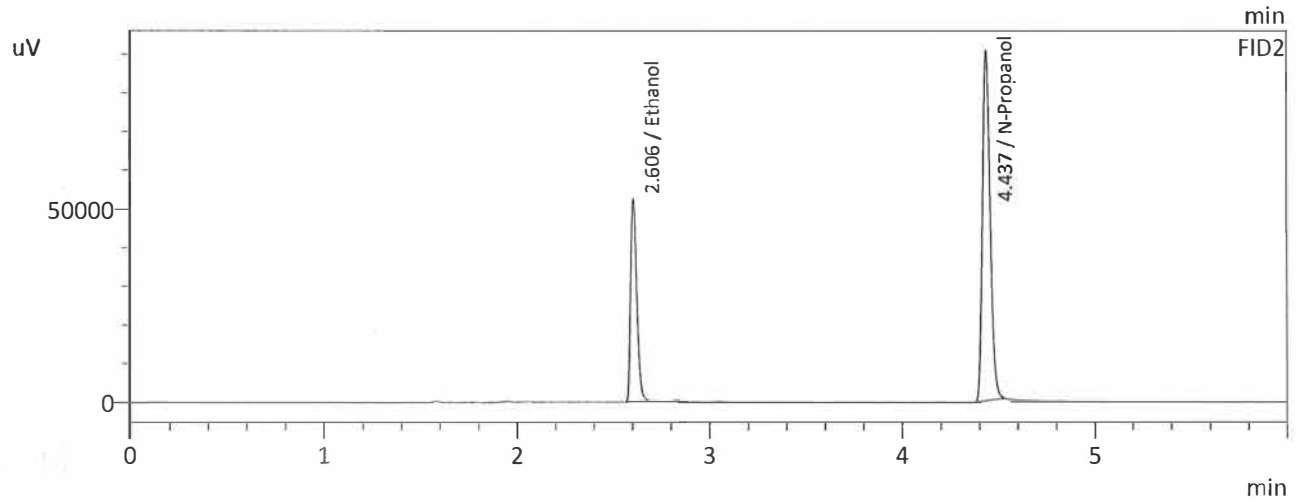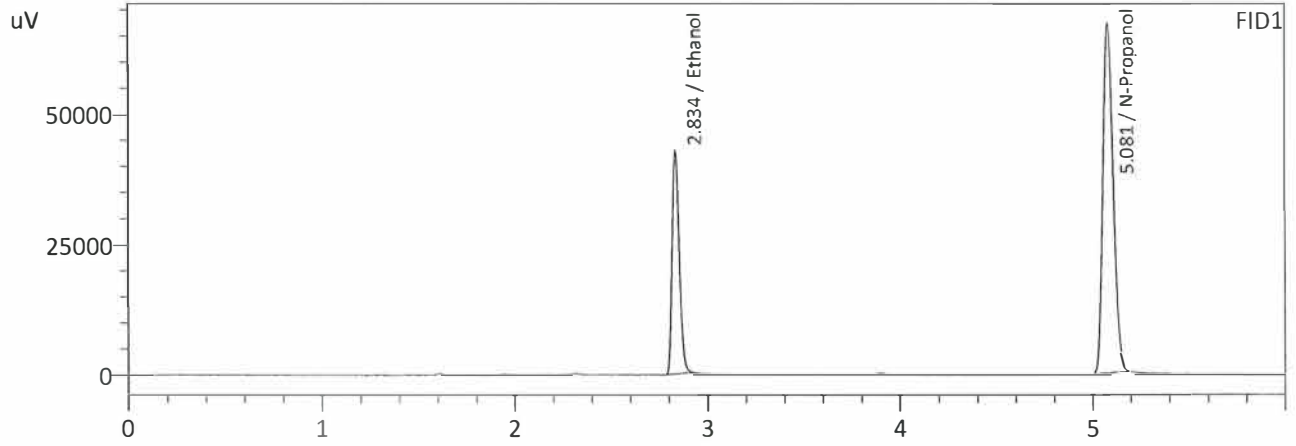

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1988	109494	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	249629	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1980	115208	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	257807	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : QC-2-1  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 6:38:47 PM  
 Vial # : 32  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

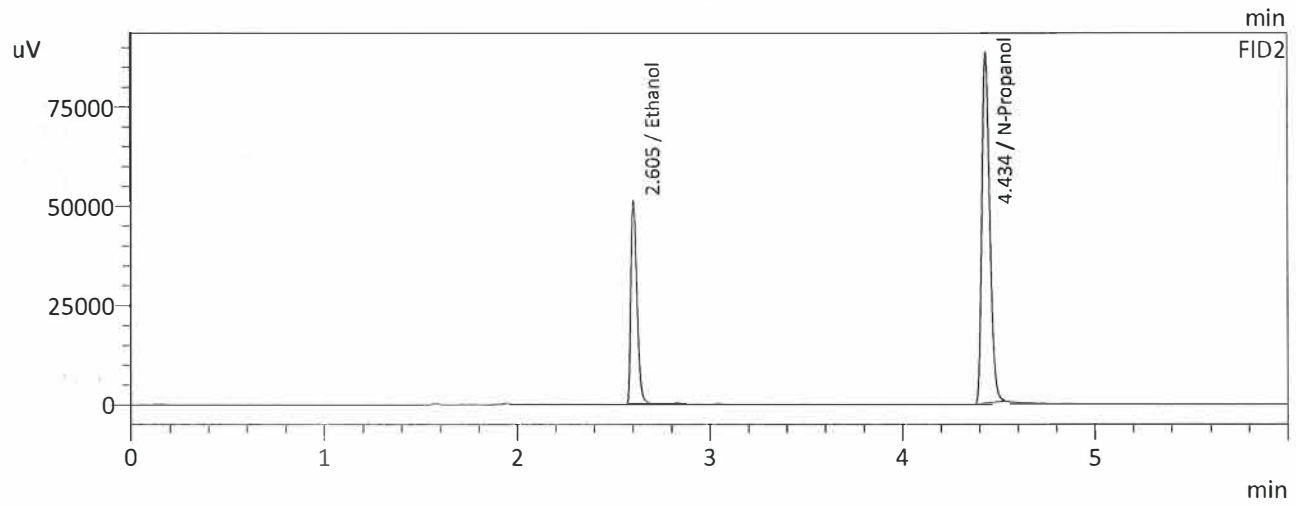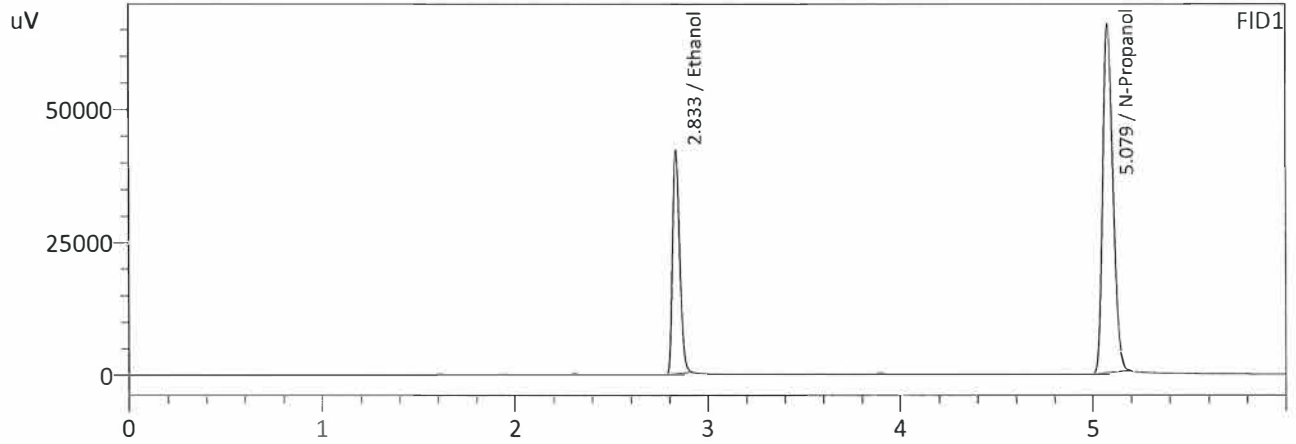

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1987	106953	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	244048	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1985	112699	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	251530	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

VOLATILES DETERMINATION CASEFILE WORKSHEET

VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-2-2		Analysis Date(s): 1/23/2024 8:35:18 PM(-08:00)				
	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.1999	0.1994	0.0005	0.1996	0.0007	0.1993
(g/100cc)	0.1994	0.1985	0.0009	0.1989		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL Long.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.199	0.189	0.209	0.010

	Reported Results
	0.199

Calibration and control data are stored centrally.

99

Sample Name : QC-2-2-B  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 8:46:03 PM  
 Vial # : 45  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

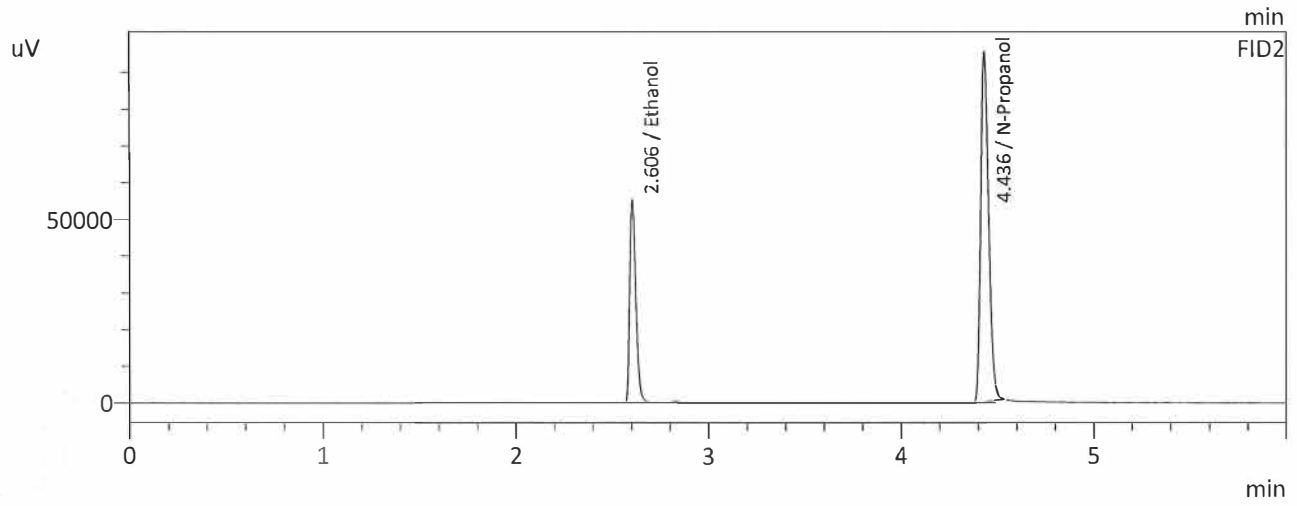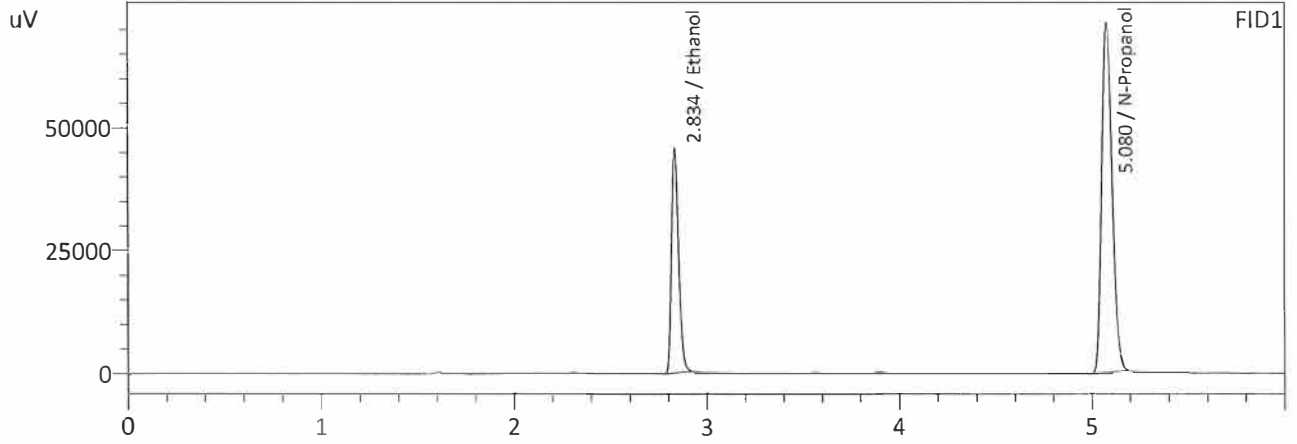

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1994	116496	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	264805	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1985	121365	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	270889	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

99

Sample Name : QC-2-2  
 Laboratory : Coeur d' Alene Lab  
 Injection Date : 1/23/2024 8:35:18 PM  
 Vial # : 44  
 Method Filename : Default Project - ALCOHOL Long.gcm  
 Instrument #GC/HS : C12255850700 / C12595700181

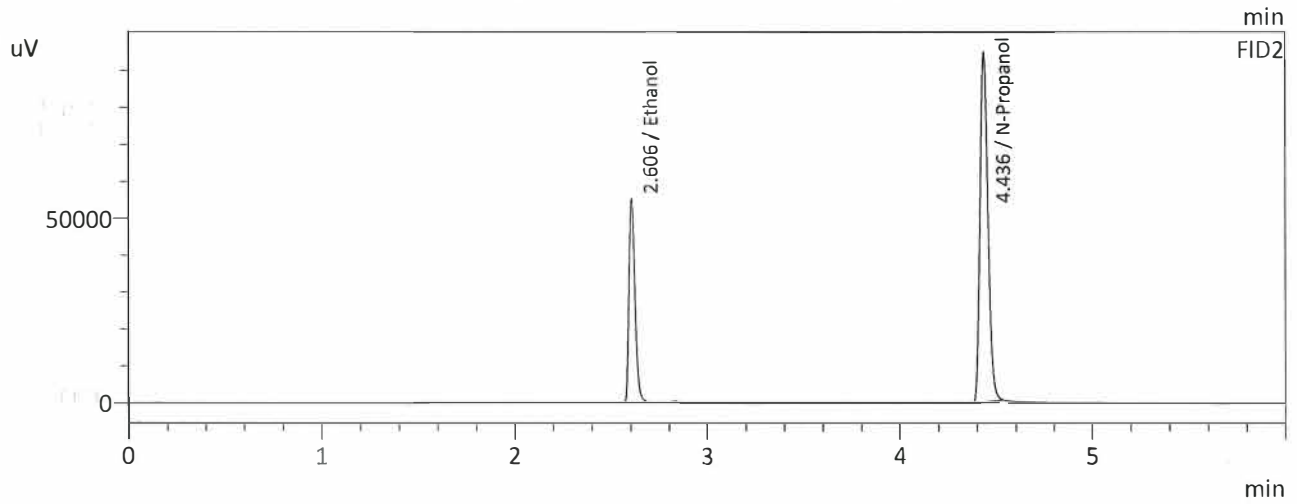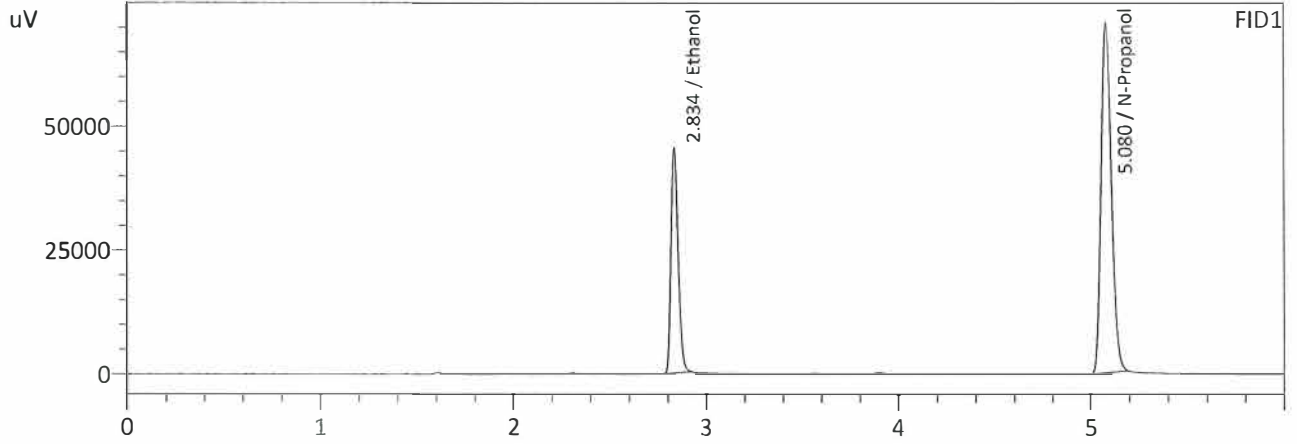

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1999	115933	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	262743	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.1994	121039	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	268832	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc