

Quantitative Analysis for Ethanol & Qualitative Analysis for Other Volatiles

Analytical Method(s): 1.0

Device: Hamilton MICROLAB Liquid Processor/Dilutor Serial Number: ML600HC11378

Volatiles Quality Assurance Controls

Run Date(s): 4/17/24

Calibration Date: (if different) 4/11/24

Worklist #: 6774

Control level	Expiration	Lot #	Target Value	Acceptable Range	Overall Results
Level 1	Feb-25	2101199	0.0808	0.0727-0.0889	0.0808 g/100cc
					0.0854 g/100cc
					g/100cc
Level 2	Mar-26	2110181	0.2030	0.1827-0.2233	0.2097 g/100cc
					g/100cc
					g/100cc
Multi-Component mixture:		Exp:	Oct. 24	Lot #	FN06041902
Curve Fit:			Column 1	0.99964	Column2 0.99959

Ethanol Calibration Reference Material

Calibrator level	Target Value	Acceptable Range	Column 1	Column 2	Precision	Mean
50	0.050	0.045 - 0.055	0.0537	0.0537	0	0.0537
100	0.100	0.090 - 0.110	0.1002	0.1005	0.0003	0.1003
200	0.200	0.180 - 0.220	0.1951	0.1946	0.0005	0.1948
300	0.300	0.270 - 0.330	0.2984	0.2983	1E-04	0.2983
400	0.400	0.360 - 0.440	N/A	N/A	#####	#DIV/0!
500	0.500	0.450 - 0.550	0.5024	0.5026	0.0002	0.5025

Aqueous Controls

Control level	Target Value	Acceptable Range	Overall Results
80	0.080	0.076 - 0.084	0.082 g/100cc

REVIEWED
 By Jeremy Johnston at 11:23 am, Apr 18, 2024

JL

Revision: 5

Issue Date: 07/05/2022

Issuing Authority: Quality Manager

Internal Standard Monitoring Worksheet

Worklist #:	6774	Run Date(s):	4/17/24
--------------------	-------------	---------------------	----------------

Internal Standard Solution:	Prep Date: 3/13/2024	Exp Date: 9/13/2024
-----------------------------	----------------------	---------------------

Sample Name	Column 1 Value	Column 2 Value
0.080	187246	200437
0.080	187800	200789
QC1	190491	204221
QC1	190068	203515
QC1	199251	214586
QC1	213180	229807
QC1		
QC1		
QC2	202387	217997
QC2	205087	220519
QC2		

	Average	(-)20%	(+)20%
Column 1	196938.8	157551.0	236326.5
Column 2	211483.9	169187.1	253780.7

JL

Revision: 5


Issue Date: 07/05/2022

Issuing Authority: Quality Manager



U6

Worklist: 6774

<u>LAB CASE</u>	<u>ITEM</u>	<u>ITEM TYPE</u>	<u>DESCRIPTION</u>
M2024-1349	2	BCK	Alcohol Analysis
M2024-1430	1	BCK	Alcohol Analysis
M2024-1431	1	BCK	Alcohol Analysis
M2024-1451	1	BCK	Alcohol Analysis
M2024-1459	1	BCK	Alcohol Analysis
M2024-1460	1	BCK	Alcohol Analysis
M2024-1489	1	BCK	Alcohol Analysis
M2024-1492	1	BCK	Alcohol Analysis
M2024-1539	1	BCK	Alcohol Analysis
M2024-1556	1	BCK	Alcohol Analysis
M2024-1563	1	BCK	Alcohol Analysis

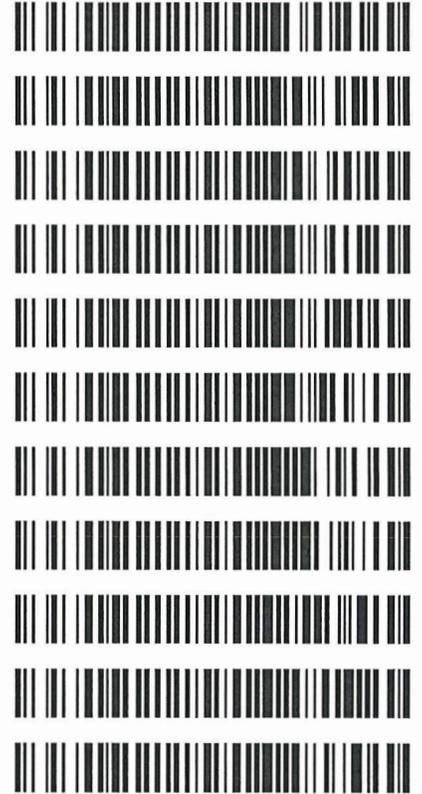

Sample Name : ISTD BLK 1
 Laboratory : Meridian
 Injection Date : 4/17/2024 2:51:09 PM
 Vial # : 1
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

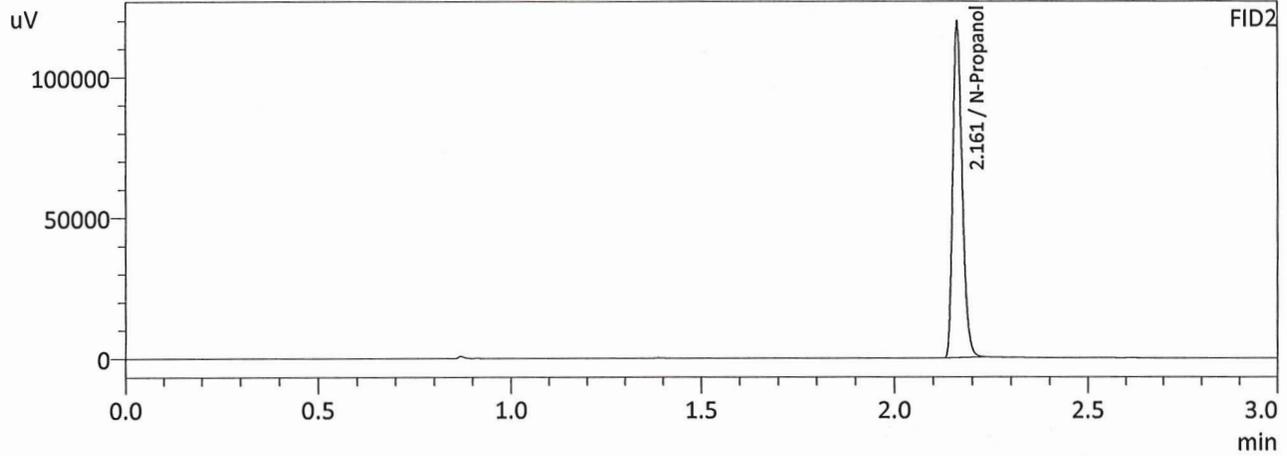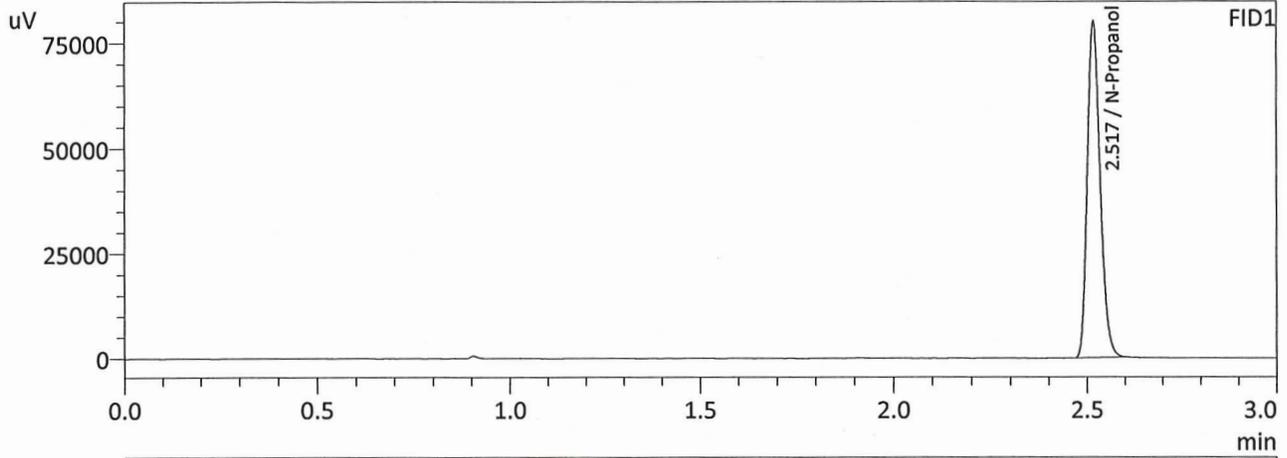

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	186696	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	199110	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

JC

Sample Name : MIXED VOLATILES FN 06041902
 Laboratory : Meridian
 Injection Date : 4/17/2024 2:58:29 PM
 Vial # : 2
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

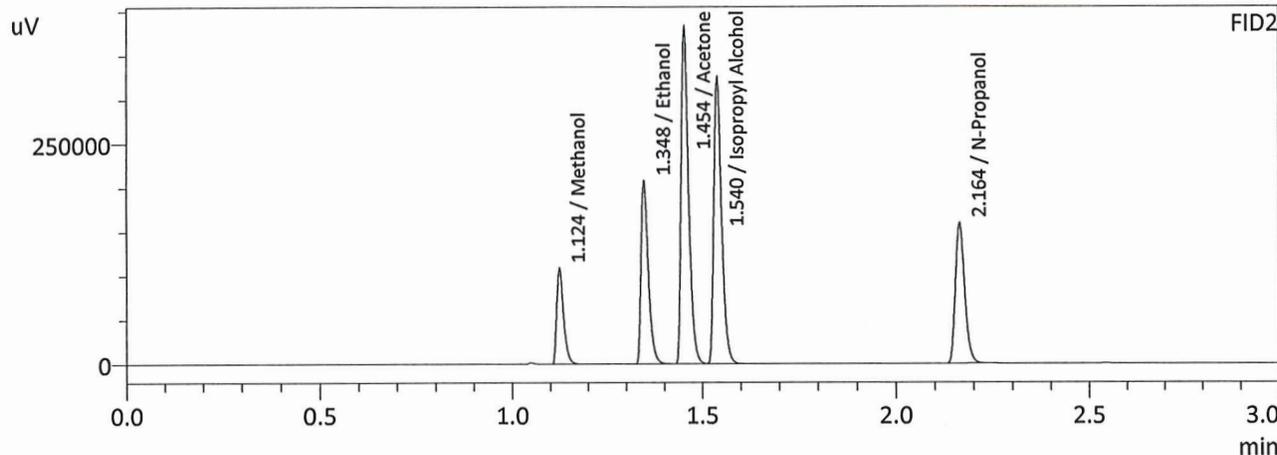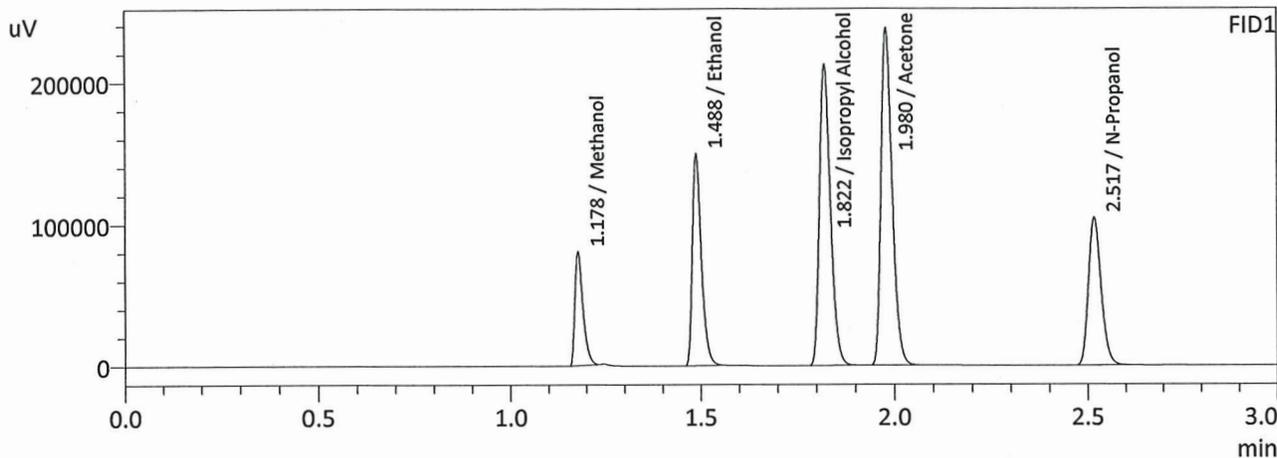

FID1

Name	Conc.	Area	Unit
Methanol	0.0000	118895	g/100cc
Ethanol	0.4710	247470	g/100cc
Isopropyl Alcohol	0.0000	412984	g/100cc
Acetone	0.0000	465694	g/100cc
N-Propanol	0.0000	242117	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	0.0000	136351	g/100cc
Ethanol	0.4747	274694	g/100cc
Acetone	0.0000	514148	g/100cc
Isopropyl Alcohol	0.0000	451976	g/100cc
N-Propanol	0.0000	264417	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

JK

VOLATILES DETERMINATION CASEFILE WORKSHEET

VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-1-1		Analysis Date(s): 4/17/2024 3:05:55 PM(-06:00)				
	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0807	0.0808	0.0001	0.0807	0.0003	0.0808
(g/100cc)	0.0810	0.0810	0.0000	0.0810		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL_240411JG.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.080	0.076	0.084	0.004

Reported Results	
0.080	

Calibration and control data are stored centrally.

JC

Sample Name : QC-1-1
 Laboratory : Meridian
 Injection Date : 4/17/2024 3:05:55 PM
 Vial # : 3
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

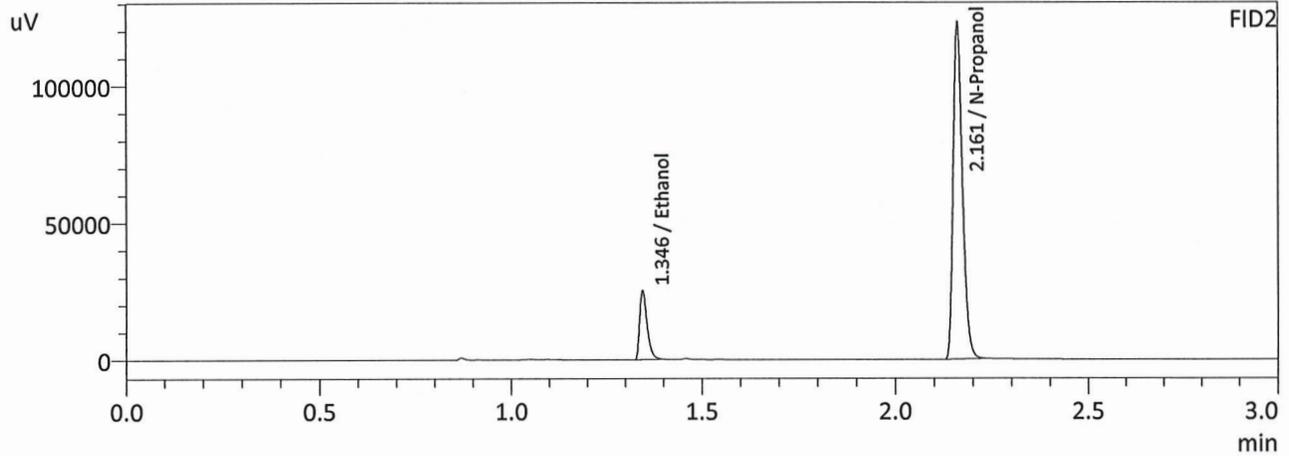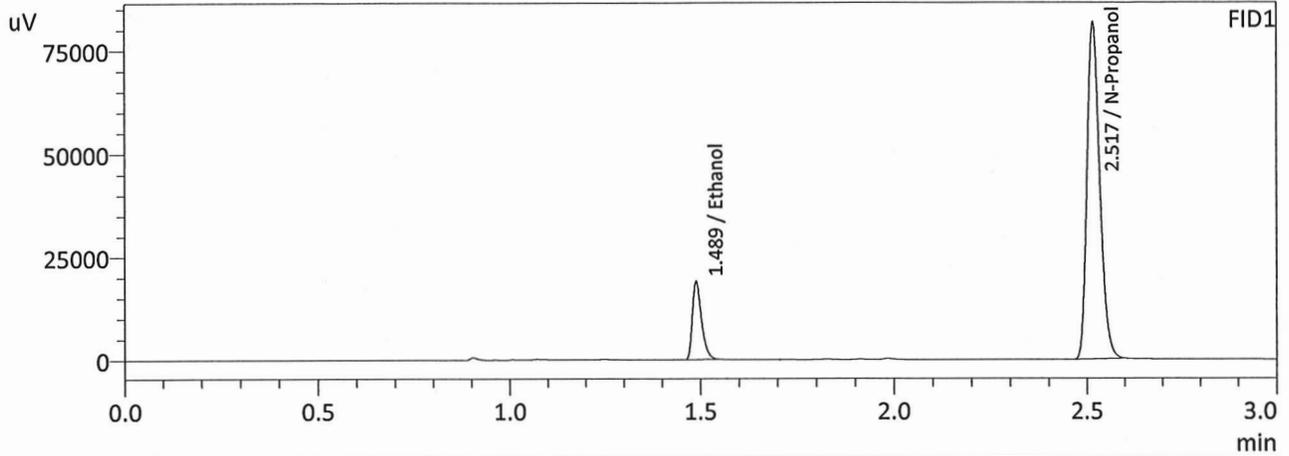

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0807	31611	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	190491	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0808	33838	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	204221	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

JG

VOLATILES DETERMINATION CASEFILE WORKSHEET

VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: 0.08 QA

Analysis Date(s): 4/17/2024 3:23:16 PM(-06:00)

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0837	0.0840	0.0003	0.0838	0.0029	0.0823
(g/100cc)	0.0804	0.0814	0.0010	0.0809		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.



Refer To Instrument Method: ALCOHOL_240411JG.gcm

Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.082	0.077	0.087	0.005

Reported Results	
0.082	

Calibration and control data are stored centrally.

Sample Name : QC-1-1-B
 Laboratory : Meridian
 Injection Date : 4/17/2024 3:14:31 PM
 Vial # : 4
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

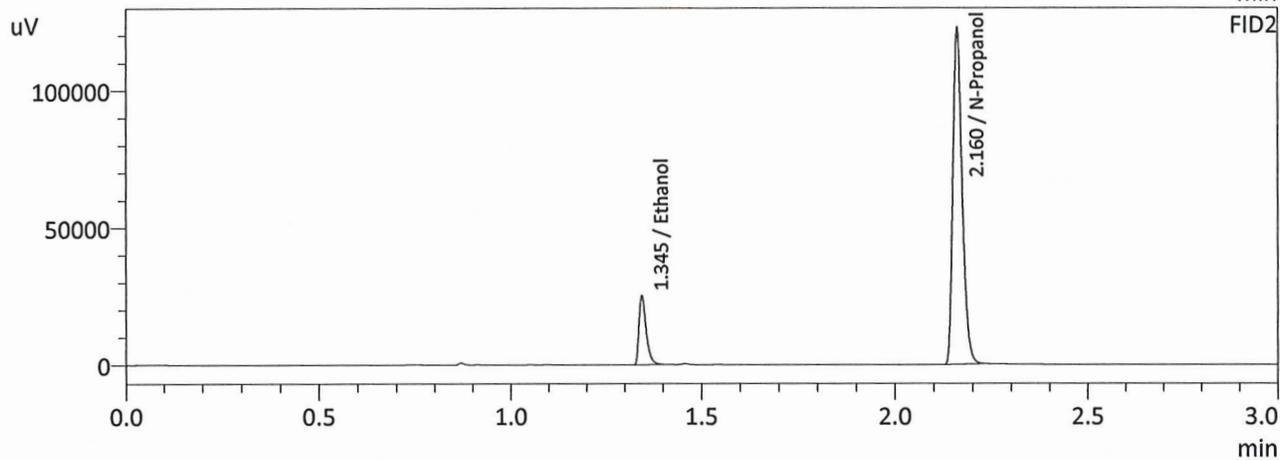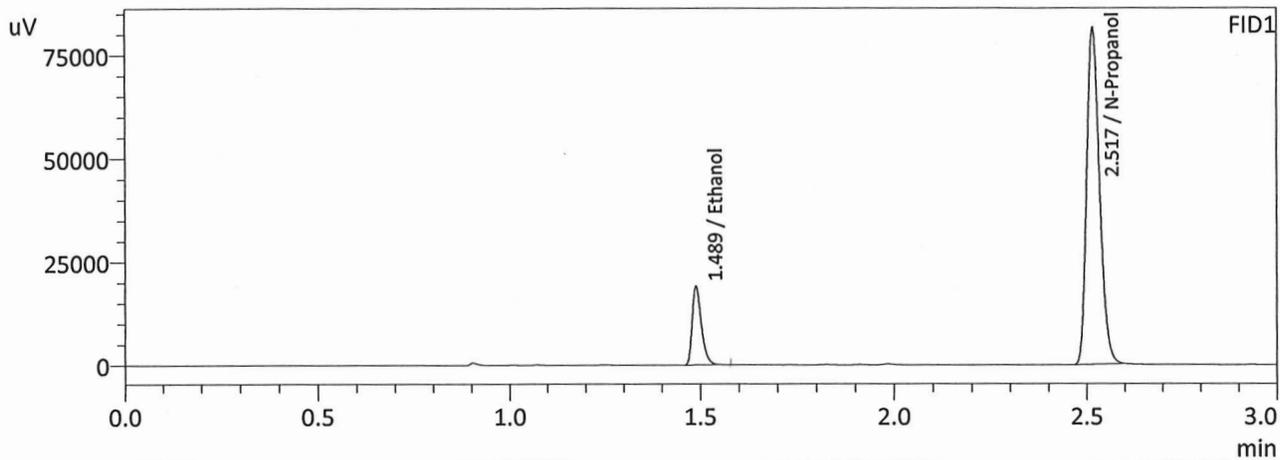



FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0810	31639	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	190068	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0810	33811	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	203515	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

Sample Name : 0.08 QA
 Laboratory : Meridian
 Injection Date : 4/17/2024 3:23:16 PM
 Vial # : 5
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

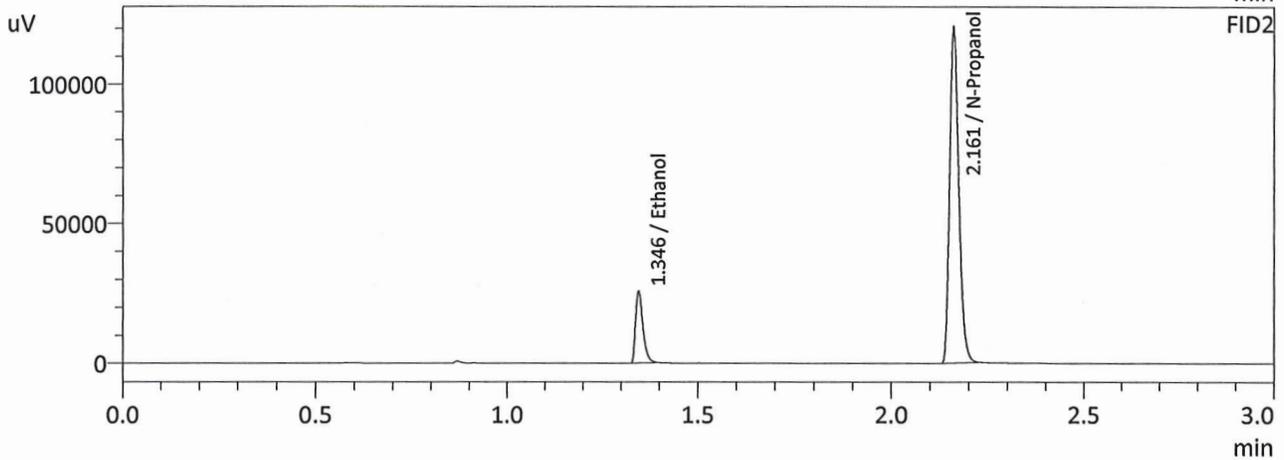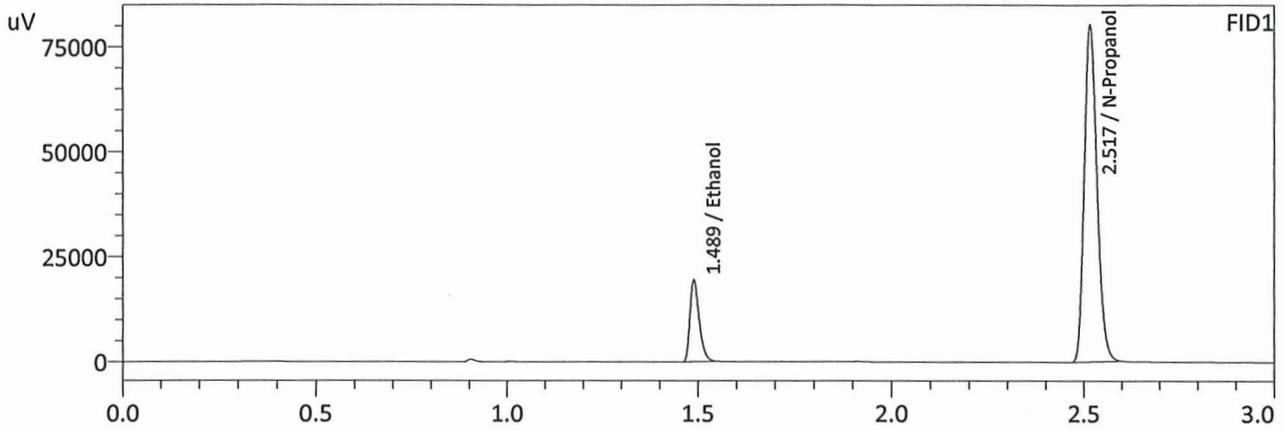

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0837	32298	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	187246	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0840	34627	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	200437	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

JG

Sample Name : 0.08 QA-B
 Laboratory : Meridian
 Injection Date : 4/17/2024 3:30:36 PM
 Vial # : 6
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

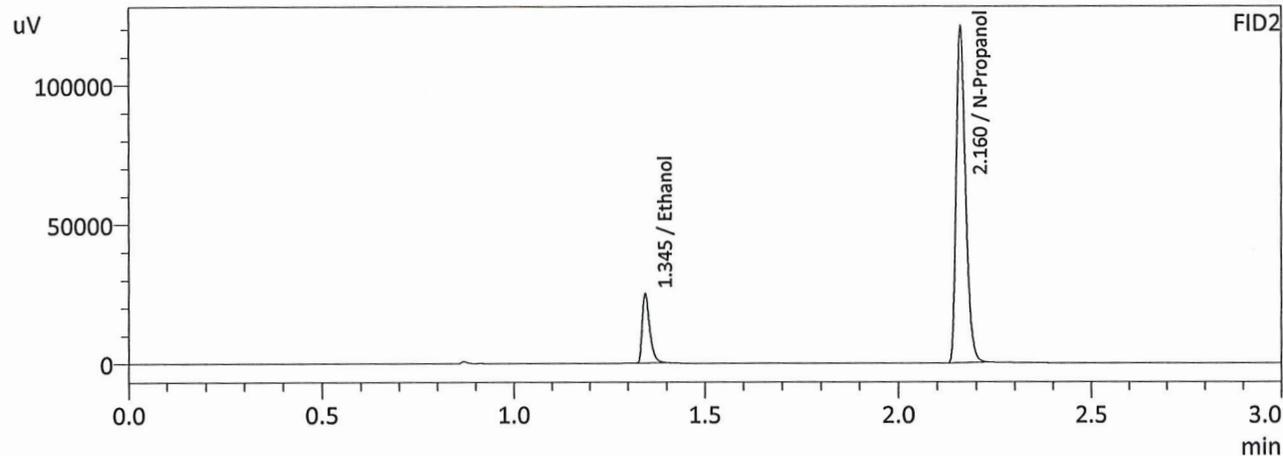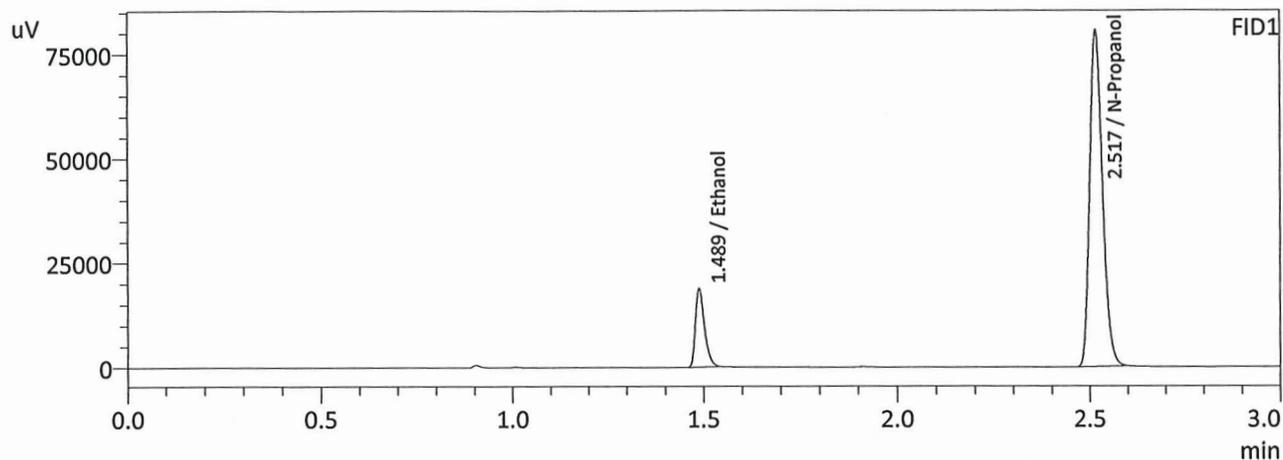

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0804	31035	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	187800	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0814	33561	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	200789	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

36

VOLATILES DETERMINATION CASEFILE WORKSHEET

VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-2-1

Analysis Date(s): 4/17/2024 6:04:27 PM(-06:00)

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.2076	0.2084	0.0008	0.2080	0.0034	0.2097
(g/100cc)	0.2110	0.2119	0.0009	0.2114		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL_240411JG.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.209	0.198	0.220	0.011

Reported Results	
0.209	

Calibration and control data are stored centrally.

JG

Sample Name : QC-2-1
 Laboratory : Meridian
 Injection Date : 4/17/2024 6:04:27 PM
 Vial # : 25
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

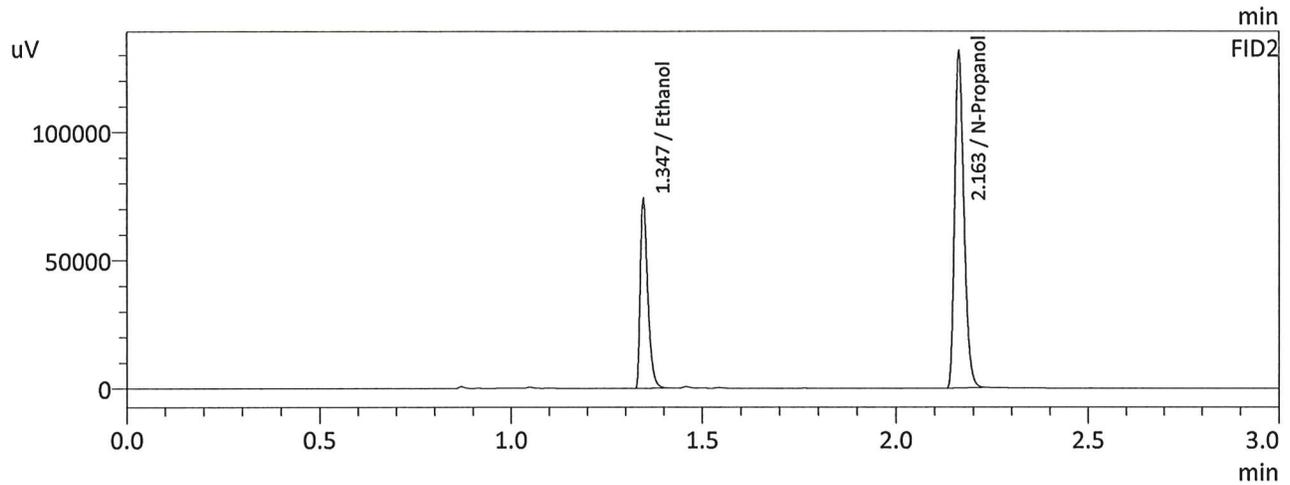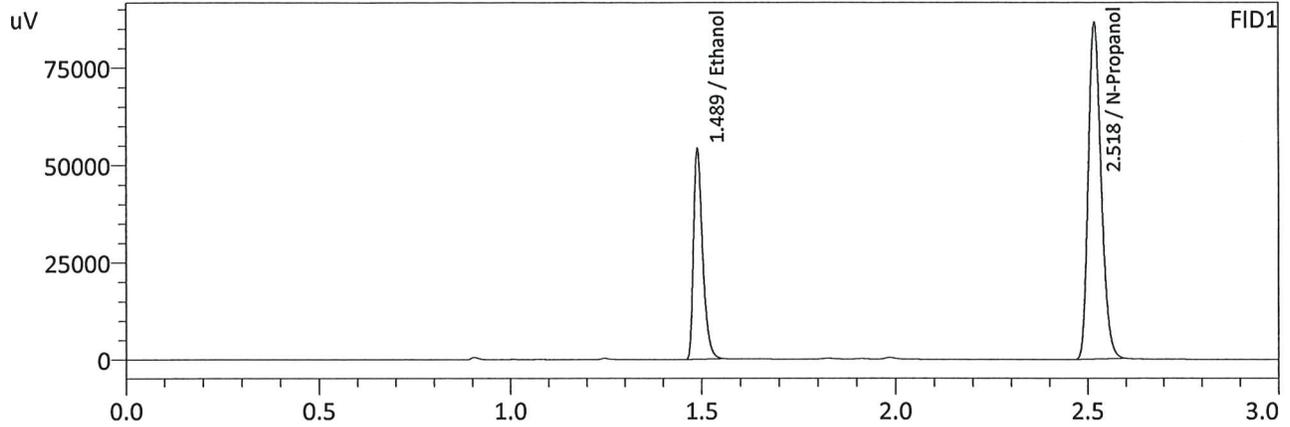



FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2076	89929	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	202387	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2084	97800	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	217997	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

Sample Name : QC-2-1-B
 Laboratory : Meridian
 Injection Date : 4/17/2024 6:11:56 PM
 Vial # : 26
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

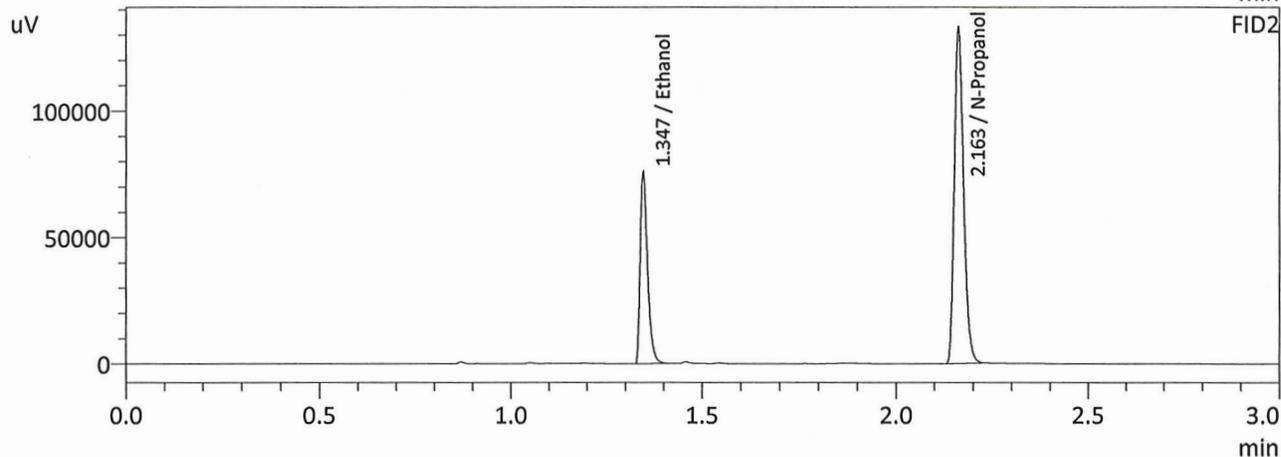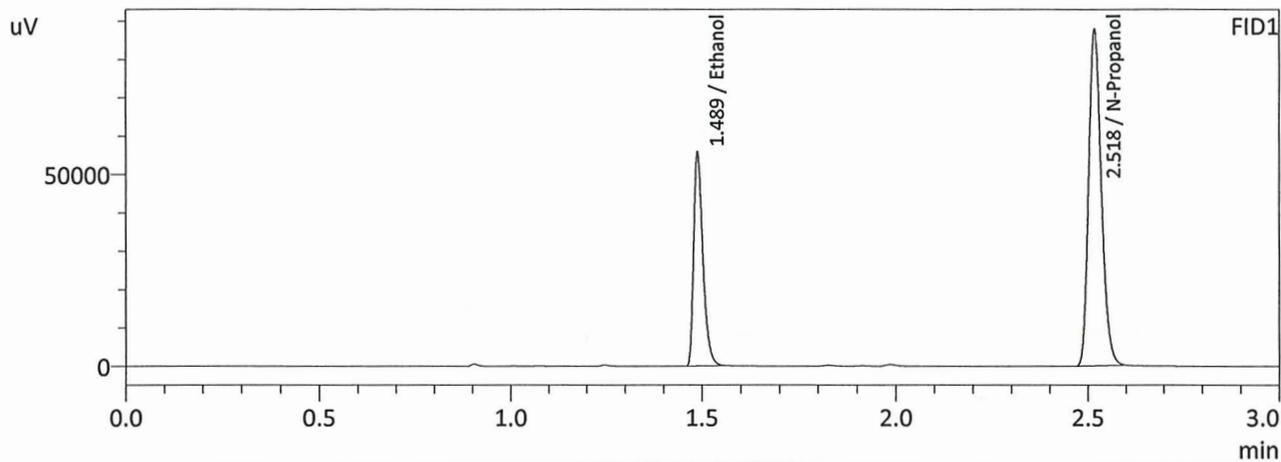

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2110	92626	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	205087	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.2119	100631	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	220519	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

26

VOLATILES DETERMINATION CASEFILE WORKSHEET

VOLATILES DETERMINATION CASEFILE WORKSHEET

Laboratory No: QC-1-2

Analysis Date(s): 4/17/2024 7:10:23 PM(-06:00)

	Column 1 FID A	Column 2 FID B	Column Precision	Mean Value	Sample A-B Difference	Over-all Mean
Sample Results	0.0846	0.0850	0.0004	0.0848	0.0013	0.0854
(g/100cc)	0.0859	0.0864	0.0005	0.0861		

Analysis Method

Refer to Blood Alcohol Method #1

Instrument Information

Instrument information is stored centrally.

Refer To Instrument Method: ALCOHOL_240411JG.gcm



Reporting of Results	Uncertainty of Measurements (UM%): 5.00%		
Overall Mean (g/100cc)	Low	High	5 % of Mean
0.085	0.080	0.090	0.005

	Reported Results
	0.085

Calibration and control data are stored centrally.

JK

Sample Name : QC-1-2
 Laboratory : Meridian
 Injection Date : 4/17/2024 7:10:23 PM
 Vial # : 33
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

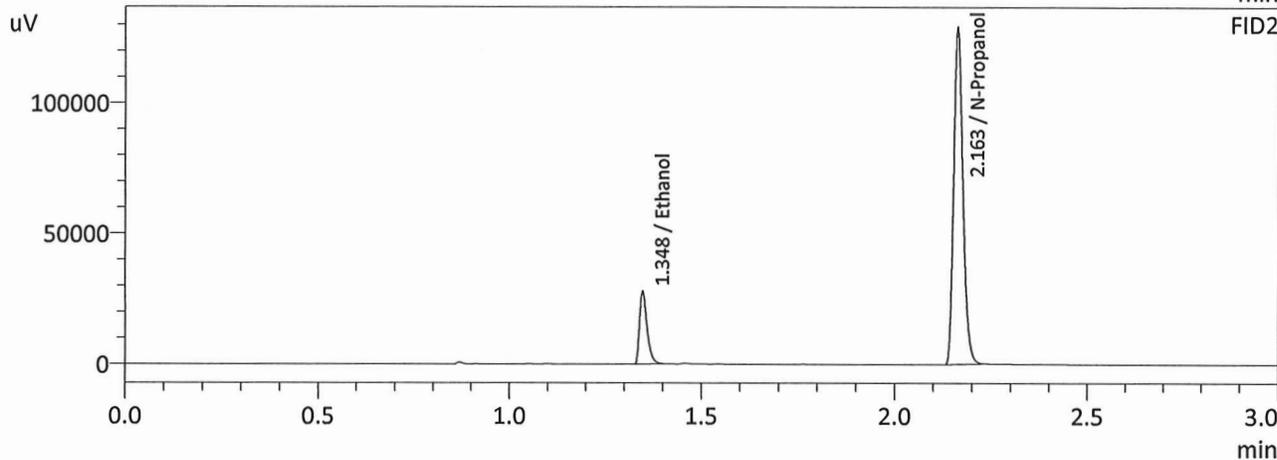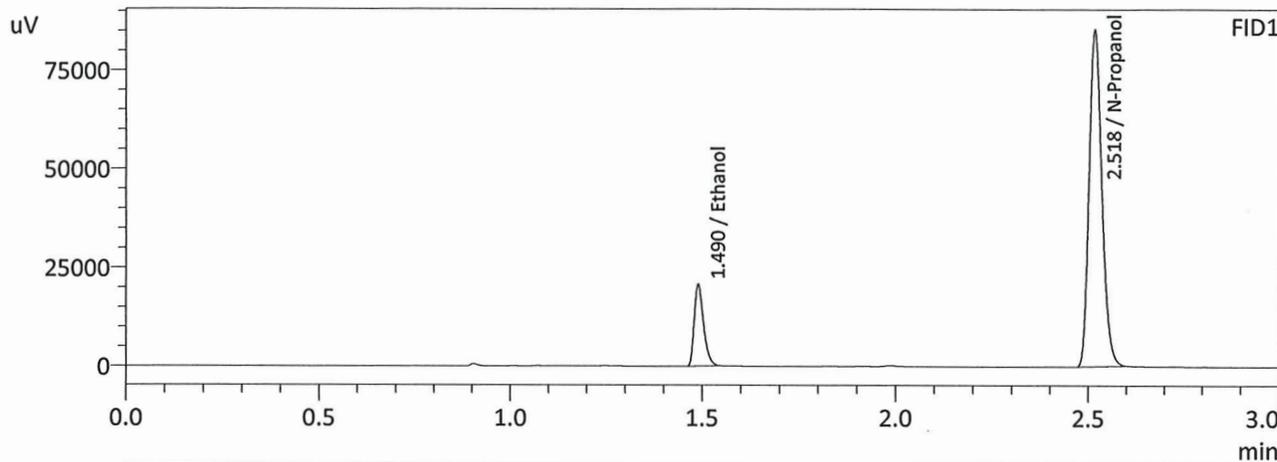

FID1



Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0846	34752	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	199251	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0850	37547	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	214586	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

CC

Sample Name : QC-1-2-B
 Laboratory : Meridian
 Injection Date : 4/17/2024 7:17:21 PM
 Vial # : 34
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

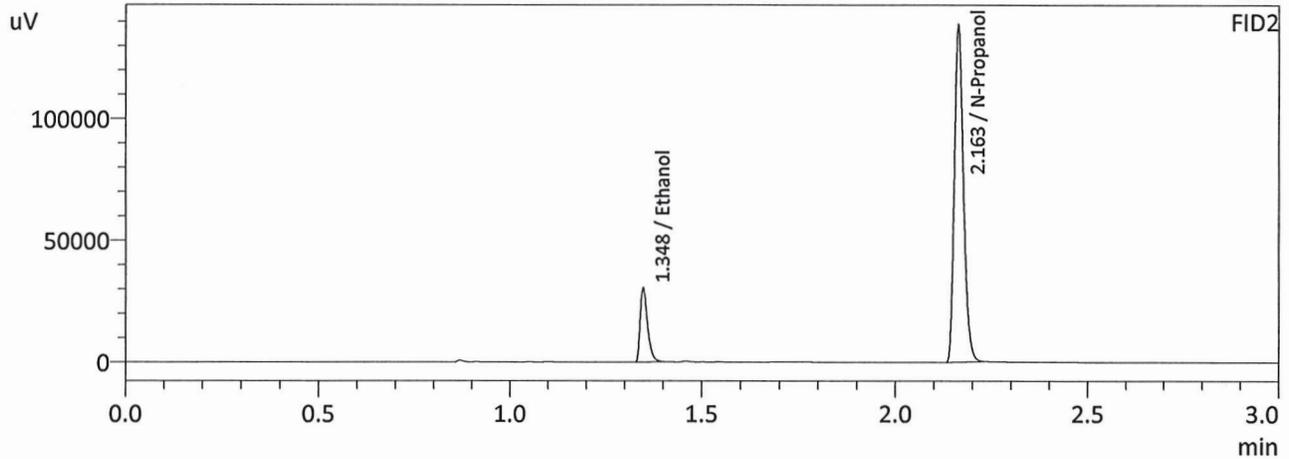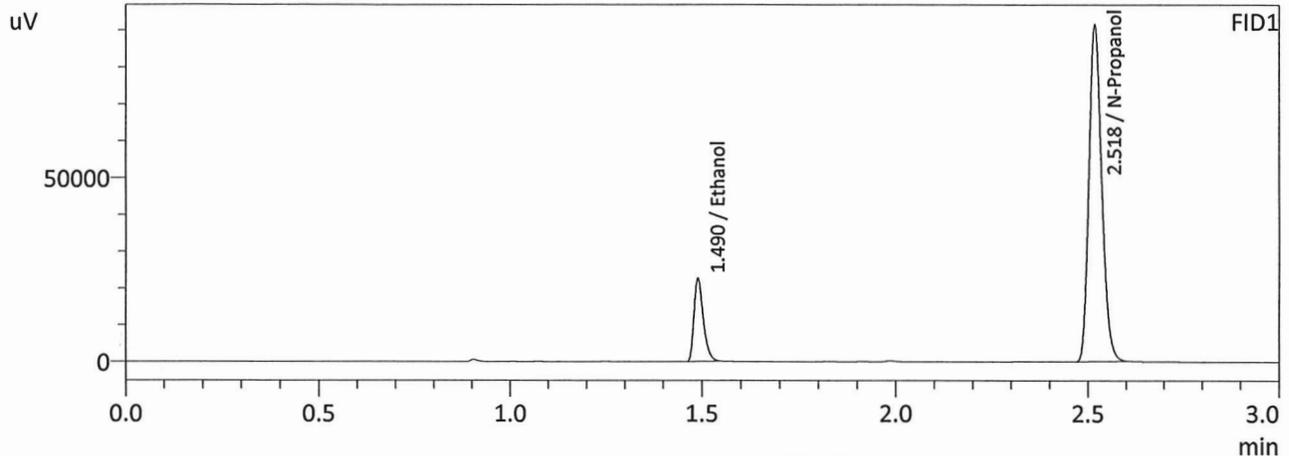



FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0859	37813	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	213180	g/100cc
Fluor. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	0.0864	40925	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	229807	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

Sample Name : ISTD BLK 2
 Laboratory : Meridian
 Injection Date : 4/17/2024 7:25:36 PM
 Vial # : 35
 Method Filename : Default Project - ALCOHOL_240411JG.gcm
 Instrument #GC/HS : C12255750548 / C12595800409

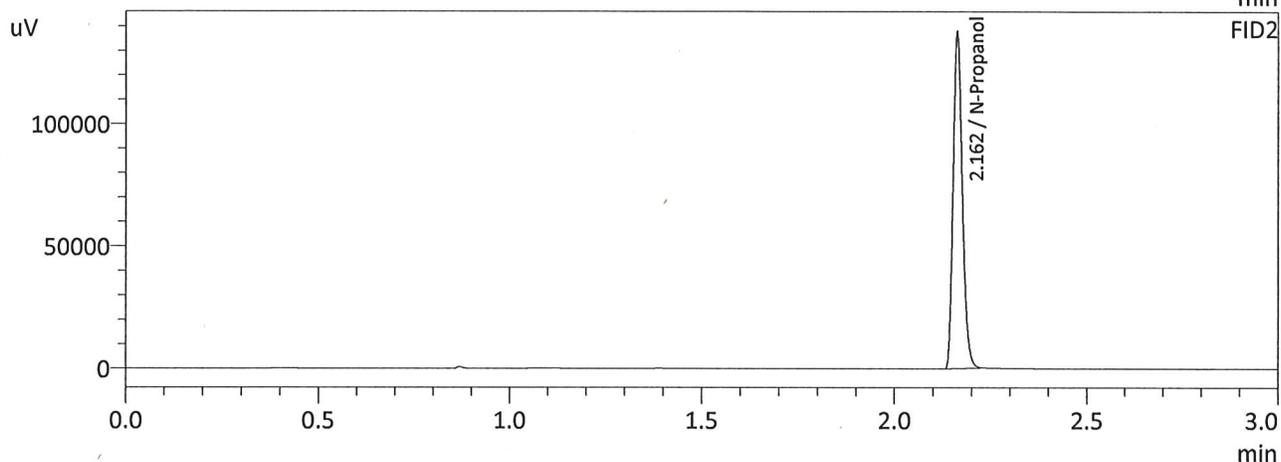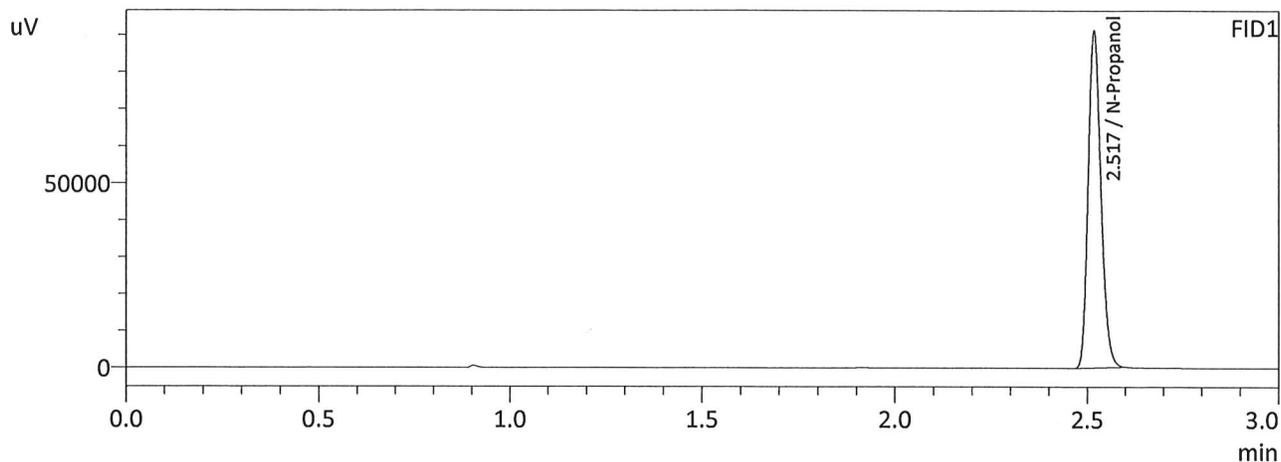

FID1

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
Acetone	--	--	g/100cc
N-Propanol	0.0000	212595	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

FID2

Name	Conc.	Area	Unit
Methanol	--	--	g/100cc
Ethanol	--	--	g/100cc
Acetone	--	--	g/100cc
Isopropyl Alcohol	--	--	g/100cc
N-Propanol	0.0000	228825	g/100cc
Flour. Hydrocarbon(s)	--	--	g/100cc

06

Meridian Blood Alcohol Analysis Batch Table

Shimadzu GC-2030 Serial #C12255750548
 Shimadzu HS-20 Serial #C12595800409
 Lab Solutions Database Software Ver. 6.111
 Copyright (C) 2008-2020 Shimadzu Corporation

Vial#	Sample Name	Sample Type	Level#	Method File
1	ISTD BLK 1	0:Unknown	0	ALCOHOL 240411JG.gcm
2	ED VOLATILES FN 0604	0:Unknown	1	ALCOHOL 240411JG.gcm
3	QC-1-1	0:Unknown	0	ALCOHOL 240411JG.gcm
4	QC-1-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
5	0.08 QA	0:Unknown	0	ALCOHOL 240411JG.gcm
6	0.08 QA-B	0:Unknown	0	ALCOHOL 240411JG.gcm
7	M2024-1349-2	0:Unknown	0	ALCOHOL 240411JG.gcm
8	M2024-1349-2-B	0:Unknown	0	ALCOHOL 240411JG.gcm
9	M2024-1430-1	0:Unknown	0	ALCOHOL 240411JG.gcm
10	M2024-1430-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
11	M2024-1431-1	0:Unknown	0	ALCOHOL 240411JG.gcm
12	M2024-1431-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
13	M2024-1450-1	0:Unknown	0	ALCOHOL 240411JG.gcm
14	M2024-1450-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
15	M2024-1451-1	0:Unknown	0	ALCOHOL 240411JG.gcm
16	M2024-1451-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
17	M2024-1459-1	0:Unknown	0	ALCOHOL 240411JG.gcm
18	M2024-1459-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
19	M2024-1460-1	0:Unknown	0	ALCOHOL 240411JG.gcm
20	M2024-1460-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
21	M2024-1489-1	0:Unknown	0	ALCOHOL 240411JG.gcm
22	M2024-1489-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
23	M2024-1492-1	0:Unknown	0	ALCOHOL 240411JG.gcm
24	M2024-1492-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
25	QC-2-1	0:Unknown	0	ALCOHOL 240411JG.gcm
26	QC-2-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
27	M2024-1539-1	0:Unknown	0	ALCOHOL 240411JG.gcm
28	M2024-1539-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
29	M2024-1556-1	0:Unknown	0	ALCOHOL 240411JG.gcm
30	M2024-1556-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
31	M2024-1563-1	0:Unknown	0	ALCOHOL 240411JG.gcm
32	M2024-1563-1-B	0:Unknown	0	ALCOHOL 240411JG.gcm
33	QC-1-2	0:Unknown	0	ALCOHOL 240411JG.gcm
34	QC-1-2-B	0:Unknown	0	ALCOHOL 240411JG.gcm
35	ISTD BLK 2	0:Unknown	0	ALCOHOL 240411JG.gcm