

IDAHO STATE POLICE FORENSIC SERVICES FOOTWEAR AND TIRE IMPRESSION ANALYTICAL METHOD

HISTORY PAGE

Revision #	Issue Date	History
1	3/29/02	Current methodology used by ISPFSS
2	3/04	Listed ESDA (6.3) as a technique. Modified 6.3 and 6.4 regarding handling of lifts/results.
3	12/04	Modified wording of conclusions. Added digital imaging guidelines.
4	11/4/05	Modified 6.5.3.1, image history; 6.3.1 regarding low resolution digital imaging; 7.0 regarding overlays.
5	1/18/07	Modified 3.1.1.1, use of digital cameras; changed "SOP" to "Analytical Method"; added 10.0, case file contents; merged separate tire and footwear AMs into one; merged reagent formula appendix with the AM.
6	6/26/07	Modified analytical method to include expanded instructions for procedures, field processing, tire impression measurements and digital imaging. Added leucomalachite green back into chemical enhancements section after previous omission.
7	6/27/08	Remove presumptive blood test work step from the following procedures; Amido Black 10.1.5.1, Leucocrystal Violet 10.6.5.1, and Leucomalachite Green 10.7.5.1. Numbering change to accommodate removal.
8	11/8/2010	Changes made to, Definitions, Appendix A, 7, 8.1.1, 8.1.5.1, 8.1.5.2, 8.1.5.6, 9.1.2.2, 9.1.3, 9.1.6, 9.1.6.5, 9.1.7.4, 9.1.8.1, 9.1.8.4, 10.4.3, 10.4.9, 10.4.12, 10.5.7.4, 10.6.7.4, 10.7.4, 10.7.5, 10.7.6.1, 10.7.7.4, 10.8.5.1, 10.9.1, 10.9.5, 11 all subsections, 12.5.1.1, 12.5.1.2, 12.7, 13.2, 13.3, 13.5.2, 14.3.1.2, 15.3, 15.5.3, 15.5.4, 15.8.10, 15.9.2.1, 15.9.2.2, 15.9.2.3, 15.9.2.4, 15.11.1.4, 15.13
9	6/17/2013	Changes made to; 2, 3 Definitions, 5.6, 10.9 Traditional Film Development (Black & Whit) Removed, 11.4.4.4, 11.5.3, 11.5.4, 11.6.4.1, 11.6.4.2, 11.6.5, 11.6.6, 12.4.5.6.7, 12.4.6.4, 14.2.3.1, 14.2.4.1.1, 14.3.1.3, 14.3.1.3.1, 14.3.1.3.2, 14.3.1.3.3, 14.3.1.3.4, 14.3.1.3.5, 14.3.1.3.6, 14.3.1.3.7, 14.4, 14.6, 15 Field Processing removed, Appendix A

**IDAHO STATE POLICE FORENSIC SERVICES
FOOTWEAR AND TIRE IMPRESSION ANALYTICAL METHOD**

TABLE OF CONTENTS

i HISTORY PAGE AND APPROVAL

ii TABLE OF CONTENTS

1 BACKGROUND/SCOPE

2 GENERAL REFERENCES

3 DEFINITIONS

4 RESPONSIBILITIES

5 EVIDENCE HANDLING PROCEDURES

6 GENERAL FOOTWEAR AND TIRE IMPRESSION PROCESSING

7 QUICK REFERENCE SEQUENTIAL PROCESSING GUIDE

8 LIGHT BASED METHODS

8.1 Alternate Light Source (ALS)

9 PHYSICAL METHODS

9.1 Iodine Fuming/ Benzoflavone

9.2 Lifting Methods

9.2.4 Handprint, Tapes and Gel Lifters

9.2.5 Casting Compounds (e.g. Mikrosil)

9.2.6 Casting – Dental Stone

9.2.7 Electrostatic Dust Lifts (EDL)

9.2.8 Electrostatic Detection Apparatus (ESDA)

9.3 Powder Detection Methods

9.4 Small Particle Reagent (SPR)

9.5 Creating Test Impressions

9.5.5 Procedure – Footwear

9.5.6 Procedure - Tires

10 CHEMICAL METHODS

10.1 Amido Black (AB)

10.2 Ammonium Thiocyanate (AT)

10.3 7,8-Benzoflavone

10.4 Cyanoacrylate Ester (CE)

- 10.5 8-Hydroxyquinoline (8HQ)
- 10.6 Leucocrystal Violet (LCV)
- 10.7 Leucomalachite Green (LMG)
- 10.8 Physical Developer (PD)

11 DIGITAL IMAGING PROCEDURE

12 FOOTWEAR AND TIRE IMPRESSION EXAMINATION AND COMPARISON

- 12.4 Procedure
- 12.5 Distribution of Footwear and Tires

13 FOOTWEAR AND TIRE SEARCH

- 13.2 Shoe Search
- 13.3 Tire Search
- 13.4 Search Procedure
- 13.5 Tire Track Measurement Vehicle Computer Search

14 CASE WORK DOCUMENTATION AND REPORT WRITING

Appendix A Instrument Maintenance

1 BACKGROUND/SCOPE

- 1.1 Footwear identification can be dated back to a homicide investigation in 1786. The initial interest in tire marks began with traffic investigation in the late 1920s. One of the earliest uses of tire impression evidence was in England in the early 1940s.
- 1.2 The discipline of Footwear and Tire Impression analysis is the process of determining whether a shoe or tire was the source of an impression.
- 1.3 It is a discipline based on the comparison of outsole or tread design, size, wear patterns and individual characteristics between an unknown or questioned impression and a known shoe or tire.
- 1.4 Footwear evidence may provide the type, make, description and approximate size of a shoe, as well as the number of suspects, sequence of events, and points of entry and exit. Tire track evidence may provide type, make, and approximate size of a tire. It may also indicate the type of vehicle used.
- 1.5 Footwear and tires may have individual characteristics produced as part of the manufacturing process. They also become worn and develop individual characteristics with use.
- 1.6 This Analytical Method defines both technical procedures for processing the majority of evidence encountered by the footwear and tire track impression analysts and comparison methodology.
 - 1.6.1 These methods will describe procedures and techniques that are routinely used in the examination of evidence.
 - 1.6.2 These methods cannot be expected to address each and every situation or type of evidence encountered.
 - 1.6.3 The individual analyst must exercise sound judgment in selecting the methods which will best suit the requirements of the evidence submitted in a specific case; therefore, the procedures are designed to accommodate the majority of evidence encountered.

2 GENERAL REFERENCES

Idaho State Police Forensic Services – Quality Manual ISO/IEC 17025:2005 Compliant.
Section 2.0 NORMATIVE REFERENCES

The Scientific Working Group on Shoeprint and Tire Tread Evidence Technology
(SWGTTREAD) - *SWGTTREAD documents are* <http://www.swgtread.org/>

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

3 DEFINITIONS

Abrasion

The progressive wearing away of a tire in service.

Adhesive Lifter

Any adhesive coated material or tape used to lift impressions.

Aging

Deterioration (physical and chemical) of rubber through oxidation over a period of time.

Alternate Light Source

Ultraviolet light source which may have variable wavelengths, used to visualize footwear or tire track impressions.

Ambient Light

The available or existing light that surrounds the object being photographed.

Attenuated Light

Supplemental light that is added when photographing luminol enhancements. The light is sufficient to allow for the exposure of the object being photographed but not to significantly interfere with the photography of luminol.

Analysis

The methodical examination of footwear and tire track impressions; separation into parts so as to determine the nature of the whole.

Aspect Ratio

A numerical term that expresses the relationship between the cross-section height and the cross-section width of a tire. An aspect ratio of 75 means that the tire is approximately 75% as high as it is wide.

Balance

Equal weight distribution around the tire.

Band Ply

The first ply inside the tire.

Beads

Layers of cords are turned around a coil or "bead" of wire at the inner edge of the sidewall, allows the casing to be fitted to the wheel rim and provides an air-tight seal in tubeless tires.

Blemish

Surface irregularity or imperfection.

Blister

A separation within the sidewall or liner stock, or between the sidewall carcass or liner and the carcass.

Body

The tire structure except the tread and sidewall rubber. Also known as carcass.

Brand Name

The company name, which appears on the sidewall, not necessarily the name of the company that manufactured it.

Buffing Rib

Raised rib at the junction of black and white rubber on a white sidewall tire. Also called a scuffing rib.

Buttress

Portion of tread running around the shoulder and blending into the sidewall.

Camber

Tilt of front wheel; inward tilt at the top of the wheels is negative camber (excessive wear on outer tire tread shoulder edge) and inward tilt at the bottom of the wheels is positive camber (excessive wear on inner tire tread shoulder edge).

Casing

The whole tire.

Cast

The result of the filling of a three-dimensional impression with an appropriate material.

Casting

The filling of a three-dimensional impression with material that takes on and retains the characteristics, which were left in that impression.

Casting Material

Dental stone, sulfur, or other suitable materials specifically used for and generally accepted to accurately reproduce impressions.

Checking

Minute cracking in the surface of the rubber caused by aging and oxidation.

Circumferential Cracks

Cracks in the tire running parallel to the beads, usually consists of cracks in the circumferential tread grooves.

Class Characteristics

Manufactured features that are shared by two or more shoes or tires.

Combined Class Characteristics

The combination of two or more independent class characteristics.

Coaxial Light

Illumination from the precise direction of the imaging lens, either through the lens or with a beam-splitter in front of the lens.

Comparison

The observation of two areas of footwear or tire track impressions for finding similarities and/or differences.

Cords

The tire casing or carcass is composed of flexible filaments of natural textile, synthetic polymer, glass fiber or steel embedded in or bonded to a matrix of natural or synthetic rubber.

Correspond

A word used to describe agreement of characteristics in the context that reflects the shoe or tire is capable of having produced a certain feature such as design, general wear, etc.

Cracking Tread or Groove

Splitting in grooves caused by excessive strain.

Cross-Section Size

External sidewall-to-sidewall measurement of inflated tire exclusive of elevations due to raised lettering or protector bands.

Crown

Section between shoulders of the tire.

Deformable Impression

An impression that causes the surface to deform, either permanently or temporarily. Permanent deformable impressions would include those impressions in sand, soil, and snow whereas temporarily deformed impressions would include those on skin or carpet.

Dental Stone

A gypsum product, similar to plaster of paris, but with different properties due to the manufacturing process.

Design

The manufactured pattern of a shoe outsole or tire tread.

Difference

A characteristic or feature that is so strong and reliable that it, in itself, indicates non-identity. Usually a difference will be a different class characteristic, such as the specific design or specific physical size of the design. Normal variations in the impression process, the absence of cuts in a questioned impression that appear on the shoe or tire, or the normal advancement of wear with time do not necessarily constitute a comparative difference.

Discrepancy/Dissimilarity

A characteristic that has the appearance of being different but lacks sufficient detail to confirm as such.

Distortion

Variances in the reproduction of impressions caused by pressure, movement, force, contact surface, etc. Distortion is not a discrepancy and is not a basis for exclusion.

DOT Number

Department of Transportation serial number assigned to every tire sold in the United States which gives information regarding the manufacturer of the tire, tire size, and date of manufacture of the tire.

Dry Origin Impression

Formed under dry conditions such as dry dust and dry residue impressions.

Elasticity

The ability of skin to recover from stretching, compression, or distortion.

Electrostatic Detection Device (EDD) or Analyzer (ESDA)

An instrument used primarily to detect indented writing on documents, which can also be used to detect footwear impressions on paper items.

Electrostatic Dust Lifting Device (EDL)

An instrument that utilizes electrostatic charges as a means of transferring dry origin impressions from a surface to a film.

Element

Unique and separate tire tread design areas distinguished by grooves, slots and sipes, that are the actual surface contact portion of the tread design. Also known as tread blocks. May also refer to design characteristics on footwear outsoles.

Elimination/Exclusion Prints

Exemplars of footwear or tires belonging to persons known to have had access to the item or area examined for impressions.

Enhancement

Rendering an impression clearer or more visible through physical, electronic, photographic or chemical means.

Erroneous Identification

The incorrect determination that two footwear or tire track impressions originated from the same source.

Evaluation

The determination of the significance, value, or clarity of an impression by careful observation and study.

Examination Quality Photographs

High quality photographs taken with a scale for use in the physical comparison of footwear and tire impressions with known footwear and tires.

Exclusion

An opinion by an examiner that is the highest degree of non-association expressed in footwear and tire impression examinations that the particular known footwear or tire was not the source of and did not make the impression.

Feather Edge

Sharp knife-like feathering along one edge of tire tread rib or blocks, due to scrubbing action of misaligned wheels.

Feathering or Schallamach Pattern

A wear pattern that is the result of frictional abrasive forces which have crossed the shoe sole perpendicular to the direction of travel. This results in a unique pattern that bears some resemblance to ridge characteristics and bifurcations of a fingerprint pattern and which continuously evolves and changes.

Fixative

A spray or powder applied to a three-dimensional impression prior to casting. It assists with the release of the substrate material from the cast. It may also refer to the application of a reagent to an impression that assists with the adherence of the impression to the substrate prior to chemical enhancement.

Gelatin Lifter

Gelatin laid on a pliable backing that can be used to lift impressions.

General Wear

The overall condition of a shoe outsole or tire tread, ranging from new to extremely worn, related to its degree of use. General wear is a class characteristic that may be used to include or exclude shoe outsoles and tire treads based on similar or different degrees and positions of wear.

Gouging

Chisel-like action of rocks or stubble on tires.

Grooves

Furrows that separate ribs and run around the tire circumference, enhancing tread performance.

Heat Breaks

Failure of fabric due to excessive heat generation in tires.

Identification

An opinion by an examiner that is the highest degree of association expressed in footwear and tire impression examinations that the particular known footwear or tire was the source of and made the impression.

Identifying Characteristic or Individual Characteristic

A particular characteristic resulting from something randomly added or removed from a surface (shoe sole or tire tread) that contributes to making that surface unique. A particular feature individual to a specific shoe or tire that resulted from a random occurrence. Examples would include cuts, burns and scratches.

Impression

The product of direct physical contact of an item resulting in a transfer and retention of characteristics of that item.

Inconclusive

The inability to either individualize or exclude an area of impression detail.

Individual Characteristic

A particular feature individual to a specific shoe or tire that resulted from a random occurrence.

Insufficient Detail

Features which fall short of allowing the confirmation of certain class or individual characteristics.

Intermediate Rib

A rib located between the center and outer ribs.

Knock-off

A copy of a popular brand of footwear.

Known Print (Footwear or Tire)/Exemplar/Test Impression

A recording of the outsole or tread pattern with black ink, electronic imaging, photography, or other medium on a contrasting background. An object of known origin that is compared to a questioned impression.

Lift

An adhesive or other medium on which recovered impression detail is preserved.

Lug

Lateral blocks or discontinuous rows of tread rubber usually located on the outer edges of the tread design.

Manufacturer

The company that made the shoe or tire.

Matrix

The substance that is deposited by the footwear or tire.

Misalignment

Maladjustment of one or more parts of steering mechanism of vehicle causing unsatisfactory operation and tire wear.

Missed Identification

The failure to make an identification (individualization) when, in fact, both impressions are from the same source.

Mold

The cavity in which shoes or tires are made that contains the design for the outsole or tread.

Mold Characteristics

Variations in the design pattern of an individual mold that is transferred to all shoes outsoles or tires made in that mold.

Non-Porous

Non-absorbent.

Non-Skid Depth

Depth of the tire tread design providing skid protection that is measured from the top of the ribs to the base of the grooves.

Oblique Light

Illumination from a light source that is at a low angle of incidence, or even parallel, to the surface of the item. (also known as side lighting)

Observable Characteristic

Class and individual features that can be visually confirmed or otherwise detected.

Patterns

The designation of footwear and tire track impression designs as well as outsoles and tire treads into basic categories of general shapes.

Photo Log

A written record of photographs taken at the crime scene.

Pitch Sequence

The elements on the tire tread that are varied to reduce the noise of the tire and improve tire performance under the rated conditions. (Also known as **noise treatment**)

Plies

Cords are fabricated in a parallel fashion into several multifilament layers or plies.

Polarized Lighting

Lighting consisting of light rays with a single propagation direction and a single vibration direction. Polarized light is produced by the use of a polarizing filter.

Porous

Absorbent, a material or surface having pores.

Physical Size and Shape

The size, shape, spacing and relative positions of the outsole design components and tire tread blocks (not the same as the manufacturer's shoe or tire size). The physical size of design is a class characteristic.

Preserved

Casting, photography, lifting, or other method used to capture impressions for further examination.

Profile

Shoulder to shoulder arc of tire tread cross-section.

Qualified Analyst

An individual who has completed the internal training program, passed competency testing, and been approved to perform case work.

Questioned Impression

An impression made by a shoe or tire whose origin is unknown. (Also known as an unknown impression).

Randomly Acquired Characteristic

Features on a footwear outsole or tire tread resulting from random events including, but not limited to: cuts, scratches, tears, holes, stone holds, abrasions and the acquisition of debris. The position, orientation, size and shape of these characteristics contribute to the uniqueness of a shoe outsole or tire tread. Randomly acquired characteristics are essential for an identification of a particular shoe or tire as the source of an impression.

Reagent

Substance used in a chemical reaction to detect, examine, measure, or produce other substances.

Recapping

A process known as “top capping” in which rubber is only applied to the tire tread surface.

Relative Position

Proximity of characteristics to each other.

Retreading

Renewing the tread on a tire by buffing the old surface and applying a new tread. It differs from recapping in that new rubber is applied to both the tread surface and shoulder area of the tire.

Ribs

Sections of tread rubber that run in rows circumferentially around the tire. The number of ribs across the width of the tread is often used when describing the tread design.

Rim

The metal support for the tire or tire and tube assembly on the wheel.

Sequential Processing

Use of a series of development methods in a specific order to maximize development of impression detail.

Shoulder

The outer edges of the tire connecting with the sidewall, may be molded round or square.

Sidewalls

The side of the tire between the rim and the tread. Sidewalls may be white or black and have information molded on them in the form of raised letters and numbers and other markings.

Similar Class Characteristics

Conclusion statement meaning that an impression and known shoe share some degree of resemblance without inferring total correspondence.

Sipes

Narrow slits incorporated into the tread design to provide improved traction and channel water off the surface of the ribs and into the grooves. Sipes that do not extend to the edge of a rib or element are called “contained sipes”. Sipes that extend to the edge of a rib or element are known as “open sipes”. Sipes are also known as blades or kerfs.

Slots

Narrow grooves (lateral slots or cross-grooves) that run in a lateral or diagonal direction across the tread, usually from one circumferential groove to another.

Spot Repair

A repair made by vulcanizing rubber to a tire without replacing cord.

Stock Solution

Concentrated solution diluted to prepare a working solution.

Stone Drilling

A condition in which a stone has penetrated the tread rubber or fabric of a tire, caused by failure of tires to eject stones picked up by the grooves.

Stud

A hard metal or plastic rivet inserted into a pinned tire for the purpose of improving winter traction.

Stud Pattern

Pattern of holes in the tread design which are strategically placed to hold small metal studs added to the tread area to improve traction.

Substrate

Surface upon which an impression is deposited.

Textured Sidewall

A tire sidewall which has been cured into a web-like pattern to decrease weather cracking.

Three-Dimensional Impression

An impression with the dimensions of length, width and depth.

Tie-Bars

Small connecting bars of rubber recessed in the grooves and sipes to add support to the elements, lugs and ribs which comprise the tread design.

Tire Balancing

Adding external weights to compensate for unequal distribution of tire weight.

Tire Growth

Slight increase in tire size over time due to inflation pressure and operating circumstances.

Tire Name

Name given to a particular design which distinguishes it from other tires of the same brand, appears in large letters on the sidewall of the tire.

Toe-in

Adjustment of front wheels so they are closer together at the front than at the back. Most wheels are set this way when the vehicle leaves the plant. Too much toe-in creates heavy wear on the outer shoulder of the tire, showing as heavy wear from one direction and sharp feathered edges from the other direction across the tread.

Traction

The grip or friction, pull or bite between the tire and the ground.

Tread

The portion of the tire that covers the plies and belts from shoulder to shoulder and makes contact with the surface.

Tread Wear Indicator

Raised rubber ridges in the grooves that run laterally across the tread which are used to signal when the useful life of the non-skid depth has expired. (Also known as **wear bars**).

Two-Dimensional Impression

An impression which for all practical purposes has only the dimensions of length and width.

Vehicle Track Width

Distance between the center points of the tires from one vehicle side to the other (i.e. from the center point of the passenger side front tire to the center point of the driver side front tire). On a dual axle vehicle, this is the distance from the center points between the dual tires on one vehicle side to the other.

Vehicle Turning Diameter

The diameter of the smallest circle that is made by the outer edge of the outer most front tire in a turn.

Vehicle Wheel Base

Distance between the front and rear axles of a vehicle.

Verification

Confirmation of an analyst's conclusion by another qualified analyst.

Wear

Erosion of the shoe outsole or tire tread due to friction and abrasive forces with the ground.

Wear Characteristics

Changes in the surface of the shoe outsole or tire tread that are observable in the impression and/or the known shoe or tire and that reflect the erosion of the surface of the outsole or tread.

Wear Pattern

The position of wear on the shoe or tire.

Weather Checking

Fine hair-like cracks in the surface of the tire rubber, caused by oxidation and other atmospheric effects.

Wet Origin Impression

Formed under wet conditions, including impressions deposited by residues of blood, grease, mud, and other wet substances.

Wheel Alignment

Adjustment of wheels for proper vehicle operation.

Working Solution

Solution at the proper dilution for processing.

4 RESPONSIBILITIES

- 4.1 The Latent/Impression Evidence Section Supervisor is responsible for ensuring that personnel adhere to established analytical methods, safety practices, and laboratory policies and procedures.
- 4.2 The Latent/Impression Evidence Section Supervisor shall ensure that analyst's training records are on file in the impression evidence section.
- 4.3 Individual analysts are responsible for adherence to established analytical methods, safety guidelines, and laboratory policies and procedures.
- 4.4 Footwear and Tire Impression analyst duties include, but are not limited to:
 - 4.4.1 Development of impressions;
 - 4.4.2 Documentation of visible or developed impressions;
 - 4.4.3 Obtaining known test impressions from footwear and tires;
 - 4.4.4 Analysis and comparison of impressions;
 - 4.4.5 Verification of compared impressions;
 - 4.4.6 Issuing reports of examination activities;
 - 4.4.7 Performing technical and administrative casework reviews;
 - 4.4.8 Responding to crime scenes to the extent to which they are trained;
 - 4.4.9 Satisfactorily completing annual proficiency tests;
 - 4.4.10 Presenting expert testimony in court.

Property of Idaho State Police Forensic Services
Uncontrolled Interim Copy
OBSOLETE DOCUMENT 11/14/2016

5 EVIDENCE HANDLING PROCEDURES

- 5.1 Evidence handling will be in accordance to ISPFS Quality/Procedure Manual Section 5.8 HANDLING ITEMS OF EVIDENCE.
- 5.2 Analysts will place unique case identifiers and initials on items where possible.
- 5.3 Wet or damp evidence will be dried in a secure area.
- 5.4 Footwear and tire impression evidence will be stored at room temperature before and after the examination process unless there are special circumstances.
- 5.4 Casts and other items should be photographed for file documentation, unless a better method of documentation can be utilized.
- 5.5 Footwear and tire impression processing has the potential to irreparably damage items of evidence.
 - 5.5.1 If an item is suspected to have great value (monetary or sentimental), the analyst should contact the submitting agency to explain potential damage and gain verbal approval prior to processing.
- 5.6 RETAINED EVIDENCE
 - 5.6.1 Evidence generated by the Impression Evidence section may be retained for future reference with the approval from the Discipline Leader.
 - 5.6.2 Impression Evidence Section will return evidence and sub items of evidence as (EDL lifts or powder lifts) to the submitting agency. (See ISPFS Quality Manual section 15.8.2.4.1)

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

6 GENERAL FOOTWEAR AND TIRE IMPRESSION PROCESSING

- 6.1 Impression print evidence is processed according to the nature of the substrate (surface) to be processed.
 - 6.1.1 Substrate types include porous and non-porous.
 - 6.1.2 Processing is generally carried out in a sequential manner employing methods appropriate to the substrate type.
- 6.2 Impression evidence is also processed with regards to what the impression matrix may consist of.
- 6.3 For the purpose of this manual, footwear and tire impression methods are divided into three categories; light based methods, physical methods, and chemical methods.
 - 6.3.1 LIGHT BASED METHODS
 - 6.3.1.1 Impressions may be visualized through the use of various angles and wavelengths of light.
 - 6.3.1.2 Visualization of impressions through the use of forensic lighting methods is non-destructive and should be attempted prior to other processing methods.
 - 6.3.2 PHYSICAL METHODS
 - 6.3.2.1 The development of impressions through the use of physical methods does not involve a chemical reaction between the impression and the method used.
 - 6.3.2.2 Physical methods encompass dusting and other discoloration methods often relying on the adhesive quality of certain impressions.
 - 6.3.2.3 The taking of known test impressions from shoes or tires shall be considered a physical method for the purposes of this manual.
 - 6.3.3 CHEMICAL METHODS
 - 6.3.3.1 The development of footwear and tire track impressions through the use of chemical methods occurs because of a chemical reaction between the impression residue components and the reagent.
 - 6.3.3.2 Reagents shall be tested after they are prepared and prior to use. If the same lot of working solution is used multiple times in the same day, only the initial control tests must be noted on the Reagent Preparation Log sheet. Subsequent use of the reagent on the same day shall revert to the prior tests. Control test results shall also be recorded in the case documentation. Reagents made up for one-time use may be documented in the case notes only.

7

QUICK REFERENCE SEQUENTIAL PROCESSING GUIDE

(Methods or combinations of methods should be chosen according to the substrate and/or matrix of the item to be examined that will generate the best possible results.)

7.1 GENERAL EVIDENCE:

7.1.1 POROUS:

1. Visual: White light/Alternate Light Source (ALS)
2. Iodine Fuming
3. Visual: White light
4. 7,8 Benzoflavone
5. Visual: White light
6. 8-Hydroxyquinoline
7. Visual: Ultraviolet light 254nm-365nm
8. Ammonium Thiocyanate
9. Visual: White light
10. Physical Developer
11. Visual: White light

7.1.2 NON-POROUS:

1. Visual: White light
2. Cyanoacrylate Fuming
3. Visual: White light/ALS
4. Powders: Non-luminescent
5. Visual: White light/ALS

7.2 BLOOD EVIDENCE:

7.2.1 POROUS:

1. Visual: White light UV
2. Amido Black, Leucocrystal Violet, or Leucomalachite Green
3. Visual: White light

7.2.2 NON-POROUS:

1. Visual: White light
2. Amido Black, Leucocrystal Violet or Leucomalachite Green
3. Visual: White light

7.3 GLASS:

1. Visual: White light
2. Electrostatic Dust Lift (EDL)
3. Cyanoacrylate Fuming
4. Visual: White light
5. Powders: Non-luminescent
6. Visual: White light

7.4 PAPER/ CARDBOARD:

1. Visual: White light
2. EDL
3. Iodine fuming
4. Visual: White light
5. 7,8-Benzoflavone
6. Visual: White light
7. 8-Hydroxyquinoline
8. Visual: Ultraviolet 254nm-365nm
9. Ammonium Thiocyanate
10. Visual: White light
11. Physical Develop
12. Visual: White light

7.5 HUMAN SKIN:

7.5.1 Living

1. Visual with white light, ultraviolet and infrared light

7.5.2 Deceased

1. Visual with white light, ultraviolet and infrared light
2. Electrostatic dust lift (EDL)

7.6 METAL:

1. Visual: White light
2. EDL
3. Visual: White light
4. Cyanoacrylate fuming
5. Visual: White light
6. Powders: Non-luminescent
7. Visual: White light

7.7 WET SURFACES:

7.7.1 POROUS:

1. Visual: White light
2. Dry to room temperature
3. Visual: White light
4. Physical developer

7.7.2 NON-POROUS:

1. Visual: White light
2. Small Particle Reagent (SPR)
3. Visual: White light

8 LIGHT BASED METHODS

8.1 ALTERNATE LIGHT SOURCE

8.1.1 BACKGROUND:

Alternate light sources (ALS) are portable, multi-waveband, and tunable light sources that are used to enhance or visualize potential items of evidence. Impressions may be composed of various substances such as blood, soil, chemicals or other organic substances that react differently to different wavelengths of light. When a luminescent deposit is excited with a particular wavelength of light, the deposit absorbs the light and re-emits it at a different wavelength. The short-lived light being re-emitted is termed fluorescence. There are several alternate light sources available to analysts that adequately meet the needs described in this manual.

8.1.2 SCOPE:

8.1.2.1 The ALS is used to attempt to create contrast between an impression and the substrate it is on.

8.1.2.2 Fluorescence may occur due to a naturally occurring substance within the impression residue itself (inherent luminescence), may be transferred to the shoe or tire via contamination and re-deposited, or may be chemically induced in impression residue with certain reagents known to exhibit fluorescent properties. Fluorescence of the substrate may also occur.

8.1.3 EQUIPMENT AND MATERIALS:

Alternate light source
Filtered goggles

8.1.4 REAGENTS:

Not applicable

8.1.5 PROCEDURE:

8.1.5.1 Turn the power switch on. Make sure the unit comes to full operating speed.

8.1.5.2 Turn on the lamp switch.

8.1.5.3 Choose the band-width you wish to use. Some models have a variable power dial that may need to be adjusted.

8.1.5.4 Observe evidence with the appropriate wavelength/goggle combination:

< 400nm	yellow or clear UV safe
400-450nm	yellow
450-540nm	orange
540->700nm	red

8.1.5.5 Push the lamp switch to off.

8.1.5.6 Wait for the unit to cool down

8.1.6 ADDITIONAL INFORMATION:

8.1.6.1 Maintenance shall consist of cleaning the exterior of the ALS with a soft cloth dampened with a mild detergent solution and using a cotton swab moistened with glass cleaner to clean the optical filters. Bulbs should be replaced as needed.

8.1.6.2 If an ALS malfunctions, it will be taken out of service until it can be repaired. The ALS shall be tagged indicating that it is out of service. Maintenance, service, etc. will be recorded in the maintenance log.

8.1.6.3 No calibration is required of these units.

8.1.6.4 The manufacturer's operator manuals shall be read prior to using the equipment.

8.1.7 CONTROLS:

Not applicable

8.1.8 SAFETY:

8.1.8.1 As with other electrical appliances, guard against electrical shock. This can be accomplished by ensuring that all connections are proper and that no loose, damaged, or frayed wires exist. Make sure the ALS is unplugged before attempting any maintenance and do not use outdoors if wet conditions exist.

8.1.8.2 The eyes are generally more vulnerable than the skin, and appropriate eye protection must be used to protect them. Permanent eye damage can occur from reflected, refracted, or direct illumination to the eye. Most of the light emitted by an ALS is not absorbed, but is reflected and scattered off the surface being examined. Extreme care should be taken around highly reflective surfaces. Never look directly into the light or allow beams to bounce off the surface into your eyes or the eyes of another person in the vicinity. Filtered goggles or shields shall be utilized when using this equipment as they provide protection from potentially harmful rays and provide additional enhancement for viewing impressions.

8.1.8.3 The nature and extent of all potential hazards are not yet known because in-depth assessments have not been made on most of the high intensity light sources used in forensic work.

8.1.9 REFERENCES:

Omnichrome Evidence Detection with Forensic Laser Technology, (1989).

Omniprint 1000A Operating Instructions, Omnichrome.

Mini-CrimeScope Tunable Forensic Light Source Model MCS-400W
Operation and Maintenance Instructions (2003).

UltraLite ALS Operator's Manual, Lynn Peavey Co., (2004).

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

9 PHYSICAL METHODS

9.1 IODINE FUMING

9.1.1 BACKGROUND:

Iodine fuming is one of the oldest impression methods currently employed in the examination processes for the visualization of impressions. Iodine vapors react with some residues in footwear or tire track impression deposits and turn the impression a yellow/brown color. Iodine fuming has an advantage over powders and chemical treatment because it does not cause the destruction of minute details in the impressions.

9.1.2 SCOPE:

9.1.2.1 Use when attempting to develop impressions that are composed of fatty, oily, or organic residues.

9.1.2.2 Other methods may cause interference. Therefore, if iodine fuming is to be used, it must be used prior to other development processes.

9.1.2.3 Iodine is not suitable for metals or dark surfaces.

9.1.3 EQUIPMENT AND MATERIALS:

Fume hood

Chamber or a heavy-duty sealable plastic bag

Glass container such as beaker or petri dish

9.1.4 REAGENT:

Iodine crystals

9.1.5 PROCEDURE 1 - CHAMBER METHOD:

9.1.5.1 In a fume hood, break open a glass ampoule of iodine crystals to reveal the iodine crystals.

9.1.5.2 Place the crystals in an airtight chamber (ex. sealable heavy plastic bag, commercial fuming chamber, etc.).

9.1.5.3 Apply heat if necessary. The application of heat may be accomplished in various ways including transfer of body heat, contained hot water, or an electric heater. Iodine crystals will start to sublime, (go from a solid to a gas), resulting in purplish fumes with the application of heat (approximately 100° F).

9.1.5.4 Place the control test and the questioned surface in the chamber and proceed with fuming.

9.1.5.5 The control test and evidence are monitored by viewing through the chamber to determine when processing is complete.

9.1.5.5.1 Impressions, if developed, will turn a yellow-brown color.

9.1.5.5.2 The process should be carefully monitored so that over-development does not occur.

9.1.5.6 Developed impressions are evaluated to determine their suitability for comparison.

9.1.5.7 Visible impressions are photographed or documented as soon as possible, and notes are taken.

9.1.6 ADDITIONAL INFORMATION:

9.1.6.1 The resulting yellow-brown impressions can vanish and must be preserved.

9.1.6.2 It is suggested that the camera be set up prior to iodine processing.

9.1.6.3 Iodine enhanced impressions that have faded, or are completely gone, can sometimes be redeveloped by processing with iodine fuming again. Iodine reprocessing cannot be done if other methods have been used or if too long of a time span has elapsed.

9.1.6.4 Shelf life of sealed iodine is indefinite.

9.1.6.5 Iodine crystals originating from glass ampoules shall be disposed of in the hazardous waste containers.

9.1.6.6 A method of fixing iodine impressions is to apply soluble starch spray or powder (shaking off the excess powder), and pass the exhibit through steam. This method should not be used if further processing is going to be done.

9.1.6.7 The iodine-enhanced impression may be further visualized by the application of 7,8-benzoflavone dip or spray.

9.1.7 CONTROLS:

9.1.7.1 Testing of iodine crystals is performed prior to each use.

9.1.7.2 This test involves the making of an oil-based mark on a test surface similar to the evidence being examined.

9.1.7.3 The test is exposed to the fumes in the same manner as the questioned surface would be.

9.1.7.4 When using the chamber method, testing of the iodine crystals and processing may be conducted at the same time.

9.1.7.5 The area surrounding the intentionally deposited mark shall serve as a negative control.

9.1.8 SAFETY:

9.1.8.1 Safety is a serious concern when using the iodine fuming method. *Iodine is toxic in any form. ALWAYS AVOID INHALING IODINE FUMES.*

9.1.8.2 Iodine fumes may irritate the skin and damage the respiratory tract. Headaches that can last for several days may result from exposure to iodine. Long-term effects to the thyroid gland may result from exposure.

9.1.8.3 Adequate ventilation when using the method is mandatory as the fumes are corrosive to metals and may discolor other surfaces that they come in contact with.

9.1.8.4 Iodine shall be purchased in glass ampoules. The ampoules shall stay sealed until use.

9.1.9 REFERENCES:

Cassidy, M.J., *Footwear Identification*, Lightning Powder Co., Inc., Salem, Oregon, 1980, page 53.

Keedwell, E., et al, "Chemical Methods for Enhancement of Footwear Marks", Metropolitan Police Forensic Services Lab, Report No. 73 (1988), page 18.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

9.2 LIFTING METHODS

9.2.1 BACKGROUND:

Lifting methods are effective for the preservation of some impressions because the adhesive on the lifting medium is stickier than the surface on which the impression deposit resides. It is a good idea to have a variety of lifting mediums as they vary in clarity, adhesion, and flexibility. The electrostatic dust lifting device (EDL) is effective on dry deposit impressions because an electrical charge induced in the particles of the deposit causes the impression to be attracted to and lifted onto an appropriate material. The electrostatic detection analyzer (ESDA) may be useful for dusty impressions on cardboard or paper because of a similar principle as the EDL. Indented impressions may be lifted from paper and cardboard by the ESDA because the cascade toner will adhere to slight indentations in the surface of the item causing visible black impression detail.

9.2.2 SCOPE:

9.2.2.1 Lifting methods are applicable to impressions that have first been developed utilizing other methods such as powders, SPR, and prints deposited in or composed of dust.

9.2.2.2 Some lifting methods may be effective on impressions made with blood, mud or other residues.

9.2.2.3 Lifts are inexpensive, easy, and a quick method of preserving developed impressions for future comparison.

9.2.2.5 Impression lifting is one method of preserving impressions at a crime scene.

9.2.2.5 Lifting may not be the most effective method of preserving a particular impression.

9.2.3 EQUIPMENT AND MATERIALS:

Exhaust vent or fume hood

Electrostatic Dust Lifting Device (EDL)

Electrostatic Detection Analyzer (ESDA)

Lifting film such as Mylar

Various sizes and types of standard lifting tapes

Adhesive lifters such as Handiprint

Gel lifters

Casting compounds

9.2.4 HANDIPRINT, TAPES, AND GEL LIFTERS:

9.2.4.1 Ensure that the surface has been prepared for lifting by removing excess powder, if powder was used first.

9.2.4.2 Lifting mediums should be removed from their backing in a smooth, continuous motion without hesitation to avoid lines in the adhesive.

- 9.2.4.3 The lifting medium is then applied to the impression bearing surface in a smooth continuous motion taking care to avoid air pockets and creases. It may be necessary to firmly rub the lifting medium onto the surface using a fair amount of pressure. If the gel lifter is applied to a residue impression, allow the lifter to remain on the impression for approximately 10 minutes to let the impression completely transfer.
- 9.2.4.4 Removal of the lifting medium from the impression bearing surface should also be performed in a smooth continuous motion and reapplied to a suitable backing material in the same manner as noted above.
- 9.2.4.5 Impression lifts shall be labeled with the following information:
Unique case identifier;
Date and initials.
Other useful information may include the impression source (description or source identifier) and significant information about the orientation and/or position of the impression on the object through description and/or diagram(s). One should be able to pinpoint the area and orientation of an impression on the object.
- 9.2.4.6 Lifts from multiple areas (different impressions) shall be placed on individual backing or paper.
- 9.2.5 CASTING COMPOUNDS (non-dental stone such as Mikrosil):
- 9.2.5.1 Ensure that the surface has been prepared for lifting by removing excess powder, if powder was used first.
- 9.2.5.2 Casting material is mixed either by hand or through the use of an extruder gun.
- 9.2.5.3 Casting material is applied to the impression-bearing surface in a manner that precludes air pockets. It may be necessary to place the casting material to the side of the impression and then smooth it across the surface.
- 9.2.5.4 The casting material is left in place until solidified.
- 9.2.5.5 The cast is then removed from the surface and attached to a suitable backing. The appropriate documentation is noted as detailed in 9.2.4.5.
- 9.2.6 CASTING – DENTAL STONE:
- 9.2.6.1 Dental stone is mixed by following the manufacturer's powder/water ratio supplied with the dental stone powder (usually 8 to 9 ounces of water for 2 pounds of dental stone). Sufficient material is mixed to pour one three-dimensional impression at a time.
- 9.2.6.2 The cast is allowed to completely harden before it is loosened and lifted.
- 9.2.6.3 The back of the cast is marked with the unique case identifier (if known), impression identifier (if applicable), initials and date.

- 9.2.6.4 The cast is loosely wrapped in paper to protect it and allowed to dry 24 to 48 hours before cleaning.
- 9.2.6.5 The cast may be cleaned in a cleaning solution of saturated potassium sulfate (K₂SO₄) for approximately 1 hour, rinsed and dried. Only dental stone casts should be cleaned in this manner.
- 9.2.6.6 Casts are positives and represent the footwear outsole or tire tread that made the impression.
- 9.2.6.7 Casts may be made at crime scenes or utilized as a test impression method for lab comparisons.

9.2.7 ELECTROSTATIC DUST LIFTS:

- 9.2.7.1 Place the item (dust print up) to be tested on the ground plate with the ground clip connected or set the antenna ground next to the area with the dust impression.
- 9.2.7.2 Cover the evidence area with a piece of lifting film (metal-coated side away from the possible impression).
- 9.2.7.3 With the Voltage Adjust Knob set at "LOW", touch the high voltage probe to the surface of the lifting film, turn on the power for 5 to 15 seconds or until the lifting film has been "pulled down" to the surface, and roll out any air bubbles with the roller. The voltage may have to be increased.
- 9.2.7.4 Turn off the power and let the charge dissipate for several seconds. Remove the lifting film from the surface and turn over away from the evidence area. View the lifted impression with oblique light.
- 9.2.7.5 To preserve the lifted impression, tape the piece of lifting film into a manila file folder or fine cardboard box, impression side up, and cover to protect it.
- 9.2.7.6 Because the dust lift is not permanently attached to the lifting film, it can be easily damaged or destroyed. The dust lift should be photographed and any analysis should be done using the photograph.
- 9.2.7.7 Dust lifts are positives (see section 9.2.6.6 above) and are the reverse of the impression as viewed on the original item or area. Photographic negatives should be reversed to produce a photograph that is the same as the original impression or tests impressions. An electrostatic dust lift can be compared directly to the shoe or tire, if extreme caution is taken to prevent damage to the lift.
- 9.2.7.8 For metal surfaces, the lifting film should be attached to a piece of heavy duty plastic film and the impression should be lifted onto the plastic film. The black mylar can be left attached to the back of the plastic film to make the impression easier to see.

9.2.8 ADDITIONAL INFORMATION:

- 9.2.8.1 Caution should be exercised in using general-purpose tapes (those not developed for lifting impressions) as they may cause migration

of some detail or may have striations or other imperfections making it hard to do comparisons.

9.2.8.2 Lifting should be performed after any necessary photography. The analyst's training and experience will determine the use and/or sequence of the lifting and photographic processes.

9.2.8.3 Store lifting mediums and casting compounds in a cool dry place.

9.2.8.4 Dispose of lifting mediums and casting compounds in the trash.

9.2.9 CONTROLS/SAFETY:

Controls are not applicable. There are no known health hazards associated with the use of lifting mediums, casting compounds or ESDA toner or Cascade Developer.

9.2.9.1 Persons using toner powders or Cascade Developer should monitor for allergic reactions (if any) to the powders.

9.2.9.2 Because the EDL uses high voltages, this device should not be used by persons who have pace makers and it should not be operated in wet conditions.

9.2.9.3 The corona unit on the ESDA must not be operated without the earth connection between the instrument body and the corona housing.

9.2.9.4 The EDL should not be used to lift dust impressions off living persons.

9.2.9.5 The lifting film or mylar should be isolated from metal conducting surfaces to prevent extreme sparking. Use of plastic barriers work well on metal surfaces.

9.2.10 REFERENCES:

Abbott, John R., *Footwear Evidence*, Charles C. Thomas, Springfield, 1964.

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990.

Cassidy, M.J., *Footwear Identification*, Lightening Powder Co., Inc., Salem, Oregon, 1995 (revised).

Electrostatic Dust Print Lifter Directions for Use, Kinderprint Co., Inc

ESDA a new instrument for the detection of indented writing in documents, Foster & Freeman Ltd., (no date).

ESDA – Applications of the Instrument for the Detection of Indented Writing In Documents, Foster & Freeman Ltd., April 1986.

ESDA Operating Instructions, Foster & Freeman, (no date).

Hamm, Ernest D., "Footprints in the Dust of Time", *Journal of Forensic Identification*, Vol. 43, No. 4, Jul/Aug 1993, pp357-361.

Hildebrand, Dwayne S., *Footwear – The Missed Evidence*, Staggs Publishing, Temecula, California, 1999.

Keedwell, E., Birkett, J., and Davis, R.J., "Chemical Methods for the Enhancement of Footwear Marks", MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

Majamaa, Heikki and Hamm, Ernest D., "Using the ESDA to Detect Dust Shoeprints on Paper", *Journal of Forensic Identification*, Vol. 41, No. 6, Nov/Dec 1991, pp 421-425.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

9.3 POWDER DETECTION METHODS

9.3.1 BACKGROUND:

Many commercially produced fingerprint powders are available and no powder is universally applicable to all types of non-porous surfaces. Most analysts stock a variety of different types and colors of powders as well as a variety of brushes for specialized applications. Powder particles physically adhere to impression residue allowing the impression to be visualized. This coloring of the residue occurs because the residue has greater adhesion properties than the substrate.

9.3.2 SCOPE:

9.3.2.1 Fingerprint powders are used to develop invisible detail, improve contrast of visible detail, and to facilitate lifting and preservation of impression evidence from non-porous surfaces.

9.3.2.2 The type of powder that is selected is dependent upon:

9.3.2.2.1 Whether resulting impression will be photographed. If so, a powder color that contrasts with the surface is often desirable.

9.3.2.2.2 The nature of the surface to be processed. Traditional powders are often most effective on non-textured surfaces while magnetic powders are often most effective on plastics and textured surfaces. The use of magnetic powders and wands should generally be avoided on substrates that contain iron. Fluorescent powders tend to have limited use. They are useful on multicolored surfaces or surfaces with a light texture that doesn't accept magnetic powder well.

9.3.2.3 The type of applicator selected is dependent upon:

9.3.2.3.1 The size of area to be dusted. Larger brushes are ordinarily used for large areas and smaller brushes on concentrated work or small partial impressions. Fiberglass brushes are often used for both instances.

9.3.2.3.2 The type of powder to be used. Magnetic wands are used in conjunction with magnetic powders while traditional powders are used with a variety of brushes. Traditional fluorescent powders are applied with a feather brush and their application requires the use of an ALS.

9.3.2.4 The prior use of cyanoacrylate ester often increases the adhesion of powders to impression residue.

9.3.2.5 Powder processing is not suitable for surfaces that are wet, tacky, or exceptionally rough and is generally the last step in the impression processing sequence.

9.3.3 EQUIPMENT AND MATERIALS:

Hood/exhaust vents/particulate filters
Traditional, magnetic, and fluorescent powders
Magnetic wand, feather brush, fiberglass brush, animal hair, etc.
Alternate light source
Filtered goggles

9.3.4 PROCEDURE 1 - TRADITIONAL POWDERS:

9.3.4.1 A variety of brushes or applicators may be utilized with the exception of magnetic wands.

9.3.4.2 Apply a small amount of powder to the brush and remove excess powder.

9.3.4.3 Powder should generally be applied to the surface in a smooth circular motion with only the tips of the brush touching the surface until optimal development is achieved

9.3.4.4 If too much powder has been applied, it may be possible to remove excess powder by tapping the object, blowing air over the surface, or by brushing it out.

9.3.4.5 Developed impressions are evaluated to determine their suitability for comparison.

9.3.4.6 Impressions deemed to be of value are marked and may be photographed or lifted.

9.3.5 PROCEDURE 2 - MAGNETIC POWDERS:

9.3.5.1 Magnetic powders generally utilize a magnetic wand in their application.

9.3.5.2 The wand is dipped into the magnetic powder where the powder is picked up by the tip of the wand. The powder actually forms a bristle-less brush that is then applied directly to the surface. The actual wand should not come in contact with the surface.

9.3.5.3 The application of magnetic powders is similar to the dusting method described in 9.3.4.3 & 9.3.4.4.

9.3.5.4 The plunger located at the end of the brush is pulled to its fully extended position to release the powder from the tip of the brush.

9.3.5.5 Excess powder may be removed by passing a wand over the surface without making contact.

9.3.5.6 Developed impressions are evaluated to determine their suitability for comparison.

9.3.5.7 Impressions deemed to be of value are marked and may be photographed or lifted.

9.3.6 PROCEDURE 3 - FLUORESCENT POWDERS:

9.3.6.1 A variety of brushes or applicators may be utilized.

9.3.6.2 Lightly dip the brush into the powder. Remove excess powder. A very small amount of fluorescent powder goes a long way.

9.3.6.3 If possible, it is best to use an ALS while applying the powder. This will prevent over powdering and loss of detail. The

application of fluorescent powders is similar to the dusting methods described in 9.3.4.3 & 9.3.4.4.

9.3.6.4 Developed impressions are evaluated to determine their suitability for comparison.

9.3.6.5 Impressions deemed to be of value are marked and may be photographed or lifted. When photographing impressions developed with fluorescent powders, it is necessary to use an ALS and a camera filter that corresponds to the color of viewing goggle utilized with the ALS. It is necessary to use black backing material with fluorescent powders.

9.3.7 ADDITIONAL INFORMATION:

9.3.7.1 Occasionally, impression quality may be enhanced by repeated powdering and lifting of the same area.

9.3.7.2 An ample number of appropriate brushes will help preclude cross-contamination of powders and brushes.

9.3.7.3 When powder-processing evidence that may be biologically contaminated, every effort shall be made to avoid cross contamination by utilizing previously unused brushes and powder. Brushes and powder will be discarded after use on contaminated items. Magnetic wands will be disinfected.

9.3.7.4 Powders stored in a cool dry place have an indefinite shelf life.

9.3.7.5 Dispose of powders in the trash.

9.3.8 CONTROLS:

Control impressions are not needed for powder processing. However, when there is doubt as to the suitability of a powder for processing a particular surface, the powder should be tested on a similar surface if available. If a similar surface is not available, then the powder may be tested on an edge of the suspected surface, carefully avoiding the area that may contain impressions.

9.3.9 SAFETY:

9.3.9.1 Safety concerns when using commercial fingerprint powders are minimal.

9.3.9.2 Analysts should use the hoods or exhaust vents when performing powdering and lifting in the laboratory.

9.3.9.3 When fingerprint powders are to be used for an extended period of time, a dust mask or half face respirator with dust filters should be worn to minimize the inhalation of the powder particles.

9.3.9.4 Persons using fingerprint powders should monitor for allergic reactions (if any) to the fingerprint powders.

9.3.10 REFERENCES:

Abbott, John R., *Footwear Evidence*, Charles C. Thomas, Springfield, 1964

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990.

Cassidy, M.J., *Footwear Identification*, Lightning Powder Co., Inc., Salem, Oregon, 1995 (revised)

Hildebrand, Dwayne S., *Footwear – The Missed Evidence*, Staggs Publishing, Temecula, California, 1999.

Keedwell, E., Birkett, J., and Davis, R.J., “Chemical Methods for the Enhancement of Footwear Marks”, MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

9.4 SMALL PARTICLE REAGENT

9.4.1 BACKGROUND:

Two types of small particle reagents (SPR) are available for use, traditional SPR that consists of a suspension of fine molybdenum disulfide (MoS₂) particles in a detergent solution and commercially available white SPR. These solutions work like a liquid fingerprint powder by adhering to the fatty or oily portion of the impression residue resulting in a gray or white colored impression.

9.4.2 SCOPE:

9.4.2.1 Small particle reagent is used to develop impressions from a variety of surfaces including adhesives and non-porous items that are or have been wet.

9.4.2.2 The color of SPR should be chosen to contrast with the background.

9.4.2.3 SPR may be used by dipping or spraying. Dipping is the preferred method as spraying is less sensitive. It is, however, difficult to prevent damage to impressions located on the lower side of an article in a dish and spraying is a valid alternative when processing large items, vehicles, or responding to crime scenes.

9.4.2.4 Surfaces that need other forensic examinations such as biology, questioned document, or trace examinations should be carefully evaluated prior to processing to determine if the SPR procedure will have an impact on subsequent examinations.

9.4.3 EQUIPMENT AND MATERIALS:

Beaker
Balance
Magnetic stirrer/stirring bar
Spray bottles
Processing tray

9.4.4 REAGENTS:

Commercially available white SPR
Molybdenum Disulfide (MoS₂)
Photo Flo 200
Distilled water

Small Particle Reagent Working Solution:

1. Place a 1500 ml beaker on magnetic stirrer base.
2. Add 1000 ml of distilled water to the beaker.
3. Place a magnetic stirring bar in the beaker.
4. Dissolve 30g of MoS₂ in the water (MoS₂ comes in 30g bottles).
5. Add three to four drops of Photo Flo 200 to the solution.

9.4.5 PROCEDURE 1 - DIPPING METHOD:

- 9.4.5.1 Shake or stir the SPR thoroughly and pour the solution into a tray.
- 9.4.5.2 Add the item to be processed to the solution. The item should be submerged.
- 9.4.5.3 Agitate the solution in the tray for 2-3 minutes, remove the item from the SPR and gently rinse with tap water.
- 9.4.5.4 Allow the surface to dry (if feasible).
- 9.4.5.5 Developed impressions are evaluated to determine their suitability for comparison.
- 9.4.5.6 Impressions deemed to be of value are marked and may be photographed or lifted. Depending on the circumstances, the item may or may not be dried prior to lifting.

9.4.6 PROCEDURE 2 - SPRAY METHOD:

- 9.4.6.1 Place the SPR into a spray bottle and shake thoroughly. The bottle should be shaken often to keep the MoS₂ in suspension.
- 9.4.6.2 Spray the SPR onto the item being examined. If the location of the impressions are known, spray the area above the impressions and allow the SPR to flow over them. Otherwise, spray the area to be examined starting at the top and working downwards.
- 9.4.6.3 Gently rinse the processed area with tap water and allow it to dry (if feasible).
- 9.4.6.4 Developed impressions are evaluated to determine their suitability for comparison.
- 9.4.6.5 Impressions deemed to be of value are marked and may be photographed or lifted. Depending on the circumstances, the item may or may not be dried prior to lifting.

9.4.7 ADDITIONAL INFORMATION:

- 9.4.7.1 Pre-mixed molybdenum has an indefinite shelf life. The shelf life the SPR working solutions is at least six months, but shall be tested prior to each use.

9.4.7.2 Excess reagent shall be collected and placed in the hazardous waste container located in the fume hood.

9.4.8 CONTROLS:

9.4.8.1 Testing of SPR is performed each day prior to use.

9.4.8.2 This test involves the making of oily mark on a test surface similar to the one being examined.

9.4.8.3 The test mark is exposed to the SPR in the same manner as the questioned surface.

9.4.8.4 An analyst shall not proceed with the processing of the evidence until a control test bearing positive results (development of a gray colored mark with traditional SPR or a white colored mark with white SPR) has been carried out and documented on the control tests work sheet.

9.4.8.5 The area surrounding the intentionally deposited mark shall serve as a negative control.

9.4.9 SAFETY:

There does not appear to be any health hazards associated with small particle reagent, but the process should be monitored to see if there are any allergies. Lab coats, gloves, and safety glasses should be worn.

9.4.10 REFERENCES:

Keedwell, E., Bickett, J., and Davis, R.J., "Chemical Methods for the Enhancement of Footwear Marks", MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

Lee, H.C. and Gaensslen, R.E., *Advances in Fingerprint Technology*, Elsevier, 1991, p 82.

Navarro, R.L., "Chemical Enhancement of Questioned Footwear Impressions", N.C. Bureau of Investigation, 1992, p 53.

9.5 CREATING TEST IMPRESSIONS

9.5.1 BACKGROUND:

Test impressions (reference or known exemplars) is a term used to describe impressions that are purposely made. These impressions may be made using a number of techniques, including, but not limited to, ink, photocopies, photography, three-dimensional impressions in sand or other material, casting, and powder/adhesive lift methods. The goal of the process is to produce legible impressions that are suitable for comparison.

9.5.2 SCOPE:

9.5.2.1 The footwear methods can be used on all footwear. If a shoe is heavily damaged, Biofoam or a casting method may be used to demonstrate what type of impression the shoe would leave.

9.5.2.2 Tires should remain mounted on the vehicle. Test impressions are made by pushing or driving the vehicle over the test impression medium.

9.5.2.3 The spare tire should not be overlooked when doing tire test impressions if there is an indication that it could have been on the vehicle at the time of the incident.

9.5.2.4 Additional test impressions may be necessary when the unknown impression shows distortion in the design elements from the way it was made.

9.5.3 EQUIPMENT AND MATERIALS:

Black printer's ink
Ink pad and roller
Kodak roller transport film
Black fingerprint powder
Inkless Shoeprint Kit
Vaseline or other grease such as WD40
Sand box
Dental stone
Fiber glass brush
Adhesive lifters
Butcher paper
Poster board (44"x28", may be cut in half lengthwise)
Biofoam
Silicone spray
Camera, 35mm or digital

9.5.4 REAGENTS:

None required

9.5.5 PROCEDURE - FOOTWEAR

9.5.5.1 Conduct a visual examination of the shoe for the presence of trace and biological evidence. If any probative evidence is present, documentation and collection must be completed before proceeding. Do test impressions of the shoes as received, then clean them (saving debris from each in separate packages) and repeat the test impression process.

9.5.5.2 Inkless Method

9.5.5.2.1 Press shoe to inkless pad.

9.5.5.2.2 Press shoe to piece of inkless paper.

9.5.5.2.3 Label test impression with "test impression", unique identifier, initials or signature, date, and which shoe is represented.

9.5.5.2.4 Repeat the method with the other shoe.

9.5.5.2.5 These inkless impressions can fade with exposure to light and a photocopy may be taken as well.

9.5.5.2.6 Transparent overlays can be made by using overhead projection transparency material and a copier.

9.5.5.3 Ink Method

9.5.5.3.1 Press shoe to inkpad.

9.5.5.3.2 Press inked shoe sole to plain white paper or damp roller transport film.

9.5.5.3.3 Label as in 9.5.5.2.3 above.

9.5.5.3.4 Repeat the method with the other shoe.

9.5.5.4 Powdered Test Impression Method

9.5.5.4.1 Cover the sole of the shoe with a heavy coating of black fingerprint powder and tap off the excess.

9.5.5.4.2 Press the shoe sole onto an adhesive lifter such as Handiprint.

9.5.5.4.3 Carefully position the transparent cover onto the test impression avoiding air bubbles.

9.5.5.4.4 Label as in 9.5.5.2.3 above.

9.5.5.4.5 Repeat the method with the other shoe.

9.5.5.4.6 By using transparent lifters, overlays can be made using this method. Otherwise, use overhead transparency material and photocopy as in 9.5.5.2.6.

9.5.5.5 Biofoam 3-D Method

9.5.5.5.1 Press the shoe into the Biofoam material.

9.5.5.5.2 Add a scale properly positioned next to the impression.

9.5.5.5.3 Photograph the impression using oblique lighting.

9.5.5.6 Sand 3-D Method

9.5.5.6.1 Press the shoe into an area of smoothed sand.

- 9.5.5.6.2 Add a scale properly positioned next to the impression.
- 9.5.5.6.3 Photograph the impression using oblique lighting.
- 9.5.5.6.4 The photograph may serve as documentation for the case file or a 1:1 enlargement of this photograph may be used for comparison.
- 9.5.5.6.5 For additional comparison to an unknown cast, the 3-D test impression may be cast with dental stone.


9.5.5.7 Other methods may be used such as potters clay or direct photography of outsole.

9.5.6 PROCEDURE - TIRES:

9.5.6.1 Ideally, tires are left in position on the vehicle. If mounted tires are received, they should be placed on a vehicle in the same position. Test impressions may be made of an un-mounted tire by rolling the tire along the test media. Conduct a visual examination of the tire for the presence of trace and biological evidence. If any probative evidence is present, documentation and collection must be completed before proceeding. Tire treads should be brushed clean of dirt, but all embedded stones should be left in place.

9.5.6.2 Inked Method

9.5.6.2.1 Outer sidewall of tire is divided into segments (usually four to seven) with tire crayon or chalk and each segment is numbered (1 up to 7 depending on number of segments)

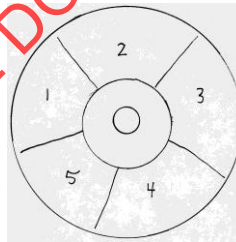

9.5.6.2.2 An inkpad is prepared by taping pieces of poster board together end to end until there is enough length for a complete tire rotation. Printer's ink is placed onto the inkpad and rolled into a thin, smooth layer from edge to edge. Position the completed inkpad in front of the selected tire.

9.5.6.2.3 The test impression board is prepared by taping pieces of poster board together end to end until there is enough length for a complete tire rotation plus 6" to 12" extra. The test impression made is made directly on the board or on plastic film taped to the board.

9.5.6.2.4 Position the test impression board in front of the inkpad. Roll the vehicle (push or drive) so the tire rotates over the inkpad. When the tire is coated with ink, continue the rotation of the tire along the test impression board.

9.5.6.2.5 As the tire moves along the test impression board, mark the segment bars and segment numbers on the test impression board beside the test impression.

9.5.6.2.6 Once the test impression is complete, the following information should be placed on the test impression:

- Unique identifier
- Initials or signature
- Date
- Which tire is represented
- Direction of rotation
- Inside and outside of impression

Other information that must be documented in notes:

- Complete vehicle information including, make, model, color, VIN, and license plate number.
- Complete tire information from sidewall of tire, including make, model, size, type of tire, etc.
- DOT number.
- Mold number, if possible.

9.5.6.2.7 Repeat the method for the other tires.

9.5.6.2.8 Allow the inked impressions to completely dry in a secure location. Package using appropriate evidence packaging and sealing techniques.

9.5.6.3 Grease Method

9.5.6.3.1 Follow the same instructions for marking the sidewalls as in 9.5.6.2.1.

9.5.6.3.2 Using proper safety precautions; jack up the tire until it "spins free". Coat the tire tread lightly with grease.

9.5.6.3.3 The test impression board is prepared by taping pieces of poster board together end to end until there is enough length for a complete tire rotation plus 6" to 12" extra. The test impression may be made directly on the board or on butcher paper (paper side up) taped to the board.

9.5.6.3.4 The board is placed in front of the tire, and the vehicle is rolled along the length of the board.

9.5.6.3.5 As the tire moves along the test impression board, mark the segment bars and segment numbers on the test impression board beside the test impression.

9.5.6.3.5 On the completed impression, add/collect the information as directed in 9.5.6.2.6.

9.5.6.3.6 Repeat the method for the other tires.

9.5.6.3.7 The "greased" impressions may be developed using black magnetic fingerprint powder.

9.5.6.3.8 Package using appropriate evidence packaging and sealing techniques.

9.5.6.4 The vehicle should be kept impounded in a secure area so that additional test impressions can be made and/or the tires can be directly compared to the impressions, if needed.

9.5.7 CONTROLS:
Not applicable

9.5.8 SAFETY:

All human biological material shall be treated as if infectious.

9.5.8.1 Gloves, eye protection, lab coat, and/or a protective disposable apron shall be worn at all times when working with footwear or tires contaminated with biological evidence.

9.5.8.2 Utensils shall be disposed of or cleaned and disinfected after use and surfaces will be disinfected.

9.5.8.3 Although there aren't any health effects associated with the black printer's ink used for tire test impressions, this ink is very viscous and difficult to remove from clothing and skin. The analyst should wear protective clothing and disposable gloves when working with this material.

9.5.8.4 Precautions should be taken when jacking up the vehicle.

9.5.9 REFERENCES:

Abbott, John R., *Footwear Evidence*, Charles C. Thomas, Springfield, 1964

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990.

Bodziak, William J., "Some Methods for Taking Two-Dimensional Comparison Standards of Tires", *Journal of Forensic Identification*, Vol. 46, No. 6, 1996, pp 689-701.

Cassidy, M.J., *Footwear Identification*, Lightning Powder Co., Inc., Salem, Oregon, 1995 (revised)

10 CHEMICAL METHODS

10.1 AMIDO BLACK BLOOD IMPRESSION PROCESSING

10.1.1 BACKGROUND:

Amido Black is also known as Amido Black 10B, Amido Black 12B, Naphthol Blue Black, or Naphthalene Black. Amido black is a dye that stains the protein portion of blood a blue-black color.

10.1.2 SCOPE:

10.1.2.1 Blood contaminated impressions may be processed with amido black to detect faint deposits of impression detail. It is generally used on dried blood stains on non-porous surfaces, but has been successful in developing impressions on some semi-porous and porous surfaces as well.

10.1.2.2 Amido black will not detect the normal constituents of impressions and therefore must be used in the proper sequence with other impression processing methods.

10.1.2.3 The amido black process utilizes a working solution, a rinse solution, and a wash solution (distilled water). Blood must be fixed prior to the application of amido black (unless using methanol in the amido black working solution as a fixative) to prevent the liquid solutions used in the process from washing away some or all of the blood deposits.

10.1.2.4 Bloodstains must be carefully examined and evaluated to preclude destruction of potentially valuable evidence. Any samples to be used for the biological examination of blood deposits or trace analysis should be collected prior to enhancement. It is often necessary to coordinate with investigators and/or other laboratory sections (biology for example) to determine which procedures may provide the most valuable findings.

10.1.3 EQUIPMENT AND MATERIALS:

Balance, magnetic stirrer/stirring bar
Pipettes
Beakers
Graduated cylinder
Storage bottles
Spray bottles

10.1.4 REAGENTS:

Amido Black
Glacial acetic acid
Methanol
Distilled water

5-sulfosalicylic acid

Amido Black Working Solution:

1. Weigh out 0.2g of amido black and place in a beaker.
2. Measure out 10mL of glacial acetic acid and add to the amido black.
3. Measure out 90mL of methanol and add to the beaker containing the amido black and glacial acetic acid.
4. Stir the solution with a magnetic stirrer for thirty minutes and transfer the solution to a clean spray bottle.

Rinse Solution 1 (de-stain):

1. Measure out 10mL of glacial acetic acid and pour into beaker.
2. Measure 90mL of methanol and add to the beaker.
3. Stir the solution for two to three minutes and transfer the solution to a clean, dry spray bottle.

Rinse Solution 2 (de-stain)

1. Measure out 5mL of glacial acetic acid and pour into a beaker.
2. Measure 98mL of methanol and add to the beaker.
3. Stir the solution for two to three minutes and pour into a spray bottle.

Fixative

1. Measure 20g of 5-sulfosalicylic acid and place into a 1500mL or larger beaker.
2. Add 1 liter of distilled water, mix well and pour into a storage bottle.

10.1.5 PROCEDURE:

10.1.5.1 Determine if samples for biology should be taken prior to processing

10.1.5.2 Conduct control tests.

10.1.5.3 Fix impressions using heat, methanol, 5-sulfosalicylic acid, or super-glue.

Blood can be fixed to an object by heating in a 100° centigrade oven for thirty minutes (restricted to non-heat sensitive objects).

Methanol may be sprayed or pipetted over the item. The item can be immersed in the sulfosalicylic acid solution for 10 minutes, then in distilled water for 5 minutes. Large areas such as floors may be sprayed with the sulfosalicylic acid solution.

The first amido black rinse that contains methanol may suffice for this "fixing" rinse. Super-glue is an effective method for non-porous evidence as it will fix all possible impressions not just those contaminated with blood.

10.1.5.4 Immerse the item in the amido black working solution for two to three minutes. Alternatively, the item may be sprayed or irrigated with the amido black working solution.

- 10.1.5.5 Immerse or irrigate the item with the de-stain rinse solution(s) to remove the excess dye.
- 10.1.5.6 Resulting impressions are a dark blue-black. The above process may be repeated to improve contrast.
- 10.1.5.7 Immerse or irrigate the surface with a distilled water wash (optional).
- 10.1.5.8 Allow the item to dry thoroughly.
- 10.1.5.9 Developed impressions are evaluated to determine their suitability for comparison.
- 10.1.5.10 Impressions deemed to be of value are marked and photographed.

10.1.6 ADDITIONAL INFORMATION:

- 10.1.6.1 Shelf life of the pre-mixed amido black working solution, de-stains and fixative is indefinite.
- 10.1.6.2 Excess reagent shall be collected, when possible, and placed in the hazardous waste container located in the fume hood.

10.1.7 CONTROLS:

- 10.1.7.1 Testing of amido black is performed each day prior to use.
- 10.1.7.2 Control tests are performed by the application of the reagent to a slide prepared with known blood. For safety reasons, analysts *will not* prepare impressions made with blood. A smear will be applied to the slide instead.
- 10.1.7.3 An analyst shall not proceed with the processing of the evidence until a control test bearing positive results (known blood staining a blue-black color) has been carried out and documented on the Reagent Preparation Log sheet or control tests work sheet.
- 10.1.7.4 The area surrounding the intentionally deposited blood smear shall serve as a negative control.

10.1.8 SAFETY:

- 10.1.8.1 Gloves, lab coats, goggles, and respirators, (if there is a chance of the reagents becoming airborne) are worn when mixing or using Amido Black.
- 10.1.8.2 Glacial acetic acid is corrosive and extremely irritating to the eyes and respiratory system. Avoid breathing the vapors and use in a fume hood, with a respirator, or with adequate ventilation. Glacial Acetic Acid will cause burns if it comes in contact with skin.
- 10.1.8.3 Methanol is *flammable*. It needs to be handled carefully and non-permeable gloves worn during the mixing and use of Amido Black. Methanol is toxic in quantities as small as 30 ml and should not be allowed to come in contact with the skin, eyes, or mouth. It is possible for methanol to be absorbed through the skin. If methanol comes into contact with the eyes or mouth, the

area should be flushed with generous amounts of water and a doctor may be consulted. Inhalation of methanol vapors should be kept at a minimum and the solution should be used in a well-ventilated area.

10.1.8.4 In addition, analysts must be aware of the biological hazards associated with blood and other body fluids and take extra precautions to protect themselves.

10.1.9 REFERENCES:

Bodziak, William J., *Footwear Impression Evidence*, CRC Press, 1999.

Journal of Forensic Identification, Vol. 45, No. 5 Sept/Oct 1995, "Superglue of Latent Shoe Prints in Blood Prior to Processing", pages 498-50.

Lee, H.C. and Gaensslen, R.E., *Advances in Fingerprint Technology*, Elsevier, 1991, p 85.

Navarro, R.L., "Chemical Enhancement of Questioned Footwear Impressions", North Carolina State Bureau of Investigations, 1992

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 1/14/2016

10.2 AMMONIUM THIOCYANATE

10.2.1 BACKGROUND:

Iron is found in many residues, particularly in soil or mud. Ammonium thiocyanate reacts with ferric iron in impression residues to form a red-brown color.

10.2.2 SCOPE:

- 10.2.2.1 This procedure works well for wet residues and muddy impressions on all surfaces.
- 10.2.2.2 Further enhancement with Physical Developer may be used after ammonium thiocyanate because nitric acid is used (other acids may interfere with Physical Developer).
- 10.2.2.3 Faint impressions sprayed with this reagent may be enhanced to show highly visible detail.
- 10.2.2.4 This reagent has better solubility than potassium thiocyanate so there isn't a need to filter after preparation of the reagent.

10.2.3 EQUIPMENT AND MATERIALS:

Balance, magnetic stirrer/stirring bar
Pipettes
Beakers
Graduated cylinder
Spray bottles

10.2.4 REAGENTS:

Ammonium thiocyanate
Acetone
2N nitric acid (1.2mL nitric acid/8.8mL distilled water)

Ammonium Thiocyanate Working Solution

1. Dissolve 2g of ammonium thiocyanate in 90mL of acetone.
2. Slowly add 10 mL of 2N nitric acid.
3. Transfer to fine mist spray bottle and loosely cap so gases can escape.

10.2.5 PROCEDURE:

- 10.2.5.1 Conduct control tests.
- 10.2.5.2 Apply two light spray applications of approximately 3 seconds each.
- 10.2.5.3 Resulting impressions are a red-brown or rust-brown color.
- 10.2.5.4 Developed impressions are evaluated to determine their suitability for comparison.
- 10.2.5.5 Impressions deemed to be of value are marked and photographed.

10.2.6 ADDITIONAL INFORMATION:

- 10.2.6.1 This reagent should be freshly prepared for each application.

- 10.2.6.2 The reaction product is water soluble so caution is required with non-absorbent surfaces.
- 10.2.6.3 For faint impressions on paper or cloth, a pre-wash treatment with 0.2% ferrous sulfate should be considered. The item is immersed in the pre-wash for 1 minute followed by 2 immersions in distilled water. The item is then placed on blotter paper and allowed to dry.
- 10.2.6.4 A source of contamination is excess soil from the surface.
- 10.2.6.5 If the substrate is red-brown or rust colored, gel lifting of the impression first and treatment of the gel lift with this reagent may be indicated.
- 10.2.6.6 Avoid over-spraying of the item so that the soil residue doesn't wash off the item.

10.2.7 CONTROLS:

- 10.2.7.1 Testing of ammonium thiocyanate is performed when the reagent is prepared.
- 10.2.7.2 Control tests are performed by the application of the reagent to a mark made with a ferric salt solution.
- 10.2.7.3 An analyst shall not proceed with the processing of the evidence until a control test bearing positive results (known mark staining a red-brown color) has been carried out and documented on the Impression Evidence Reagent Prep Log or control tests work sheet.
- 10.2.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.2.8 SAFETY:

- 10.2.8.1 Ammonium thiocyanate is an eye, skin and respiratory irritant and safety glasses should be worn when measuring it.
- 10.2.8.2 Acetone is highly flammable and the reagent should be prepared in the fume hood.
- 10.2.8.3 If the spray bottle cap is tightened, a build-up of gases will cause the bottle to blow the top and release the contents in a violent manner.

10.2.9 REFERENCES:

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990.

Davis, R.J., "Notes on the use of Chemical Methods for the Enhancement of Footwear Marks", Metropolitan Police Forensic Science Lab, January, 1988.

Davis, R.J., "A Systematic Approach to the Enhancement of Footwear Marks", Canadian Society of Forensic Science Journal, Vol. 21, No. 3, 1988, pp 98-105.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

10.3 7,8-BENZOFLAVONE

10.3.1 BACKGROUND:

Benzoflavone is absorbed by a range of fatty, oily and organic materials to give a blue colored impression after enhancement by iodine fuming.

10.3.2 SCOPE:

10.3.2.1 The blue color isn't permanent and will weaken with time. Photography should closely follow the treatment.

10.3.2.2 Post-treatment of an iodine-enhanced impression with benzoflavone will often visualize additional detail because of the darker blue color.

10.3.2.3 The item may be sprayed or dipped into the reagent.

10.3.2.4 Further chemical enhancement may be utilized after this treatment.

10.3.3 EQUIPMENT AND MATERIALS:

Balance, magnetic stirrer/stirring bar

Pipettes

Beakers

Graduated cylinder

Storage bottle

Spray bottle

10.3.4 REAGENTS:

7,8-benzoflavone

Chloroform

Petroleum ether

Benzoflavone Working Solution

1. Dissolve 0.2g of benzoflavone in 2 to 3mL of chloroform.
2. Make up to 100mL with petroleum ether.
3. Place in storage bottle or spray bottle.

10.3.5 PROCEDURE:

10.3.5.1 Conduct control tests.

10.3.5.2 Dip the item in the reagent for 10 seconds or spray.

10.3.5.4 Developed impressions are evaluated to determine their suitability for comparison.

10.3.5.5 Impressions deemed to be of value are marked and photographed.

10.3.6 ADDITIONAL INFORMATION:

10.3.6.1 Benzoflavone working solution has a shelf life of up to one year, but is generally made just before using.

10.3.6.2 This reagent is used as an additional enhancement step following iodine fuming.

10.3.7 CONTROLS:

10.3.7.1 Testing of benzoflavone is performed when the reagent is prepared and before each use if reagent has been stored.

10.3.7.2 Control tests are performed by the application of the reagent to a known mark in starch.

10.3.7.3 An analyst shall not proceed with the processing of the evidence until a control test bearing positive results (known mark staining a blue color) has been carried out and documented on the Impression Evidence Reagent Prep Log or control tests work sheet.

10.3.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.3.8 SAFETY:

10.3.8.1 Benzoflavone is an irritant and caution should be taken when weighing the powder.

10.3.8.2 Chloroform is a known mutagen and carcinogen. The reagent should only be used in a fume hood or with proper respiratory protection.

10.3.9 REFERENCES:

Bodziak, William J. *Footwear Impression Evidence*, Elsevier, New York, 1990.

Davis, R.J. "Notes on the use of Chemical Methods for the Enhancement of Footwear Marks", Metropolitan Police Forensic Science Lab, January, 1988.

Davis, R.J., "A Systematic Approach to the Enhancement of Footwear Marks" Canadian Society of Forensic Science Journal, Vol. 21, No. 3, 1988, pp 98-105.

Keedwell, E., Birkett, J., and Davis, R.J., "Chemical Methods for the Enhancement of Footwear Marks", MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

10.4 CYANOACRYLATE ESTER

10.4.1 BACKGROUND:

Cyanoacrylate ester (CAE) also referred to as "superglue," is sold as a number of brands and in a number of viscosities. Items that are to be processed with CAE need to be exposed to an atmosphere rich in CAE fumes. This may be accomplished through the use of a traditional fuming chamber, superglue-fuming wand, or vacuum chamber.

10.4.2 SCOPE:

10.4.2.1 Fuming with cyanoacrylate ester (super-glue) is a process that is used to visualize impression deposits on non-porous and some semi-porous objects. CAE processing also prepares the surface for the acceptance of powders and dye-stains that may enable further visualization of the impressions.

10.4.2.2 CAE is most effective on wet-origin impressions and least effective on dry-origin impressions.

10.4.2.3 When superglue vapors contact moisture and other components of impression residue the cyanoacrylate ester polymerizes fixing the impressions to the surface. This makes them more stable and less easily damaged.

10.4.2.4 The process is temperature, humidity, and pressure sensitive.

10.4.2.5 Objects that need additional forensic examinations such as trace or questioned document examinations should be carefully evaluated prior to processing to determine if this procedure will have an impact on subsequent examinations.

10.4.3 EQUIPMENT AND MATERIALS:

Fuming chamber (computerized)

Relatively airtight container such as a tank or sealed plastic bag

Vacuum chamber

Superglue fuming wand

Cups/warm water (optional)

Low temperature heating element (optional)

10.4.4 REAGENTS:

Cyanoacrylate gel or liquid

One shot fuming kit or equivalent

Superglue cartridges

10.4.5 PROCEDURE 1 - TRADITIONAL FUMING CHAMBER:

10.4.5.1 Select the appropriately sized fuming chamber.

10.4.5.2 Place the surface to be processed in the chamber (suspend if possible).

10.4.5.3 Add control test.

- 10.4.5.4 Add humidity to the chamber via cups of hot water (larger chambers will require more cups, smaller chambers fewer cups).
- 10.4.5.5 Allow the chamber to warm (if necessary) and humidity to build (80 degrees Fahrenheit and 80 % humidity is optimal but satisfactory results may be obtained at varying temperatures and humidity levels).
- 10.4.5.6 Add the CAE source.
- 10.4.5.6.1 Hot Plate Method - plug in the hot plate and place in the chamber. Add an approximately 2-3 cm in diameter pool of liquid superglue to a disposable aluminum dish and place on the hot plate.
- 10.4.5.6.2 Gel Packet Method - open and add one or more foil CAE gel packets (dependent on size of chamber, fuming rate, and analyst's preference) to the chamber. Once the gel is exposed to the air, the CAE will begin to vaporize at a controlled rate.
- 10.4.5.6.3 "ONE-SHOT" fuming kits - place the "activator solution" in the jar provided. Add the "activator canister" to the solution. Empty the CAE on to the top of the "activator canister." This method is generally reserved for crime scene response.
- 10.4.5.7 Secure the door to the chamber.
- 10.4.5.8 Fuming times will vary by the size of the chamber, the properties of the cyanoacrylate being used, the amount of heat and humidity, and the properties of the evidence being fumed. A control test should be carefully monitored by the analyst to prevent over or under fuming. Proper development is achieved when impression details on the control turn slightly white in color and begin to show good contrast. In the event of under fuming, the item may be re-fumed.
- 10.4.5.9 When development is complete evacuate the CAE fumes and remove the CAE source from the chamber.
- 10.4.5.10 Remove the item from the chamber and examine for comparable impression detail.
- 10.4.5.11 Impressions may be marked and photographed at this point, but are more commonly further enhanced with powders or dyes prior to preservation.

10.4.6 PROCEDURE 2 – SUPER GLUE FUMING WAND METHOD:

- 10.4.6.1 In a fume hood or other well ventilated area, place a superglue cartridge over the end of the fuming wand. Select cartridge size dependent upon amount and size of evidence.
- 10.4.6.2 Set control level to high and ignite the fuming wand. Fumes should be visible once the wand is hot, approximately 1-2 minutes.
- 10.4.6.3 Lower the heat level if desired.

- 10.4.6.4 Conduct a control test.
- 10.4.6.5 Fume the item by holding the fuming wand approximately 4-8 inches away. Fumes from the wand will rise so it is best to direct the fumes below your item if possible or deflect the fumes toward your item. Do not hold the wand too close or in the same area too long as damage and/or over development may occur.
- 10.4.6.6 Turn the fuming wand off and allow the unit to cool completely prior to removing cartridges or repackaging.
- 10.4.6.7 Examine item for impression detail.
- 10.4.6.8 Impressions may be marked and photographed at this point, but are more commonly further enhanced with powders or dyes prior to preservation.

10.4.7 PROCEDURE 3 - VACUUM CHAMBER METHOD:

- 10.4.7.1 Place items of evidence and controls into the vacuum chamber. It is not necessary to unfold garbage bags or leave large amounts of space between the items. *Do not place pressurize items such as sealed cans, bottles etc. in the chamber as they may explode.*
- 10.4.7.2 Add the CAE source. Foil CAE gel packs are recommended (number is dependent on chamber size and space), but a small dish with liquid CAE may also be used.
- 10.4.7.3 Place the lid on the vacuum chamber and close the release valve.
- 10.4.7.4 Turn on the vacuum pump.
- 10.4.7.5 Open the Gas Ballast Valve about one half turn.
- 10.4.7.6 Open the Isolation Valve (up position). If necessary, press on the lid until the chamber begins to evacuate.
- 10.4.7.7 Close the Gas Ballast Valve.
- 10.4.7.8 Evacuate the chamber to approximately 25 inches of mercury as shown on the chamber gauge.
- 10.4.7.9 Close the Isolation Valve.
- 10.4.7.10 Open the Gas Ballast Valve, wait 2-3 seconds and turn off the pump.
- 10.4.7.11 Close the Gas Ballast Valve.
- 10.4.7.12 Leave the items under vacuum for at least 20 minutes. There is no danger of over fuming.
- 10.4.7.13 Evacuate the chamber by slowly opening the release valve.
- 10.4.7.14 Remove glue and evidence. Examine item for impression detail.
- 10.3.7.15 Impressions may be marked and photographed at this point, but are more commonly further enhanced with powders or dyes prior to preservation.

10.4.8 ADDITIONAL INFORMATION:

- 10.4.8.1 In the event of over-fuming, it may be possible to use an adhesive lifting technique (tape, gel lifter etc.) to lift away heavy upper deposits, revealing underlying detail.

- 10.4.8.2 The “foil packets” may be stored at room temperature and have a shelf life of six months to a year. Liquid CAE and cartridges may be stored at room temperature with an indefinite shelf life.
- 10.4.8.3 CAE may be disposed of in the trash.
- 10.4.8.4 Analysts shall read the manufactures operating instructions for the super glue fuming wand and vacuum chambers prior to operating this equipment.

10.4.9 PROCEDURE 4 – CYANOACRYLATE FUMING CHAMBER:

- 10.4.9.1 Turn on power
- 10.4.9.2 The menu screens are designed to prompt the use of action to be taken to complete a full cycle. Screen is touch operated. Do not tap on screen with any object that could damage it.
- 10.4.9.3 Upon start-up the unit will load software and self calibrate.
- 10.4.9.4 Once running, the unit will prompt the user for each activity.
- 10.4.9.5 You will need to set the desired humidity level and fuming time. The unit default is 80% RH for 15 minutes. Refer to user’s manual for operating instructions

10.4.10 CONTROLS:

- 10.4.10.1 Testing of CAE and processing are performed at the same time.
- 10.4.10.2 A test impression is applied to a non-porous surface and put into the tank in an easily-monitored position with the questioned surface
- 10.4.10.3 When the development of the control test is complete, the questioned surface is also finished. Positive results are indicated by development of a white impression.
- 10.4.10.4 The area surrounding the intentionally deposited impression shall serve as a negative control.
- 10.4.10.5 Results of control tests shall be documented on control tests work sheet.

10.4.11 SAFETY:

- 10.4.11.1 Super glue fuming should only be conducted in well-ventilated areas. Precautions should be taken to avoid inhaling or allowing the vapors to contact the eyes, as the vapors can be irritating to the eyes, nose, and throat. Persons wearing contact lenses should not open CAE chambers without taking proper precautions. Non-vented goggles should be worn.
- 10.4.11.2 Precautions include using relatively sealed CAE chambers and evacuating the fumes from the chambers prior to removal of the questioned and test surfaces.
- 10.4.11.3 Gloves should be worn to prevent the cyanoacrylate from contacting the skin. If liquid glue is allowed to contact the skin, adhesion may result. If the skin sticks together, immerse

affected areas in warm water. This will loosen the skin so that it can be gently pulled apart.

10.4.12 REFERENCES:

Coleman Vacu-Print Instructions and Notes, Lightning Powder, (1995).

Cushman, B. Q. and Simmons, N.J., "A Cyanoacrylate Fuming Method for the Development of Footwear Impressions", *Journal of Forensic Identification*, Vol. 46, No. 4, Jul/Aug, 1996, p 412.

Lee, Henry C. and Gaensslen, R. E., *Advances in Fingerprint Technology*, Elsevier, 1991, p 67.

Llewellyn, P.E. and Scott-Dinkins, L., "A New Use for an Old Friend", *Journal of Forensic Identification*, Vol. 45, No. 5 July/August, 1996; Vol. 46, No. 1 Sep/Oct, 1995, p 498.

Paine, N., "Use of Cyanoacrylate Fuming and Related Enhancement Techniques to Develop Shoe Impressions on Various Surfaces", *Journal of Forensic Identification*, Vol. 48, No. 3, Sep/Oct, 1998, p 585.

Air Science, Operating Manual, SAFEFUME Cyanoacrylate Fuming Chamber, Rev 2 July-11-2008.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 1/14/2016

10.5 8-HYDROXYQUINOLINE

10.5.1 BACKGROUND:

8- hydroxyquinoline reacts with calcium, magnesium, iron, aluminum and other metal cations (positive ions) that may be present in small amounts in impression residues. The reaction causes fluorescence that is visible under ultra-violet (UV) light.

10.5.2 SCOPE:

10.5.2.1 The fluorescence is not permanent and will fade within a day or so. Photography with color or black and white film, an appropriate filter, UV light and long exposure times should follow within a short time.

10.5.2.2 This procedure works well on impressions made with either wet or dry residues.

10.5.2.3 This reagent is very useful on raw wood surfaces.

10.5.2.4 If the substrate fluoresces with this reagent but the residue doesn't contain any metal cations and the impression prevents fluorescing where it lies, a dark impression on a fluorescing background will be produced.

10.5.2.5 If the fluorescence fades, repeat spraying will not work.

10.5.2.6 Other chemical enhancement methods may be used after this treatment.

10.5.3 EQUIPMENT AND MATERIALS:

Balance, magnetic stirrer/stirring bar

Pipettes

Beakers

Graduated cylinder

Storage bottle

Spray bottle

10.5.4 REAGENTS

8-Hydroxyquinoline

Acetone

Distilled water

8-Hydroxyquinoline Working Solution

1. Mix 90mL of acetone with 10mL of distilled water.
2. Dissolve 0.5g of 8-hydroxyquinoline into the above mixture.
3. Stir well and pour into a storage bottle or spray bottle.

10.5.5 PROCEDURE

10.5.5.1 Conduct control tests.

10.5.5.2 Lightly spray the area of interest.

10.5.5.3 View with short wave and long wave UV light.

- 10.5.5.6 Developed impressions are evaluated to determine their suitability for comparison.
- 10.5.5.7 Impressions deemed to be of value are marked and may be photographed or lifted.

10.5.6 ADDITIONAL INFORMATION

- 10.5.6.1 8-Hydroxyquinoline working reagent has a shelf life of 6 months.
- 10.5.6.2 If both the substrate and residue have cations that fluoresce, this reagent will not work. Lifting the impression from the substrate by gel lifter first and spraying the gel lift may work.

10.5.7 CONTROLS

- 10.5.7.1 Testing of 8-hydroxyquinoline is performed each day prior to use.
- 10.5.7.2 This test involves the making of a mark in chalk on a test surface and following the processing procedure.
- 10.5.7.3 An analyst cannot proceed with the processing of the evidence until a control test bearing positive results (development of a purple fluorescing mark with UV light) has been carried out and documented on the Impression Evidence Reagent Prep Log or control tests work sheet.
- 10.5.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.5.8 SAFETY

- 10.5.8.1 This reagent is both a skin and respiratory irritant and should only be used in a fume hood or with proper respirator protective equipment. Disposable gloves should also be worn.
- 10.5.8.2 Acetone is very flammable and this reagent should not be used near an open flame or ignition source.
- 10.5.8.2 The reagent should be stored in an amber bottle in the fume hood.

10.5.9 REFERENCES

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990.

Davis, R.J., "Notes on the use of Chemical Methods for the Enhancement of Footwear Marks", Metropolitan Police Forensic Science Lab, January, 1988.

Davis, R.J., "A Systematic Approach to the Enhancement of Footwear Marks", Canadian Society of Forensic Science Journal, Vol. 21, No. 3, 1988, pp 98-105.

Keedwell, E., Birkett, J., and Davis, R.J., "Chemical Methods for the Enhancement of Footwear Marks", MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

10.6 LEUCOCRYSTAL VIOLET

10.6.1 BACKGROUND:

Leucocrystal Violet is a biological stain used to dye the blood hemoglobin components of impression residues an intense purple color. Due to the toxic nature of this reagent, it should only be used in small quantities with the appropriate safety precautions observed.

10.6.2 SCOPE:

10.6.2.1 Leucocrystal violet is a dye stain used to visualize impression deposits on many types of adhesive surfaces.

10.6.2.2 Leucocrystal violet may also be used on small non-porous surfaces contaminated with grease and oils. It is not suitable for water-soluble adhesives or porous surfaces.

10.6.2.3 Surfaces that need other forensic examinations such as biology or trace should be carefully evaluated prior to processing to determine if this procedure will have an impact on subsequent examinations.

10.6.3 EQUIPMENT AND MATERIALS:

Balance
Magnetic stirrer/stirring bar
Graduated cylinder
Glass beaker
Glass tray
Storage bottles

10.6.4 REAGENTS:

Leucocrystal Violet
5-sulfosalicylic acid
3% hydrogen peroxide
Distilled water

10.6.4.1 Formula "A"

10.6.4.1.1 Dissolve 10g of 5-sulfosalicylic acid in 100ml distilled water.

10.6.4.1.2 Add 400ml 3% hydrogen peroxide to sulfasalicylic acid solution.

10.6.4.1.3 Immediately prior to use, add 0.75g leucocrystal violet to above. Stir the mixture vigorously.

10.6.4.2 Formula "B"

10.6.4.2.1 10g 5-sulfosalicylic acid dissolved in 500ml 3% hydrogen peroxide.

10.6.4.2.2 Add 3.7g sodium acetate and 1.0g leucocrystal violet. Stir the mixture vigorously.

10.6.5 PROCEDURE:

- 10.6.5.1 Determine if samples for biology should be taken prior to processing.
- 10.6.5.2 Conduct control tests.
- 10.6.5.3 Spray the impression using a fine mist sprayer. Items may also be soaked or the surface flooded with the solution.
- 10.6.5.4 Development of dark purple impressions should occur in 30 seconds.
- 10.6.5.5 Developed impressions are evaluated to determine their suitability for comparison.
- 10.6.5.6 Impressions deemed to be of value are marked and may be photographed or lifted.

10.6.6 ADDITIONAL INFORMATION:

- 10.6.6.1 Shelf life of the working solution is approximately 3 months.
- 10.6.6.2 Excess reagent shall be collected and placed in the hazardous waste container located in the fume hood.

10.6.7 CONTROLS:

- 10.6.7.1 Testing of leucocrystal violet is performed each day prior to use.
- 10.6.7.2 This test involves the making of a mark in blood on a slide and following the processing procedure.
- 10.6.7.3 An analyst cannot proceed with the processing of the evidence until a control test bearing positive results (development of a purple mark) has been carried out and documented on the Reagent Preparation Log sheet or control tests work sheet.
- 10.6.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.6.8 SAFETY:

- 10.6.8.1 Leucocrystal may be harmful by inhalation, ingestion, or skin absorption. May cause skin and eye irritation. May cause irritation to mucous membranes and upper respiratory tract.
- 10.6.8.2 Leucocrystal violet should not be used in large amounts.
- 10.6.8.3 A respirator should be used when working with the dry form. Leucocrystal violet should be prepared and used in a fume hood or well-ventilated area. The analyst should wear a lab coat, heavy-duty (non-disposable) gloves, and safety glasses.
- 10.6.8.4 In addition, analysts must be aware of the biological hazards associated with blood and other body fluids and take extra precautions to protect themselves.

10.6.9 REFERENCES:

Bodziak, William J., "Use of Leucocrystal Violet to Enhance Shoe Prints in Blood", *Forensic Science International*, Vol. 82, No. 1, September 1996.

Chemical Formulas and Processing Guide for Developing Latent Prints, US Department of Justice, 1994, pp 47-48.

Fisher, John F., "An Aqueous Leucocrystal Violet Enhancing Reagent for Blood Impressions", *Symposium on the Forensic Aspects of Footwear and Tire Impression Evidence*, FBI Academy, 1994.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

10.7 LEUCOMALACHITE GREEN (LCG)

10.7.1 BACKGROUND:

Leucomalachite green, is a biological stain used to dye the blood hemoglobin components of impression residues an intense green color. Due to the toxic nature of this reagent, it should only be used in small quantities with the appropriate safety precautions observed.

10.7.2 SCOPE:

- 10.7.2.1 Leucomalachite green is a dye stain used to visualize impression deposits on many types of surfaces without coloring the background. There is no need for destaining.
- 10.7.2.2 Leucomalachite green may also be used on small non-porous surfaces contaminated with grease and oils.
- 10.7.2.3 Surfaces that need other forensic examinations such as biology or trace should be carefully evaluated prior to processing to determine if this procedure will have an impact on subsequent examinations.

10.7.3 EQUIPMENT AND MATERIALS:

Balance
Magnetic stirrer/stirring bar
Graduated cylinder
Glass beaker
Glass tray
Spray bottle

10.7.4 REAGENTS:

Leucomalachite green
Sodium perborate
Methanol
Glacial acetic acid
1,1,1-trichloroethane)

Leucomalachite Green Working Solution

1. Mix together 0.06g of leucomalachite green, 0.2g sodium perborate, and 20mL methanol and 10mL glacial acetic acid.
2. Stopper the flask tightly and wrap the stopper with Parafilm.
3. Shake the flask vigorously for 1 minute.
4. Transfer to a spray bottle.

10.7.5 PROCEDURE:

- 10.7.5.1 Biological samples should be taken prior to processing.
- 10.7.5.2 Conduct control tests.
- 10.7.5.3 Spray the surface lightly 2 or 3 times, holding the sprayer approximately 14" or more from the surface.

- 10.7.5.4 Development of dark green impressions should occur rapidly.
- 10.7.5.5 Developed impressions are evaluated to determine their suitability for comparison.
- 10.7.5.6 Impressions deemed to be of value are marked and may be photographed or lifted.

10.7.6 ADDITIONAL INFORMATION:

- 10.7.6.1 Leucomalachite green reagent should be prepared with each application.
- 10.7.6.2 This reagent does not color the substrate and is very good for absorbent paper, cloth, carpet or other porous substrates.

10.7.7 CONTROLS

- 10.7.7.1 Testing of leucomalachite green is performed prior to use.
- 10.7.7.2 This test involves the making of a mark in blood on a slide and following the processing procedure.
- 10.7.7.3 An analyst cannot proceed with the processing of the evidence until a control test bearing positive results (development of a dark green mark) has been carried out and documented on the Reagent Preparation Log sheet or control tests work sheet.
- 10.7.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.7.8 SAFETY:

- 10.7.8.1 It can be harmful if inhaled, and is irritating to the eyes and skin.
- 10.7.8.2 Leucomalachite green should not be used in large amounts.
- 10.7.8.3 A respirator should be used when working with the dry form outside of the fume hood. Leucomalachite green should be prepared and used in a fume hood or well-ventilated area. The analyst should wear a lab coat, heavy-duty (non-disposable) gloves, and safety glasses.
- 10.7.8.4 In addition, analysts must be aware of the biological hazards associated with blood and other body fluids and take extra precautions to protect themselves.

10.7.9 REFERENCES:

Moses, Kenneth R., "Special Techniques for Latent print Development", San Francisco Police Department instructions.

Patty, John R. and Giberson, Michael W., "Bloody Fingerprint Spray Reagent", CAC Newsletter, 1984.

Stow, K.M., "Direct lift-enhancement of blood-contaminated shoe marks by leuco malachite green-impregnated membranes", Journal of the Forensic Science Society, 34(4), 1994, pp 241-244.

10.8 PHYSICAL DEVELOPER (PD)

10.8.1 BACKGROUND:

Physical developer is a silver-based aqueous reagent that reacts with lipids, fats, oils, and waxes present in the impression residue to form a silver-gray deposit.

10.8.2 SCOPE:

10.8.2.1 Physical developer is a method used for the development of impressions on porous substrates. It is not suitable for non-porous surfaces.

10.8.2.2 This method is the final step in the sequential processing of porous items.

10.8.2.3 Physical developer is the only method to show adequate results on paper that has been wet.

10.8.2.4 Surfaces that need other forensic examinations such as body fluid, trace, or questioned document examinations should be carefully evaluated prior to processing to determine if this procedure will have an impact on subsequent examinations.

10.8.3 EQUIPMENT AND MATERIALS:

Graduated cylinder

Glass trays

Plastic tongs

10.8.4 REAGENTS:

Physical Developer Kit (parts A & B)

1. Any contamination may ruin the physical developer working solution. To avoid contamination use clean glassware rinsed with tap water, then with distilled water prior to beginning.
2. Add 5 ml of solution A (20% silver nitrate solution) to 90 ml of solution B (reductant solution) in a beaker.
3. Stir the working solution for approximately one minute with a clean glass/plastic stirring rod.
4. Do not mix the working solution until you are ready to use it as it does not have a very long shelf life once mixed.

10.8.5 PROCEDURE:

10.8.5.1 Arrange glass trays in a sink so that the evidence can be moved easily from one tray to another in the proper sequence.

10.8.5.2 Add the physical developer working solution to its dedicated glass tray.

10.8.5.3 Use plastic photographic tongs or plastic forceps without serrated edges to add or remove articles from PD solutions. Do not use metal tools.

10.8.5.4 Conduct control tests.

- 10.8.5.5 Immerse the item and gently rock the tray for approximately 5-15 minutes until impression development is complete or adequate time has elapsed (analyst's discretion).
- 10.8.5.6 Remove the item from the physical developer working solution and place into a tray with running tap water. Rinse until the water runs clear.
- 10.8.5.7 Dry completely.
- 10.8.5.8 Developed impressions are evaluated to determine their suitability for comparison.
- 10.8.5.9 Impressions deemed to be of value are marked and photographed.

10.8.6 ADDITIONAL INFORMATION:

- 10.8.6.1 Cleanliness is important in the physical developer method. A good deal of the instability in the earlier solutions was a result of laboratory equipment that was not spotless. Some contaminants, especially salts, will cause the silver nitrate in the solution to come out of suspension thus spoiling the physical developer solution and perhaps ruining the item being examined. It is important to keep the glassware spotless and rinsed with distilled or de-ionized water prior to use. When washing glassware, use detergent, not abrasive cleaners.
- 10.8.6.2 Physical developer will cause dark stains on many surfaces. Care must be taken to avoid spills in the laboratory. Full strength chlorine bleach will usually remove any stains from counter tops and floors, but the bleach may cause damage to fabrics stained with physical developer.
- 10.8.6.3 Shelf life for ready to use kit (un-mixed) is six months from date of purchase. The reagent shall be mixed upon each use.
- 10.8.6.4 Excess reagent shall be collected and placed in the hazardous waste container located in the fume hood.

10.8.7 CONTROLS:

- 10.8.7.1 Testing of physical developer is performed prior to each use.
- 10.8.7.2 This test involves the making of a quality (oil based) mark on a test surface similar to the evidence being examined and following the processing procedure.
- 10.8.7.3 An analyst shall not proceed with the processing of the evidence until a control test bearing positive results (development of a silver-gray mark) has been carried out and documented on the control tests work sheet.
- 10.8.7.4 The area surrounding the intentionally deposited mark shall serve as a negative control.

10.8.8 SAFETY:

10.8.8.1 Physical developer should only be used in well-ventilated areas, as it is irritating to the respiratory tract. Standard laboratory protocol is followed for chemical handling.

10.8.9 REFERENCES:

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990, pp 144-148.

“Chemical formulas and Processing Guide for Developing Latent Prints”, U.S. Department of Justice, FBI Laboratory, 1994, p 35.

Keedwell, E., Birkett, J., and Davis, R.J., “Chemical Methods for the Enhancement of Footwear Marks”, MPFSL Report No. 73, Metropolitan Police Forensic Science Laboratory, January, 1988.

Navarro, R.L., “Chemical Enhancement of Questioned Footwear Impressions”, North Carolina State Bureau of Investigation, 1992.

“Physical Developer”, Manual of Fingerprint Development Techniques, British Home Office Scientific Research & Development Branch, 1986.

Technical Note #1-2730, Lightning Powder Co., (1993).

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy 1/14/2016
OBSOLETE DOCUMENT

11 DIGITAL IMAGING PROCEDURE

11.1 BACKGROUND:

Shoe and tire track images are frequently captured, enhanced, and stored using digital devices. The intent of image enhancement is to make details of an image that are less visible more visible. Enhancement may be used to increase the contrast between the print and the substrate, reverse the color of the ridges, etc.

11.2 SCOPE:

This sets forth the Impression Evidence Section's procedures for the capture, storage, enhancement, and output of digital images.

11.3 RESPONSIBILITIES:

11.3.1 Latent/Impression Evidence Section Supervisor

11.3.1.1 The Latent/Impression Evidence Section Supervisor shall act as the Digital Imaging System Administrator or appoint a Digital Imaging System Administrator.

11.3.1.2 The Latent/Impression Evidence Section Supervisor shall oversee and document the training of each new digital imaging system operator. This includes documenting competency testing.

11.3.1.3 The Latent/Impression Evidence Section Supervisor shall ensure access is limited to authorized users.

11.3.1.4 The Latent/Impression Evidence Section Supervisor or designee shall act as a liaison with CJIS and Foray technical staff on system maintenance, upgrades, and when technical difficulties arise.

11.3.1.5 The Latent/Impression Evidence Section Supervisor or designee shall be the only personnel authorized to delete images or cases entered into Digital Workplace or equivalent software.

11.3.2 Digital Imaging System Administrator

11.3.2.1 The Digital Imaging System Administrator shall update the Shoe/Tire track Section Digital Imaging System User's Manual.

11.3.2.2 The Digital Imaging System Administrator shall be responsible for system maintenance to include: tape back-ups, deletion of images/cases, archiving, etc.

11.3.2.3 The Digital Imaging System Administrator shall communicate system status to the supervisor and other system users.

11.3.3 Analysts

11.3.3.1 Analysts shall only use enhancement techniques that are supported by their training and/or experience.

11.3.3.2 Analysts shall maintain system security.

11.3.3.2.1 Network and/or program passwords are not to be distributed to unauthorized users. Operators may change their passwords as needed.

11.3.3.2.2 The external modem shall remain off unless technical support is being contacted.

11.3.3.3 Analysts shall fill out the CD/DVD Log when filing or retrieving archived images from the vault.

11.4 DIGITAL IMAGE CAPTURE

11.4.1 All digital evidentiary images shall be acquired through Digital Workplace or equivalent software.

11.4.2 Digital Workplace or equivalent software shall establish a chain of custody from the time of acquisition into the program.

11.4.3 Images shall be designated using a file name structure generated by Digital Workplace or equivalent software.

11.4.4 Analysts shall use one of the following digital image capture devices to acquire images of the print(s) in question.

11.4.4.1 Flat Bed Scanner

11.4.4.2 Digital Camera

11.4.4.3 Digital Media (e.g. Thumb Drive, CD/DVD, etc.)

11.4.4.4 Outside agencies may submit processed film for digital capture or digitally submit impression images.

11.4.4.4.1 Images of shoe and tire prints should contain a scale.

11.4.4.4.2 It is preferred that existing images be submitted in a loss-less format such as '.tif' and at as high a resolution as possible.

11.4.5 All original close up images captured by impression evidence section analysts shall contain a scale in centimeters.

11.4.6 Images shall be captured at the highest resolution practical for their intended use.

11.5 DIGITAL IMAGE ENHANCEMENT

11.5.1 All digital evidentiary images requiring enhancement shall be enhanced via Adobe Photoshop (using a copy of the original image) through Digital Workplace or equivalent software.

11.5.2 Enhanced images will be designated using a file name structure generated by Digital Workplace or equivalent software.

11.5.3 All enhancement history shall be recorded via Digital Workplace or equivalent software.

11.5.4 All images stored in a secured digital imaging system maintained by ISP Forensic Services shall serve as documentation for the case record.

11.6 DIGITAL IMAGE STORAGE, ARCHIVAL, AND RETRIEVAL

11.6.1 All images, both original and enhanced, shall be temporarily stored on the digital imaging system hard drive until the examination is completed.

11.6.2 Once completed, the case's originating analyst shall make an entry in the Digital Imaging Notebook on the "Cases to be Archived" form. The entry

- shall contain the date entered, analyst's initials, and complete case number to be archived.
- 11.6.3 A backup shall be completed by the LPS Supervisor or Digital Imaging System Administrator no less than once a week, more often if needed.
 - 11.6.4 Archiving of images shall be completed by the LPS Supervisor or Digital Imaging System Administrator on an as needed basis.
 - 11.6.4.1 It is recommended that images be recorded on Write-once Compact Disk DVD+R.
 - 11.6.4.2 The DVD along with a printed copy of listed contents shall be stored in the evidence vault. A printed copy of listed contents shall also be placed in the Archived Cases Log Book.
 - 11.6.5 All DVDs shall be logged in and out of the vault using the "Section DVD Log" sheet. The log shall detail the DVD title, date out/date returned, requesting analyst, and the person checking it in or out.
 - 11.6.6 Cases with no statute of limitation shall be stored on the digital imaging server hard drive indefinitely.
- 11.7 QUALITY CONTROL:
- 11.7.1 Performance checks shall be conducted on equipment as needed.
 - 11.7.2 When a problem is noted with a particular piece of equipment, software program, etc., the Digital Imaging System Administrator and/or Latent Section Supervisor shall be notified.
 - 11.7.3 If it is determined that the situation is persistent or cannot be easily rectified, an entry shall be made on the "Instrument Maintenance Log".
 - 11.7.3.1 The log shall detail the date, the person making the entry, the piece of equipment/software involved, and relevant details of the situation.
 - 11.7.4 Effected equipment/software shall be taken off line and all users notified.
 - 11.7.5 If necessary, technical support shall be sought and/or the equipment repaired before being put back into operation.
 - 11.7.6 Actions taken to repair or correct the problem shall be documented on the "Instrument Maintenance Log."
 - 11.7.7 Image calibration shall be checked, as needed by comparing the scale in the printed image with a standard metric scale.
- 11.8 TRAINING
- 11.8.1 All analysts utilizing imaging technologies shall be trained and tested for competency in the standard operating procedures and the operation of the relevant imaging technologies.
 - 11.8.2 Formal training may be modified at the discretion of the Latent Section Supervisor dependent upon previous training and/or experience.
 - 11.8.3 Recommended formal training consists of:
 - 11.8.3.1 Reviewing the ISP-FS Impression Evidence Section Digital Imaging Procedure.

- 11.8.3.2 Reviewing the Digital Workplace Quick Reference Guide or equivalent.
- 11.8.3.3 Reviewing the ISP Impression Evidence Section Digital Imaging User's Manual.
- 11.8.3.4 Review of relevant chapters of the Adobe Photoshop Users Manual and/or completion of a digital imaging course that utilizes Adobe Photoshop.
- 11.8.3.5 Satisfactory creation and digital processing of a mock-case using Digital Workplace and Adobe Photoshop software or equivalent software.
- 11.8.3.6 Satisfactory completion of a written test.
- 11.8.4 Continuing education shall be provided as courses become available through outside sources such as Foray, the FBI, etc.
- 11.8.5 Competency testing shall be repeated when significant changes in hardware or software are made.

11.9 REFERENCES:

International Association for Identification "Resolution 97-9."

Digital Workplace Quick Reference Guide.

Scientific Working Group on Imaging Technologies (SWGIT), "Best Practices for Documenting Image Enhancement" Version 1.2, March 04, 2004.

Scientific Working Group on Imaging Technologies (SWGIT), "Definitions and Guidelines for the use of Digital Image Technologies in the Criminal Justice System," Version 2.3-June 6, 2002.

Scientific Working Group on Imaging Technologies (SWGIT), "Recommendations and Guidelines for the Use of Digital Image Processing in the Criminal Justice System," Version 1.2-February 2001.

Scientific Working Group on Imaging Technologies (SWGIT), "Guidelines and Recommendations for Training in Imaging Technologies in the Criminal Justice System," Version 1.2-December 6, 2001. *Forensic Science Communications*, April 2002-Volume 4-Number 2.

Scientific Working Group on Imaging Technologies (SWGIT), "Guidelines for Field Applications of Imaging Technologies in the Criminal Justice System," Version 2.3, December 6, 2001. *Forensic Science Communications*, April 2002-Volume 4-Number 2.

12 FOOTWEAR AND TIRE IMPRESSION EXAMINATION AND COMPARISON PROCESS

12.1 BACKGROUND:

12.1.1 When shoes and tires impact with a surface, a record of that impact remains on the surface.

12.1.2 The amount of detail in the impression depends on many factors (not all-inclusive):

- Type of substrate or surface.
- Color, pattern and texture of the surface.
- Type of residue in the impression.
- Type of footwear or vehicle tire.
- Resolution of impression photographs.
- Type of casting material used.
- Tire pressure, vehicle load, wheel alignment, wheel balance, ambient temperature, and other factors.
- How the footwear impression was made, such as shoe impact on door or step on flooring
- Physical characteristics of shoe wearer such as weight, type of stride or infirmities.

12.1.3 Footwear and tires may leave unique impressions depending on the level of detail and the presence of one or more individual characteristics.

12.1.4 An impression contains class characteristics.

12.1.5 An impression containing a sufficient quantity and quality of detail may be individualized to or excluded from a particular source.

12.1.6 No scientific basis exists for requiring a pre-determined minimum number of individual characteristics to be present in two impressions in order to establish a positive identification.

12.1.7 Identification/Exclusion is supported by the theories of the random nature of individual characteristics, probability modeling, and empirical data gained through more than one hundred years of operational experience.

12.2 SCOPE:

After footwear and tire impressions are documented and collected in the field, analysts will use a combination of visual inspection, documentation, photography, test impression preparation, and other procedures to compare the known exemplar to the questioned impression to reach a decision about the source of the questioned impression. Impressions from one crime scene or area of a crime scene may be compared to impressions in another crime scene or area of the same crime scene to determine if they have a common origin.

12.3 EQUIPMENT AND MATERIALS

Magnifiers

Dividers

Finely divided rulers and scales

Stereomicroscope

Digital and/or 35mm cameras
Digital imaging system

12.4 PROCEDURE:

12.4.1 The case report is examined to determine if other testing is required on any items of evidence prior to the impression comparison process. This determination may also be made by discussion with other analysts and/or the investigating officer.

12.4.2 The impression evidence is documented in the notes and is properly labeled with the unique case identifier and initials of the analyst if possible. Photographic documentation may also be done.

12.4.3 If casts are included, they are cleaned, either by hand or in the cleaning tank if necessary.

12.4.4 The footwear or tire impression is assessed to determine suitability for comparison.

12.4.4.1 The analyst should determine if the impression contains sufficient detail to proceed with the comparison.

12.4.4.2 If “no”, the notes should reflect that the impression was not suitable for comparison.

12.4.4.3 If “yes”, the analysis continues.

12.4.4.4 At any stage of the process, exclusion is reached when the impressions being compared are in disagreement or contain an unexplainable discrepancy.

12.4.4.5 Inconclusive findings result from the absence of sufficient details (individual characteristics) to effect a conclusion of individualization or exclusion, although class characteristics may be comparable.

12.4.5 Comparison Process For Footwear and Tires

12.4.5.1 If the design is different, document, discontinue this process and report accordingly. This is an exclusion.

12.4.5.1 If the design is similar and there is sufficient detail, continue with the process.

12.4.5.2 Size and Shape: The analyst compares the physical size and shape of the impression design with the test impressions and/or shoes and tires.

12.4.5.2.1 Tire Noise Treatment Analysis

12.4.5.2.1.1 Noise treatment analysis can help determine if the suspect tire is consistent in size and design with the crime scene impression, if the impression has sufficient size and detail for this process.

12.4.5.2.1.2 Noise treatment analysis can define the location(s) on the circumference of the suspect tire that could have made the crime scene impression and once the location(s) is determined, the

examination for general wear and individual characteristics can proceed.

12.4.5.2.1.3 The process utilizes a high quality crime scene impression photograph (or cast) and a full circumference transparent tire test impression, or a crime scene impression transparency and a full circumference tire test impression on poster board or paper. In circumstances where a transparency can't be made, noise treatment can be determined by comparing the test impression to the manufacturer's mold drawings.

12.4.5.2.1.4 Any wear bars that are visible in the impression are marked. Wear bars are marked on the tire test impression, bases on the markings on the sidewall and the segment mark and segment number on the test.

12.4.5.2.1.5 The transparency is moved along the item until the noise treatments coincide. This area is marked.

12.4.5.2.1.5 Once the location(s) has been found, the comparison can proceed. If noise treatment doesn't coincide and a cast is available, repeat using the cast. If only a photograph is available, consider scaling, perspective and other issues that may account for the differences.

12.4.5.2.2 If the physical size and shape are different, document, and consider scaling, perspective and other issues that may account for differences. If the differences can't be explained, discontinue this process and report accordingly. This is an exclusion.

12.4.5.2.3 If the physical size and shape, to include noise treatment of tires, are similar and there is sufficient detail, continue with the process.

12.4.5.3 Wear Pattern: The analyst compares the wear pattern in the impression with the test impressions and/or shoes and tires.

12.4.5.3.1 If the position and degree of general wear are different, document, and evaluate possible wear changes between the date of the crime (when impression was made) and the date that the shoes or tires were recovered. If the differences can't be explained, discontinue this process and report accordingly. This is an exclusion.

12.4.5.3.2 If the position and degree of general wear are similar and there is sufficient detail, continue with the process.

- 12.4.5.4 Randomly Acquired Characteristics: The analyst examines the impression for the presence of individual characteristics.
 - 12.4.5.4.1 Randomly Acquired characteristics can include bubbles, tears and other flaws imparted during the manufacturing process.
 - 12.4.5.4.2 Randomly Acquired characteristics can include cuts, nicks, tears, burns, feathering and chunk outs occurring during use of the shoe or tire.
 - 12.4.5.4.3 Randomly Acquired characteristics can include additions to the outsole or tire tread such as rocks, glass, wads of gum and thorns.
 - 12.4.5.4.4 These characteristics should be evaluated according to their position, size, shape, orientation, and clarity.
 - 12.4.5.4.5 Tire impression casts (or photos, if that is all that's available) are compared directly to the known unmounted tires. The impression/tire are compared for randomly acquired characteristics.
 - 12.4.5.4.6 When sufficient randomly acquired characteristics are present in the questioned impression and correspond with respective features in the known shoe or tire, an identification or match can be effected. Evaluate, document and report accordingly.
 - 12.4.5.4.7 Due to varying circumstances, not all randomly acquired characteristics will reproduce in every impression. Therefore, the absence of an individual characteristic is not a basis for exclusion and does not preclude identification.

12.4.6 Verification

- 12.4.6.1 Verification is the confirmation of an analyst's conclusion by another qualified analyst.
- 12.4.6.2 A qualified analyst shall verify all footwear and tire impression comparisons.
- 12.4.6.3 The verifier will complete the Footwear and Tire Impression Verification Worksheet and place it in the case file.
- 12.4.6.3 Analysts shall not verify any conclusions with which they are not comfortable. Comfort level is a function of training and experience.
- 12.4.6.4 Analysts are encouraged to work out differing conclusions through collaboration. If the differing conclusion(s) cannot be resolved, the ISP Quality Manual Section 15.9.4.3 "conflict resolution" policy will be followed.

12.5 DISTRIBUTION OF FOOTWEAR AND TIRES:

12.5.1 Background: Footwear and tires are produced, marketed, distributed and sold in stores for the purpose of generating profits at various levels. This information is tracked by various manufacturers, intermediaries and retail stores and can be obtained by forensic laboratories to assist with criminal cases.

12.5.1.1 Tire manufacturers can be contacted through *Who Makes it and Where Tire Directory* printed yearly, and *Who Retreads Tire*

12.5.1.2 When contacting the manufacturer, the request should be addressed to the public relations branch, sales department or engineering department. They may be able to supply a considerable amount of information regarding the shoe or tire in question:

- If the manufacturer has an exclusive contract for the molds producing the outsole or tire tread which is only sold with one brand name on it (may be several models with same outsole under the brand name) or if the mold produces an outsole or tire tread which may be distributed under several brand names.
- Number of shoe or tires produced from that mold in the year in question (production runs).
- Retail stores handling that design in the area of interest.
- Number of shoes or tires sold during the year in question.
- Samples of outsoles from the mold.
- Mold drawings of a particular tire.

12.5.1.3 Requests usually have to be followed up with a formal letter on letterhead stationary listing the request and general reasons for the request. The companies will not require specific case information, but may require information about the severity of the crime for an “expedite” request.

12.6 CONTROLS/SAFETY

12.6.1 Controls do not apply.

12.6.2 Contamination of footwear and tires by biological material is a possibility. Proper precautions should be taken such as gloves, lab coats and eye protection.

12.6.3 Shoe and tires may contain inclusions such as sharp rocks and glass, which cannot be removed prior to the comparison process. The analyst should exercise caution to prevent cuts.

12.6.4 Mounted tires may present a lifting hazard. The analyst should seek assistance with moving these, if necessary.

12.7 REFERENCES:

Bodziak, William J., *Footwear Impression Evidence*, Elsevier, New York, 1990, pp 144-148.

Nause, S/Sgt. Lawren, *Forensic Tire Impression Identification*, Canadian Police Research Centre, 2001.

The Scientific Working Group on Shoeprint and Tire Tread Evidence Technology (SWGTTREAD) - *SWGTTREAD documents are officially published in the Journal of Forensic Identification*, 766/55(6), 2005.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

13 FOOTWEAR AND TIRE SEARCH

- 13.1 The footwear and tire impression analysts may be asked to provide information about the type of shoe or tire that left the impression at a crime scene in cases where a suspect hasn't been identified yet. Generally, the analyst will work with a high quality photograph and/or cast of an impression that shows the finer design details. A complete shoe impression from heel to toe and a complete tire impression at least 12" long showing the design from inner to outer edge will make this task easier.
- 13.2 Shoe searches may be accomplished by:
- Using a commercially available shoe database software program,
 - Hand searching in stores,
 - Contacting manufacturers,
 - Internet searches of producer's collections of outsoles and manufacturer's and producer's images of outsoles.
 - Sending a copy of the impression to another lab such as the FBI Lab.
- 13.3 Tire searches may be accomplished by:
- Using a commercially available tire database software program such as Tread Assistant.
 - Hand searching the *Tread Design Guide*, available in the reference library.
 - Contacting local tire stores.
 - Contacting manufacturers
 - Sending a copy of the impression to another lab such as the FBI Lab.
- 13.4 SEARCH PROCEDURE
- 13.4.1 The analyst may work from a photograph or a sketch of the impression showing design details.
- 13.4.2 The actual evidence must not be taken from the lab to conduct the search.
- 13.4.3 Unless the impression is in the form of a cast, the analyst will have to do a negative to positive conversion to find the outsole or tread design.
- 13.4.4 Based on the footwear impression design elements, the analyst may be able to determine a category of shoe, such as hiking, work boot, cowboy boot, or athletic shoe to narrow the search.
- 13.4.5 Examination of a tire impression can give the analyst the number of ribs and grooves, whether there's a center rib and its appearance, and appearance of grooves, slots and sipes. This information will assist in the search for a possible tire(s).
- 13.4.6 When the analyst finds potential candidates, they should be printed for re-examination later.
- 13.4.7 The search candidates should account for all the design elements in the impression. If a design element cannot be accounted for, a knock-off shoe or manufacturer's design change may be possible. The analyst should consider contacting the manufacturer of the "close" candidate for assistance.

13.4.8 The analyst should be aware that new shoe or tire brands and models may be produced with the same design after the search is finalized.

13.5 **TIRE TRACK MEASUREMENT VEHICLE COMPUTER SEARCH:**

13.5.1 When measurements of wheel base and/or vehicle track width have been collected at the scene, the analyst can attempt to use resources such as the internet, tire manufacturers, or other forensic laboratories to develop a candidate list.

13.6 **DOCUMENTATION AND REPORTING:**

13.6.1 The search results should be documented in the notes and copies of any pictures of the shoes or tires should also be kept with the notes.

13.6.2 When a report is generated with the search result(s), a disclaimer should always be added that the list is not all-inclusive, unless the manufacturer has indicated that the list is complete.

13.6.3 Because of 13.4.8, a recommended start to the search report would be a general disclaimer "As of _____ date".

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016

14 CASE WORK DOCUMENTATION AND REPORT WRITING

- 14.1 Case work documentation and report writing will be according to ISPFS Quality Manual, Section 5.10 Reporting the Results.
- 14.2 Documentation shall be to the extent that another qualified analyst would be able to determine each examination activity conducted, their sequence, results of the activities, and any conclusions reached.
- 14.2.1 As each development method is completed it is noted in sequence on the EXAMINATION WORKSHEET and the evidence is visually examined for the presence of comparable impression detail.
- 14.2.2 When comparable impression detail is observed, it should be preserved prior to additional processing.
- 14.2.2.1 Comparable impression detail may be photographed upon initial examination, as additional detail develops, after a specific method, and/or prior to a subsequent method.
- 14.2.2.1.1 Impression photographs/images and/or case documentation shall include a scale, unique case identifier, analyst's initials.
- 14.2.2.2 Impressions developed via powder processing may be lifted in lieu of photography.
- 14.2.2.2.1 Impression lifts shall contain the unique case identifier, date, and analyst's initials.
- 14.2.3 Impression examination documentation shall include which impressions were analyzed, compared, evaluated, and the conclusions reached.
- 14.2.3.1 Each impression analyzed shall have an individualizing numeric designation.
- 14.2.3.2 The analyst will note any impressions deemed to not be of comparison quality.
- 14.2.3.3 Impression documentation may also include, at the analyst's discretion, the source of the impression and pattern if discernable.
- 14.2.3.4 Analysts shall document what the impressions were compared to and the results of those comparisons.
- 14.2.3.4.2 The verifying analyst shall complete the Footwear and Tire Impression Verification Worksheet.
- 14.2.3.4.3 Documentation of exclusions shall include, at a minimum, which impression(s) was excluded and the known shoes or tires the impression was excluded to.
- 14.2.3.4.4 Documentation of inconclusive findings shall include the reason(s) for the inconclusive finding. If the quality of the test impression(s) is at fault, additional test impressions should be made to improve the quality. If quality test impressions cannot be produced because the shoes or tires are severely damaged (for example, in a

fire), this problem would be listed as a reason for the inconclusive finding.

14.2.4 A reproduction suitable for comparison of both the compared impressions and the known test impressions must be retained as part of the case record.

14.2.4.1 When the laboratory cannot ensure that the original impression or test impressions used and relied upon in the examination, will be maintained by the contributing agencies, the laboratory must maintain an image of the actual data.

14.2.4.1.1 The analyst will produce shoe test impressions (except exclusion prints/photographs) from footwear submitted by the agency. Footwear that can be excluded based on outsole design differences does not require test impressions.

14.2.4.1.2 Case documentation shall contain copies of all impressions submitted by the customer. All impressions deemed of value for comparison shall be preserved in photographic or digital form.

14.2.4.1.3 Impression lifts produced and retained in the laboratory do not need to be copied for the case file.

14.2.4.1.4 Case documentation shall contain originals or machine copies of all known test impressions used in the comparison. Known test impressions submitted by the customer agency shall be copied prior to being returned. Copies must be suitable for comparison.

14.3 The report shall be as clear and concise as possible, convey the analytical findings and conclusions, and will be supported by scientific procedures.

14.3.1 The report is generally divided into three sections: EVIDENCE DESCRIPTION, EXAMINATION, and CONCLUSION. The following are some basic report wording guidelines categorized according to where they would appear in the report. There may be situations that do not fit the examples given. Unique wording for these situations will be developed as the need arises.

14.3.1.1 EVIDENCE DESCRIPTION:

Forensic Services Exhibit#, (Agency Exhibit #), a brief description of the evidence with an explanation if any unique identifiers were used to differentiate more than one item in the exhibit.

14.3.1.2 EXAMINATION:

14.3.1.2.1 Use for enhancement and development of impressions: Describe in general terms how the item(s) were examined (physical and/or chemical methods).

14.3.1.2.2 Use for retained evidence: Describe any impression lifts or photographic evidence that was marked and sealed for preservation and is being retained by ISP Forensic Services.

14.3.1.2.3 Use for footwear or tire search and vehicle track computer program: Describe what resources were used to do the search for make and model information.

14.3.1.2.4 Use for comparison: Describe which shoes and/or tires were compared to which impressions.

14.3.1.3 RANGE OF CONCLUSION (this exact wording need not be used):

14.3.1.3.1 Lacks sufficient detail – In the opinion of the examiner, insufficient detail was present in the questioned impression for a meaningful comparison with the known footwear or tire.

14.3.1.3.2 Exclusion – This is the highest degree of non-association expressed in footwear and tire impression examinations. Sufficient differences were noted in the comparison between characteristics in the questioned impression and the known footwear or tire. In the opinion of the examiner, the particular known footwear or tire was not the source of and did not make the impression.

14.3.1.3.3 Indications of non-association – In the opinion of the examiner, the questioned impression exhibit dissimilarities when compared to the known footwear or tire; however, certain details or features were not sufficiently clear to permit an exclusion.

14.3.1.3.4 Limited association of class characteristics - Some similar class characteristics were present; however, there were significant limiting factors in the questioned impression that did not permit a stronger association between the questioned impression and the known footwear or tire. These factors may include but were not limited to insufficient detail, lack of scale, improper position of scale, improper photographic techniques, distortion or significant lengths of time between the date of the occurrence and when the footwear or tires were recovered that

could account for a different degree of general wear. No confirmable differences were observed that could exclude the footwear or tire.

In the opinion of the examiner, factors (such as those listed above) have limited the conclusion to a general association of some class characteristics. Other shoes/tires with the same class characteristics are included in the population of possible sources of the impression.

14.3.1.3.5 Association of class characteristics - The class characteristics of both design and physical size must correspond between the questioned impression and the known footwear or tire. Correspondence of general wear may also be present.

In the opinion of the examiner, the known footwear or tire is a possible source of the questioned impression and therefore could have produced the impression. Other footwear or tires with the same class characteristics are included as possible sources.

14.3.1.3.6 High degree of association - The questioned impression and known footwear or tire must correspond in the observable class characteristics of the design, physical size, and general wear. For this degree of association there must also exist: (1) wear that, by virtue of its specific location, degree and orientation make it unusual and/or (2) one or more randomly acquired characteristics.

In the opinion of the examiner, the characteristics observed exhibit strong associations between the questioned impression and known footwear or tire; however, the quality and/or quantity were insufficient for an identification. Other footwear or tires with the same class characteristics are included as possible sources only if they display the same wear and/or randomly acquired characteristics observed in the questioned impression.

14.3.1.3.7 Identification – This is the highest degree of association expressed by a footwear and tire impression examiner. The questioned impression and the known footwear or tire share agreement of class

and randomly acquired characteristics of sufficient quality and quantity.

In the opinion of the examiner, the particular known footwear or tire was the source of, and made, the questioned impression and another tire or item of footwear being the source of the impression is considered negligible.

14.3.1.3.8 Footwear and tire search results, if done by ISP Forensic Services.

- 14.4 The case file should contain the following items if generated:
Report;
Examination documentation packet (notes, exam worksheets, copies/forms, digital imaging reports, photographs, test impressions);
Administrative documentation (copies of submission forms, communication logs, agency reports, etc.).
- 14.5 The verifying analyst shall perform the technical and administrative review on the case. In the event that the verifying analyst is unavailable, another approved analyst may perform the reviews.
- 14.6 Report writing and case documentation will be modified with the implementation of a Laboratory Information Management System (LIMS). Report writing and case documentation will follow ISPFS procedures.
- 14.7 REFERENCES

The Scientific Working Group on Shoeprint and Tire Tread Evidence Technology (SWG TREAD) - SWG TREAD documents are officially published in the *Journal of Forensic Identification*, 766/55(6), 2005.

Appendix A

Footwear and Tire Impression Instrument Maintenance

Manufacturer's instrument manuals

ALS

For ALS maintenance see method 8.1.

Balances

Footwear and tire impression analysts will use the balances located in the drug lab following the Controlled Substances Discipline "Balance Calibration Verification Analytical Method".

CAE Fuming Chambers

Maintenance shall consist of cleaning the tanks as needed.

Cameras

General maintenance consists of wiping camera bodies with a soft cloth, blowing off lenses and mirrors to remove dust or dirt, and then cleaning with a soft cloth or eyeglass cleaner.

If a camera malfunctions, it will be taken out of service until it can be repaired. The camera shall be tagged indicating that it is out of service. Maintenance, service, etc. will be recorded in the maintenance log.

ESDA

The bronze plate may become contaminated by fingerprints or small quantities of toner during use. The best method of cleaning this plate is to brush vigorously with tissue paper. If the plate becomes seriously contaminated or damaged such that the vacuum is impaired, the plate should be replaced. Additional maintenance shall be conducted as needed and all maintenance will be recorded in the maintenance log.

Electrostatic Dust Print Lifter

Recharge the battery overnight approximately every two months or as needed. Maintenance shall be conducted as needed and all maintenance will be recorded in the maintenance log.

Eyewashes and showers

Eyewash function shall be checked monthly and documented on the Chemistry Section QC worksheet. The showers are also checked annually by building maintenance.

Fume Hoods

All hoods are equipped with continuous flow monitoring devices. Capture velocity at the open face of the commercially purchased hoods is at least 100 feet/minute. Capture velocity at the open face of the sink hoods ranges between 75-100 feet/minute.

If a hood fails a monthly check, the check will be repeated. If the hood still fails, it will be taken out of service until it can be repaired. The hood shall be tagged indicating that it is out of service.

General maintenance consists of cleaning. Filters are changed regularly by building maintenance staff. Additional maintenance shall be conducted as needed and will be recorded in the maintenance log.

Fume hoods are certified yearly by an outside source.

Magnifying Glasses

Magnifying glasses should be cleaned regularly with a high quality lens cleaner and soft cloth. No caustic chemicals should be applied to the lens.

Microscopes

All microscopes have year servicing performed by an outside source.

Powder Station Exhaust Vents

Filters are washed periodically and replaced annually by the building maintenance staff. Additional maintenance is performed as needed. Maintenance, service calls, etc. will be recorded in the maintenance log.

Property of Idaho State Police Forensic Services
Uncontrolled Internet Copy
OBSOLETE DOCUMENT 11/4/2016